



## SARS-CoV-2 Antigen Rapid Test (COVID-19 Antigen Rapid Test) (Swab)

### Package Insert

REF INCP-502	English
--------------	---------

SARS-CoV-2 Antigen Rapid Test (Swab) is a rapid chromatographic immunoassay for the qualitative detection of SARS-CoV-2 Nucleocapsid Protein antigens present in swab specimen.

For professional *in vitro* diagnostic use only.

#### 【INTENDED USE】

The SARS-CoV-2 Antigen Rapid Test (Swab) is a rapid chromatographic immunoassay for the qualitative detection of SARS-CoV-2 Nucleocapsid protein antigens in swab specimen from individuals with suspected SARS-CoV-2 infection in conjunction with clinical presentation and the results of other laboratory tests.

Results are for the detection of SARS-CoV-2 Nucleocapsid protein Antigens. An antigen is generally detectable in upper respiratory specimens during the acute phase of infection. Positive results indicate the presence of viral antigens, but clinical correlation with patient history and other diagnostic information is necessary to determine infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. The agent detected may not be the definite cause of disease.

Negative results do not preclude SARS-CoV-2 infection and should not be used as the sole basis for treatment or patient management decisions. Negative results should be treated as presumptive and confirmed with a molecular assay, if necessary for patient management. Negative results should be considered in the context of a patient's recent exposures, history and the presence of clinical signs and symptoms consistent with COVID-19.

#### 【SUMMARY】

The novel coronaviruses belong to the  $\beta$  genus. COVID-19 is an acute respiratory infectious disease. People are generally susceptible. Currently, the patients infected by the novel coronavirus are the main source of infection; asymptomatic infected people can also be an infectious source. Based on the current epidemiological investigation, the incubation period is 1 to 14 days, mostly 3 to 7 days. The main manifestations include fever, fatigue and dry cough. Nasal congestion, runny nose, sore throat, myalgia and diarrhea are found in a few cases.

#### 【PRINCIPLE】

The SARS-CoV-2 Antigen Rapid Test (Swab) is a qualitative membrane-based immunoassay for the detection of SARS-CoV-2 Nucleocapsid protein Antigens in swab specimen. SARS-CoV-2 Nucleocapsid protein antibody is coated in the test line region. During testing, the specimen reacts with SARS-CoV-2 Nucleocapsid protein antibody-coated particles in the test. The mixture then migrates upward on the membrane by capillary action and reacts with the SARS-CoV-2 Nucleocapsid protein antibody in test line region. If the specimen contains SARS-CoV-2 Antigens, a colored line will appear in test line region as a result of this. If the specimen does not contain antigens to SARS-CoV-2, no colored line will appear in the test line region, indicating a negative result. To serve as a procedural control, a colored line will always appear in the control line region, indicating that the proper volume of specimen has been added and membrane wicking has occurred.

#### 【REAGENTS】

The test contains anti-SARS-CoV-2 Nucleocapsid protein antibody as the capture reagent and anti-SARS-CoV-2 Nucleocapsid protein antibody as the detection reagent.

#### 【PRECAUTIONS】

- This package insert must be read completely before performing the test. Failure to follow directions in package insert may yield inaccurate test results.
- For professional *in vitro* diagnostic use only. Do not use after expiration date.
- Do not eat, drink or smoke in the area where the specimens or kits are handled.
- Do not use test if pouch is damaged.
- Handle all specimens as if they contain infectious agents. Observe established precautions against microbiological hazards throughout in the collection, handling, storage and disposal of patient samples and used kit contents.
- Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are assayed.

- Wash hands thoroughly after handling.
- Please ensure that appropriate amounts of samples are used for testing. Too much or too little sample size may lead to deviation of results.
- Sterile Swabs for the collection of Nasopharyngeal specimen and Nasal specimen are different. Do not mix the using of the two types of sampling swabs.
- Extracted specimens for PCR tests cannot be used for the test.
- The used test should be discarded according to local regulations.
- Humidity and temperature can adversely affect results.

#### 【MATERIALS】

Material Provided		
• Test cassettes	• Sterile swabs	• Package insert
• Extraction buffer	• Extraction tubes and tips (Optional)	• Workstation
• Procedure card		

#### Materials Required But Not Provided

- Timer

#### 【STORAGE AND STABILITY】

Store as packaged in the sealed pouch at room temperature or refrigerated (2-30 °C). The test is stable through the expiration date printed on the sealed pouch. The test must remain in the sealed pouch until use. **DO NOT FREEZE.** Do not use beyond the expiration date.

#### 【SPECIMEN COLLECTION, TRANSPORT AND STORAGE】

##### Nasal Swab Specimen Collection

- Insert a sterilized swab less than one inch (about 2 cm) into a nostril (until resistance is met at the turbinates).
- Rotate the swab 5-10 times against the nasal wall. Using the same swab repeat the collection procedure with the second nostril.
- Withdraw the sterile swab; avoid excess volume and high-viscous nasal discharge.



##### Nasopharyngeal Swab Specimen Collection

- Insert a sterile swab into the nostril of the patient, reaching the surface of the posterior nasopharynx.
- Swab over the surface of the posterior nasopharynx 5-10 times.
- Withdraw the sterile swab from the nasal cavity and avoid excess volume and highly-viscous nasopharyngeal discharge.



Caution: If the swab stick breaks during specimen collection, repeat specimen collection with a new swab.

##### Specimen Transport and Storage

Specimens should be tested as soon as possible after collection.

If swabs are not been processed immediately, it is highly recommended the swab sample is placed into a dry, sterile and tightly sealed plastic tube for storage. The swab specimen in dry and sterile condition is stable for up to 24 hours at 2-8 °C.

#### 【SPECIMEN PREPARATION】

Only the extraction buffer and tubes provided in the kit is to be used for swab specimen preparation.

Please refer to the Procedure card for detailed information of Specimen Extraction.


- Place the swab specimen in the Extraction tube with Extraction buffer. Rotate the swab for **10-15 seconds** while pressing the head against the inside of the tube to release the antigen in the swab.
- Remove the swab while squeezing the swab head against the inside of the Extraction tube as you remove it to expel as much liquid as possible from the swab. Discard the swab in accordance with your biohazard waste disposal protocol.

\*NOTE: The storage of the specimen after extraction is stable for 2 hours at room temperature and 24 hours at 2-8 °C.

#### 【DIRECTIONS FOR USE】

Allow the test, extracted specimen and/or controls to equilibrate to room temperature (15-30 °C) prior to testing.

- Remove the test cassette from the sealed foil pouch and use it within one hour. Best results will be obtained if the test is performed immediately after opening the foil pouch.
- Invert the specimen extraction tube and add **3 drops of extracted specimen** (approx.75-100  $\mu$ L) to the sample well (S) and then start the timer.
- Wait for the colored line(s) to appear. Read the result at **15 minutes**. Do not interpret the result after 20 minutes.

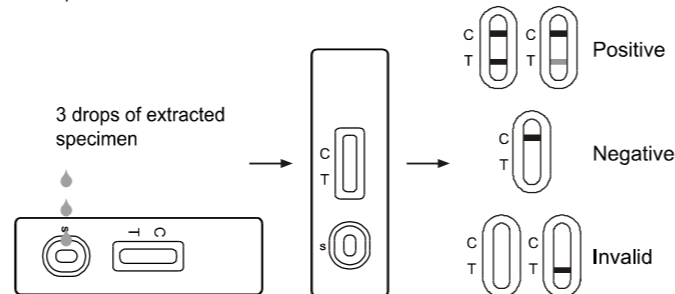

#### 【INTERPRETATION OF RESULTS】

(Please refer to the illustration above)

**POSITIVE:\*** **Two colored lines appear.** One colored line should be in the control region (C) and another colored line should be in the test region (T). Positive result in the test region indicates detection of SARS-CoV-2 antigens in the sample.

\*NOTE: The intensity of the color in the test region (T) will vary based on the amount of SARS-CoV-2 antigen present in the sample. So any shade of color in the test region (T) should be considered positive.

**NEGATIVE: One colored line appears in the control region (C).** No colored line appears in the test region (T).

**INVALID: Control line fails to appear.** Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

#### 【QUALITY CONTROL】

##### External Quality Control

Positive/negative controls are not included in this kit. However, in compliance with Good Laboratory Practice (GLP), these controls are recommended.<sup>1</sup>

##### Internal Quality Control

Internal procedural controls are included in the test. A colored line appearing in the control region (C) is an internal procedural control. It confirms sufficient specimen volume and correct procedural technique.

#### 【LIMITATIONS】

- The performance of the SARS-CoV-2 Antigen Rapid Test (Swab) was evaluated using the procedures provided in this product insert only. Modifications to these procedures may alter the performance of the test. Extracted specimens for PCR tests cannot be used for the test.

- The Direction for use and the Interpretation of test Result must be followed closely when testing for the presence of SARS-CoV-2 Nucleocapsid protein antigens in the human nasopharynx specimen from suspected individuals. For optimal test performance, proper sample collection is critical. Failure to follow the procedure may give inaccurate results.
- The SARS-CoV-2 Antigen Rapid Test (Swab) is for *in vitro* diagnostic use only. This test should be used for detection of SARS-CoV-2 Nucleocapsid protein Antigens in swab specimens as an aid in the diagnosis of patients with suspected SARS-CoV-2 infection in conjunction with clinical presentation and the results of other laboratory tests. Neither the quantitative value nor the rate of increase in the concentration of SARS-CoV-2 antigens can be determined by this qualitative test.
- The SARS-CoV-2 Antigen Rapid Test (Swab) will only indicate the presence of SARS-CoV-2 Antigens in the specimen and should not be used as the sole criteria for the diagnosis of SARS-CoV-2 infections.
- The results obtained with the test should be considered with other clinical findings from other laboratory tests and evaluations.
- If the test result is negative or non-reactive and clinical symptoms persist. It is recommended to re-sample the patient and test again or test with a molecular diagnostic device to rule out infection in these individuals.
- The test will show negative results under the following conditions:
  - The concentration of the novel coronavirus antigens in the sample is lower than the minimum detection limit of the test.
  - The optimal sampling time (peak virus concentration) after infection has not been verified, so collecting samples at different times for the same patient may avoid false negatives.
- Negative results do not rule out SARS-CoV-2 infection, particularly in those who have been in contact with the virus. Follow-up testing with a molecular diagnostic should be considered to rule out infection in these individuals.
- The accuracy of the test depends on the quality of the swab sample. False negatives may result from improper sample collection or storage.
- Positive results of SARS-CoV-2 may be due to infection with non-SARS-CoV-2 coronavirus strains or other interference factors.

#### 【PERFORMANCE CHARACTERISTICS】

##### Sensitivity, Specificity and Accuracy

The SARS-CoV-2 Antigen Rapid Test (Swab) has been evaluated with swab specimens obtained from the patients. RT-PCR (Nasopharyngeal Swab) is used as the reference method for the SARS-CoV-2 Antigen Rapid Test (Swab). Specimens were considered positive if RT-PCR (Nasopharyngeal Swab) indicated a positive result. Specimens were considered negative if RT-PCR (Nasopharyngeal Swab) indicated a negative result.

##### Nasopharyngeal Swab Specimen

SARS-CoV-2 Antigen Rapid Test (Swab)	RT-PCR (Nasopharyngeal Swab)		Total
	Positive	Negative	
SARS-CoV-2 Antigen	Positive	124	125
	Negative	3	318
<b>Total</b>	<b>127</b>	<b>316</b>	<b>443</b>
<b>Relative Sensitivity</b>	<b>97.6% (95%CI*: 93.3%~99.5%)</b>		
<b>Relative Specificity</b>	<b>99.7% (95%CI*: 98.3%~99.9%)</b>		
<b>Accuracy</b>	<b>99.1% (95%CI*: 97.7%~99.8%)</b>		

\*Confidence Intervals

##### Nasal Swab Specimen

SARS-CoV-2 Antigen Rapid Test (swab)	RT-PCR (Nasopharyngeal Swab)		Total
	Positive	Negative	
SARS-CoV-2 Antigen	Positive	604	605
	Negative	16	1092
<b>Total</b>	<b>620</b>	<b>1077</b>	<b>1697</b>
<b>Relative Sensitivity</b>	<b>97.4% (95%CI*: 95.8%~98.5%)</b>		
<b>Relative Specificity</b>	<b>99.9% (95%CI*: 99.5%~100%)</b>		
<b>Accuracy</b>	<b>99.0% (95%CI*: 98.4%~99.4%)</b>		

\*Confidence Intervals

#### Limitation of Detection

The SARS-CoV-2 Antigen Rapid Test (Swab) can detect out SARS-CoV-2 heat-inactivated virus strain as low as 1X10<sup>2</sup> TCID<sub>50</sub>/mL.

#### Specificity Testing with Various Viral Strains

The SARS-CoV-2 Antigen Rapid Test (Swab) was tested with the following viral strains. No discernible line at either of the test-line regions was observed at these concentrations:

Description	Test Level
Human coronavirus 229E	5 x 10 <sup>5</sup> TCID <sub>50</sub> /mL
Human coronavirus NL63	1 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
Human coronavirus OC43	1 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
MERS coronavirus Florida	1.17 x 10 <sup>4</sup> TCID <sub>50</sub> /mL
Human coronavirus HKU1	1 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
Influenza A H1N1	3.16 x 10 <sup>5</sup> TCID <sub>50</sub> /mL
Influenza A H3N2	1 x 10 <sup>5</sup> TCID <sub>50</sub> /mL
Influenza B	3.16 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
Parainfluenza virus 2	1.58 x 10 <sup>7</sup> TCID <sub>50</sub> /mL
Parainfluenza virus 3	1.58 x 10 <sup>8</sup> TCID <sub>50</sub> /mL
Respiratory syncytial virus	8.89 x 10 <sup>4</sup> TCID <sub>50</sub> /mL
Adenovirus type 3	3.16 x 10 <sup>4</sup> TCID <sub>50</sub> /mL
Adenovirus type 7	1.58 x 10 <sup>5</sup> TCID <sub>50</sub> /mL
Human Rhinovirus 2	2.81 x 10 <sup>4</sup> TCID <sub>50</sub> /mL
Human Rhinovirus 14	1.58 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
Human Rhinovirus 16	8.89 x 10 <sup>6</sup> TCID <sub>50</sub> /mL
Measles	1.58 x 10 <sup>4</sup> TCID <sub>50</sub> /mL
Mumps	1.58 x 10 <sup>4</sup> TCID <sub>50</sub> /mL

TCID<sub>50</sub> = Tissue Culture Infectious Dose is the dilution of virus that under the conditions of the assay can be expected to infect 50% of the culture vessels inoculated.

#### Specificity Testing with Various Organisms

The following organisms were tested at 1.0x10<sup>8</sup> org/mL and all found to be negative when tested with the SARS-CoV-2 Antigen Rapid Test (Swab):

<i>Arcanobacterium</i>	<i>Pseudomonas aeruginosa</i>
<i>Candida albicans</i>	<i>Staphylococcus aureus subsp. aureus</i>
<i>Corynebacterium</i>	<i>Staphylococcus epidermidis</i>
<i>Escherichia coli</i>	<i>Streptococcus pneumoniae</i>
<i>Moraxella catarrhalis</i>	<i>Streptococcus pyogenes</i>
<i>Neisseria lactamica</i>	<i>Streptococcus salivarius</i>
<i>Neisseria subflava</i>	<i>Streptococcus sp. group F</i>

#### Interfering Substances

The interfering substances below were spiked with negative, SARS-CoV-2 Antigen weak positive. No substances showed any interference with the SARS-CoV-2 Antigen Rapid Test (Swab).

Substance	Concentration
Whole Blood	20 $\mu$ L/mL
Mucin	50 $\mu$ g/mL
Budesonide Nasal Spray	200 $\mu$ L/mL
Dexamethasone	0.8 mg/mL
Flunisolide	6.8 ng/mL
Mupirocin	12 mg/mL
Oxymetazoline	0.6 mg/mL
Phenylephrine	12 mg/mL
Rebetol	4.5 $\mu$ g/mL
Relenza	282 ng/mL
Tamiflu	1.1 $\mu$ g/mL
Tobramycin	2.43 mg/mL

#### Precision Intra-Assay & Inter-Assay

Within-run and Between-run precision has been determined by using three specimens of SARS-CoV-2 standard control. Three different lots of SARS-CoV-2 Antigen Rapid Test (Swab) have been tested using negative, P1 and P5 specimens. Ten replicates of each level were tested each day for 3 consecutive days. The specimens were correctly identified >99% of the time.

#### 【BIBLIOGRAPHY】

- Westgard JO, Barry PL, Hunt MR, Groth T. A multi-rule Shewhart for quality control in clinical chemistry. *Clinical Chemistry* 1981;27:493-501

#### Index of Symbols

	For <i>in vitro</i> diagnostic use only		Tests per kit		Authorized representative in EU
	Store between 2-30°C		Use by		Do not reuse
	Do not use if package is damaged		Lot number		Catalog #
	Manufacturer		Consult Instructions For Use		

**Hangzhou AllTest Biotech Co., Ltd.**  
#550, Yinhai Street  
Hangzhou Economic & Technological Development Area  
Hangzhou, 310018 P.R. China  
Web: www.alltests.com.cn Email: info@alltests.com.cn



Statement: Information about manufacturer of sterile swab is placed on the packaging.

Number: 14602108300

Revision Date: 2024-02-07