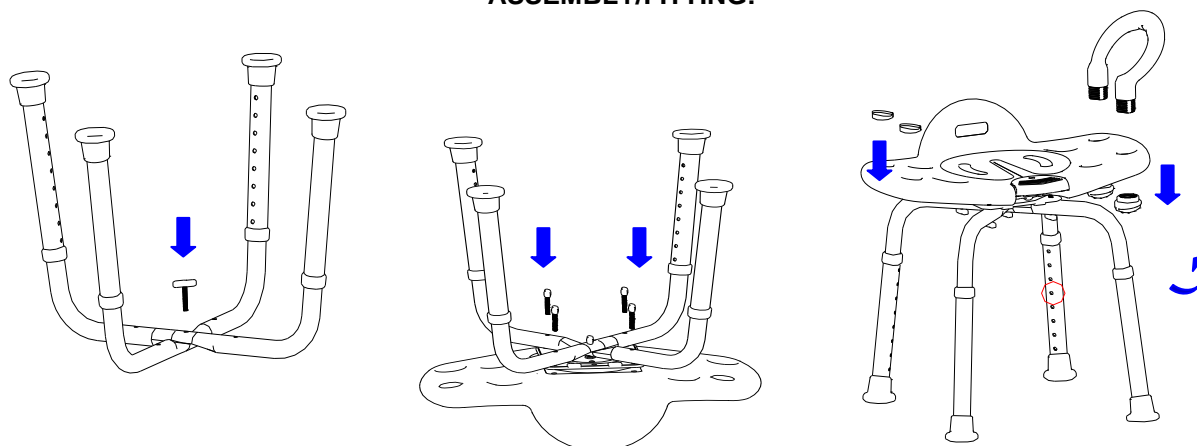


CODE 27733 SWIVEL SHOWER CHAIR USER'S MANUAL

⚠ WARNING: Be sure to read all information and directions before use.

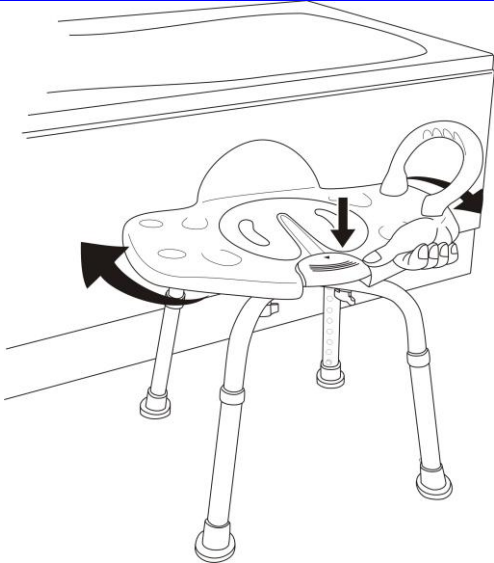
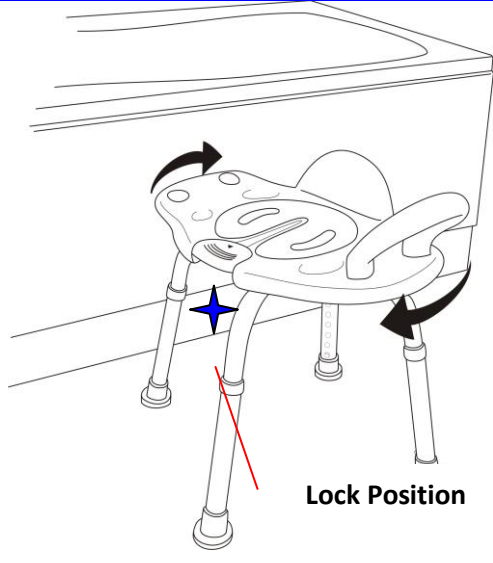
- ▲ **DO NOT** install this equipment without first reading and understanding this instruction sheet.
- ▲ **ALWAYS** verify that the device is secure prior to use.
- ▲ **ALWAYS** observe the weight limit on the labelling of your product.
- ▲ **ENSURE** to follow all instruction required while operating this device.

ASSEMBLY/FITTING:



- | | | |
|--|---|---|
| <p>1. Put two legs cross to form right angle and screw on the knob in the center hole.</p> | <p>2. To screw assembly other knobs with the swivel seat.</p> | <p>3. Adjust the seat height.
4. Select one end for handle screwed on, and covers up the small lids on the other end.</p> |
|--|---|---|

OPERATION:

- | | |
|--|---|
|  |  <p style="text-align: right;">Lock Position</p> |
| <p>1. Adjust the seat height, Sit on, Push down the front bottom to swivel the top seat.
2. Turn around right or left with one hand holding on the handle.</p> | <p>3. Hand off the front bottom.
4. Position lock at right angle location.</p> |