MINIASPEED

Battery Evo Plus


ISTRUZIONI D'USO

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INSTRUCTION MANUAL

MANUEL D'INSTRUCTIONS

MANUAL DE INSTRUCCIONES

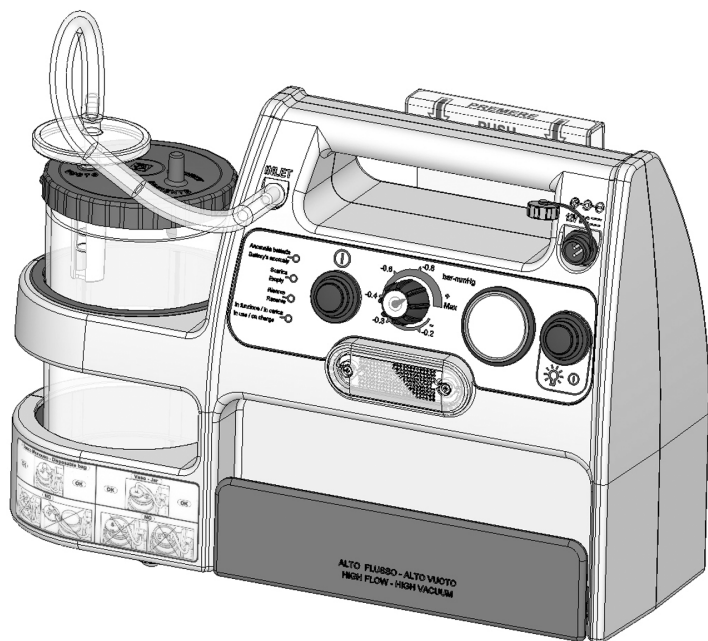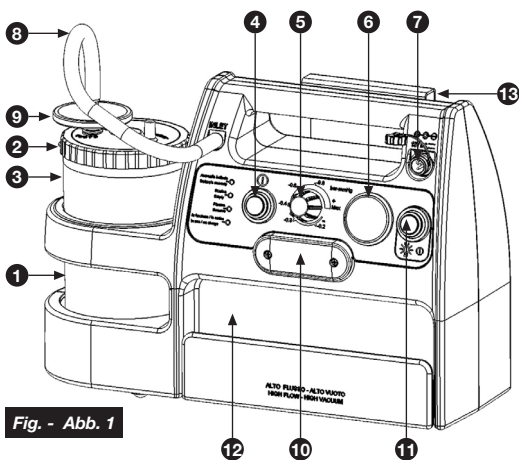
Fig. - Abb. 1

E

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

- GB**
- 1. Flashing red Led:** the battery is no longer capable of supplying maximum charge duration; replace it.
 - 2. Red Led permanently on:** the battery is flat, use the appliance for a maximum of 3 minutes and re-charge the battery
 - 3. Yellow Led permanently on:** the battery is in reserve (about 10 minutes' operating time left); recharge it if possible.
 - 4. Green Led permanently on:** appliance operating, battery charged.
Green Led flashing: battery under charge.

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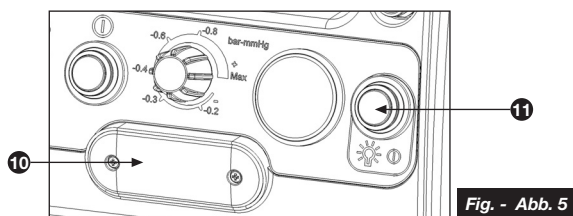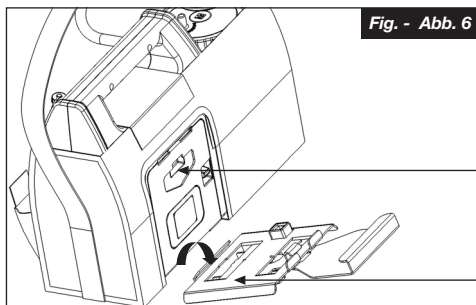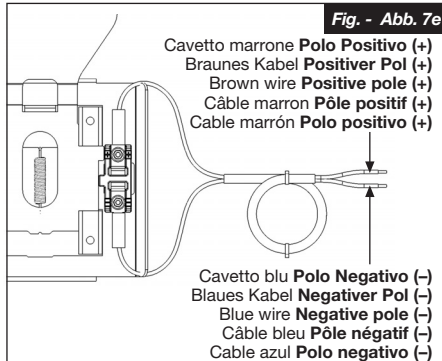
Fig. - Abb. 7e

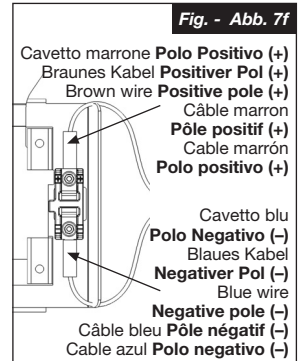
Fig. - Abb. 7f

The MINIASPEED BATTERY EVO PLUS aspirator, specific for secretion removal, is a portable professional device suitable for use in ambulances, in homes and in clinics. It has a vacuum regulator, a vacuum gauge and 1000 ml container with a protection device to prevent liquid from entering the suction pump, which cuts off the aspiration flow. It has lifetime lubrication and is easy to handle, simple to use, reliable, strong and silent. This model has a 12 V rechargeable battery; when connected to the ambulance mounting bracket it not only operates on the ambulance's battery but also recharges. The MINIASPEED BATTERY EVO PLUS aspirator has the following accessories: 1000 ml container with protection device, power supply cord with multi-voltage power supply unit, 24 cm long sterilisable silicone connection tube, 130 cm long sterilisable silicone connection tube, **sterile, disposable** cannula, **sterile, disposable** manual flow regulator, **disposable** bag and **disposable** antibacterial filter. **N.B.: Only use genuine accessories supplied by 3A Health Care; the accessories supplied with the device by the manufacturer have been tested and proven compliant with the latest current safety standards. WARNING! Correct operation of the device is not guaranteed in the event of use of accessories other than those supplied.**

INTENDED USE

Battery-powered medical / surgical aspirator intended for use in the field and/or in vehicles. It may also be used in home care and/or outpatient settings.

Medical purposes: This product is intended for use for the aspiration of body fluids.

Intended users of the product:

- Legally certified medical personnel (doctors, nurses and therapists).
- For home treatment or home care under the guidance of medical personnel.
- Users must also be able to understand the operation of the medical device, and the contents of the instruction manual, in general terms.

Intended patients for the product: People who need to remove body fluids (saliva, blood, etc.).

Environment: This product is intended for use in an ambulance. However, it may also be used in home care and/or outpatient settings.

Expected duration: Duration may vary based on the operating environment. The lifetime of the device is 5 years and that of the collection vessel and the silicone tubes 1 year or 30 sterilisation cycles. The cannula, manual flow regulator, antibacterial filter and bag are disposable devices and must therefore be replaced after each application. Frequent use of the product may shorten the duration.

Precautions for use: The warnings and cautions described in the instruction manual must be observed.

IMPORTANT WARNINGS

This is a medical device and must be used by qualified staff. It must be operated as indicated in this user instruction manual. It is important for the operator to read and understand the information on use and maintenance of the unit. If you have any questions, contact your stockist. MICROBIAL CONTAMINATION: in the case of illnesses with a risk of infection or microbial contamination, the accessories should be thoroughly cleaned and sterilised after each use.

The manufacturer has made every effort to ensure that all its products are of the highest quality and are reliable. Nevertheless, as for all electrical appliances, fundamental safety standards must be observed:

- Children and persons who are not self-sufficient may only use the unit under the strict supervision of a responsible adult who has read this manual.
- The device must always be used by specifically trained staff who have read this manual.
- Never use adapters for supply voltages different to the voltage shown on the data plate on the back of the unit. Keep the cord away from hot surfaces.
- The unit should not be used in the presence of inflammable anaesthetic mixtures with air, oxygen or nitrous oxide.
- Never handle the supply cord plug with wet hands or use the unit when taking a bath or a shower.
- Never leave the unit close to water, do not submerge it in any liquids. If by chance it falls into water, pull the plug out of the socket immediately before recovering it. Do not use the unit if the plug or supply cord are worn or wet (send it immediately to your stockist).
- Although the device is protected against sprays of water, do not use it in the rain.
- Only authorised personnel may perform maintenance and/or repair work. Unauthorised repairs annul the guarantee.
- Ensure that the connections and container closing are performed carefully to prevent suction losses.
- Do not tip the container over while connected to the unit when operating, as liquid may be sucked inside the appliance causing pump damage. If this happens, switch the aspirator off immediately and empty and clean the container. Send the appliance to your stockist.
- When the protection device intervenes aspiration is interrupted; empty the container and perform the cleaning operations.
- The cannula and the manual control of the aspirated flow are **sterile, disposable** products: They must be replaced after every application.
- Check the use-by-date on the original packaging of the cannula and check the integrity of the sterile packaging. If expired and/or deteriorated, replace it.

- The **disposable** antibacterial filter must be replaced after every application.
- Never use the battery charger with other appliances or for uses other than that established by this manual. Never use MINIASPEED BATTERY EVO PLUS with other power supply units.
- In view of their length, the power cord and connection hose could constitute a strangulation hazard.
- Use of the device in ambient conditions other than those specified in the manual may seriously impair its safety and technical characteristics.
- In the event of aspiration without the container and/or antibacterial filter, or if it is suspected that substances have entered the aspiration circuit, contact your stockist at once.
- Never modify the device's mounting bracket or the bracket installed in the ambulance as this may impair the safety of both the device and its installation environment.
- Always use the device in the vertical position on an unobstructed, stable, flat surface.

USE INSTRUCTIONS

Each time before use, ensure that all the accessories are perfectly clean according to the instructions indicated in the "CLEANING AND DISINFECTION OPERATIONS".

1. OPERATION WITH MULTI-VOLTAGE POWER SUPPLY UNIT

- 1.1 Check the charge level of the ambulance's battery before using the device
- 1.2 For use outside the ambulance (i.e. not fixed to the bracket [19]) with the internal battery flat or not fully charged, connect the multi-voltage power supply (16) to the socket on the device (7) and connect the cord (17) to a mains power socket.
- 1.3 Connect the appliance as indicated in Figure 1 (page 1).
- 1.4 Start the appliance by placing the switch to the "I" (ON) position (4) (green pilot light on).
- 1.5 The vacuum regulator (5) can be used to set the level of vacuum required (bar/KPa). Turn the knob clockwise, in the "+" direction, to obtain a higher vacuum or anticlockwise, in the "-" direction, for a lower vacuum; these values are read on the vacuum gauge (6),
Important: the vacuum values on the control decal are purely for guidance; always refer to the vacuum gauge reading.
- 1.6 After the application, switch off the device, remove the supply cord (17) from the mains socket and disconnect the power supply unit (16) from the device. Perform the cleaning operations as described in the "CLEANING AND DISINFECTION" paragraph.

2. OPERATION WITH BATTERY AND BATTERY CHARGER WITH MULTI-VOLTAGE POWER SUPPLY UNIT

- 2.1 The device is supplied with the battery partially charged. We therefore recommend charging it before use.
- 2.2 To charge the battery, with the device switched off connect the multi-voltage power supply unit (16) to the device's socket (7) and to the mains power supply by means of the cord (17). Charging time: about 4 hours. Charge duration: about 40/45 minutes at maximum suction power.
- 2.3 Operation with internal battery only:
Switch the device on by pressing the "ON-OFF button (4) (green light comes on). If the green light goes out during use, and the yellow light comes on, there is about 10/15 minutes of battery charge left (reserve level). Therefore: terminate the application if possible. If the application cannot be terminated, it can be continued until the red (battery flat) light comes on, but if this occurs do not use the appliance for more than 3 minutes to avoid damaging the battery. To continue the application with the battery flat (red light on), connect the multi-voltage power supply unit (16) (as described in point 2.2).
- 2.4 At the end of the application, switch the device off by pressing the button (4). Connect the power supply unit (16) to charge the battery (as described in point 2.2). If the device is not in use and is not connected to the bracket (19) in the ambulance, the power supply unit should be left connected (16) to ensure that the battery is always charged to the optimal level.
- 2.5 For aspiration procedures see points 1.3; 1.4; 1.5.
N.B.: when the battery power drops below a set level the aspirator shuts down to avoid damage to the battery. If absolutely necessary, the user can restart the aspirator for 1 minute by pressing the ON/OFF button again.

3. OPERATION WITH THE DEVICE INSTALLED IN THE AMBULANCE

- 3.1 The device can also operate connected to the ambulance mounting bracket (19).
- 3.2 With the device connected to the bracket (19), there is no limit to patient treatment time because the device is powered directly by the ambulance battery.
- 3.3 The device's battery is charged by leaving it connected to the ambulance mounting bracket (19).
3.3.1 The device's battery is also charged while treatment is in progress (with the device connected to the ambulance mounting bracket).
- 3.4 For aspiration procedures see points 1.3; 1.4; 1.5.

4. SECRETION COLLECTION VESSEL - 1000 ml

The 1000 ml collection bottle supplied with the aspirator can be used in two ways: as a collection vessel which can be sterilised as shown in Figure 3a or as a collection vessel with disposable bag (9) as illustrated in Figure 3b.

- 4.1 Sterilisable secretion vessel (1): the secretion vessel set consists of an overflow valve, a vessel (1) in clear material (polycarbonate) and a blue plug (2). Fit the antibacterial filter (6) straight into the plug (2); it will only

fit into the hole marked VACUUM/VUOTO, which is outlined in yellow. The antibacterial filter also protects the aspiration circuit from any contaminating agents sucked in during use. **Do not use the aspirator without the antibacterial filter, because from a bacteriological point of view, it becomes dangerous for the patient. Keep the device vertical to allow the overflow to function correctly.** All the components of the vessel can be sterilised using a conventional system in an autoclave at a temperature of 121° C, or by boiling for 10 minutes. We recommend replacing the complete vessel at every 30 sterilisation cycles. Do not overturn the vessel during use, in order to prevent the intervention of the non-return valve (3-4-5); should this occur, switch the aspirator off and detach the tube connected to the antibacterial filter. Never use the aspirator without the secretion collection vessel and/or without the antibacterial filter.

4.1.1 Connection: connect one end of the short, sterilisable silicon tube (7) to the antibacterial filter connector (6) and insert this latter into the “VACUUM” hole of the blue top (2); connect the other end to the “INLET” connector of the aspirator. Connect one end of the long sterilisable silicon tube (8) to the “PATIENT/PAZIENTE” connector of the blue top (2); to the other end, connect the disposable, sterile manual regulator (11) and connect the disposable, sterile cannula (10) to this latter.

- 4.2 Secretion collection vessel with disposable bag (9). The aspirator can be used with the re-usable 1000 ml. transparent secretion collection vessel (1) and with the disposable bag (9) supplied. In this case, the antibacterial filter is integrated in the disposable bag and therefore the antibacterial filter (6) and blue top with valve (2) are not used. The filter integrated in the bag also performs the function of preventing the return of aspirated liquids towards the aspirator. **In this case, to restore device operation, it is necessary to replace the disposable bag.** For the cleaning and disinfecting operations of the tubes (7-8) and vessel (1), sterilise the single parts in an autoclave at a maximum temperature of 121° C, or by boiling for 10 minutes. The bag is disposable and must always be replaced after every use. **The bag must be completely inserted in the vessel in order to prevent any vacuum losses.**

N.B.: maximum disposable bag usage vacuum: – 0.75 bar (75 kPa).

4.2.1 Connection: connect one end of the short sterilisable silicon tube (7) to the yellow connector (VACUUM) of the top (photograph 3b) and the other end to the “INLET” connector of the aspirator. Connect one end of the long sterilisable silicon tube to the red connector (PATIENT) and connect the disposable sterile manual regulator (11) and the disposable sterile cannula (10) to the other.

N.B.: only use the disposable bags supplied by 3A - Code 3A1687.

DESCRIPTION OF USE OF WORKING LIGHT

The MINIASPEED BATTERY EVO PLUS is fitted with a LED working light (10) which allows use of the device in dark conditions and/or at night.

Press the button (11) to switch on the light, as shown in Figure 5 of page 3; press the button (11) again to switch it off. The light switches off automatically after about 60 minutes.

The light is of LED type, so the amount of current it consumes is minimal and it can be used even during treatment; however, it should be switched off whenever it is not required and/or at the end of the treatment.

CLEANING AND DISINFECTION OPERATIONS

N.B.: If using chemical disinfectants, follow the manufacturer's instructions exactly.

- The cannula and the aspirated flow manual command are **sterile, disposable** products and must be replaced after every application.
- The antibacterial disposable filter must be replaced after every application.
- Never leave the appliance in water or submerged; clean the external casing of the appliance using only a damp cloth with detergent (non abrasive).

DESCRIPTION OF AMBULANCE MOUNTING / RELEASE SYSTEM

1. Fit the bottom of the bracket fixed to the aspirator into the recess provided in the bottom of the bracket mounted in the ambulance (as shown in Figure 6 of page 3).
2. Keeping the device engaged with the bracket, push it until it engages with the hook in the top of the yellow mobile bracket (as shown in Figure 6a):

N.B.: to ensure that the device engages correctly, press the mobile bracket, push the device fully down and release the mobile bracket (as shown in Figure 6b).

 **ALWAYS CHECK THAT THE DEVICE HAS ENGAGED WITH THE BRACKET.**

3. To release the device, press the yellow mobile bracket and lift the aspirator up and towards you (as shown in Figure 6c).

 **KEEP A FIRM GRIP ON THE HANDLE OF THE ASPIRATOR WHILE RELEASING IT.**

DESCRIPTION OF PROCEDURE FOR MOUNTING BRACKET IN AMBULANCE

The ambulance mounting bracket is supplied assembled and complete (as shown in Figure 7):

1. Wall mounting bracket.
2. Mobile bracket.
3. Mobile bracket spring.
4. Block complete with electrical contacts and cord.
5. 2 metre power cord for connecting the ambulance battery.

Mounting the bracket in the ambulance

- The bracket must be installed on the panel inside the ambulance by skilled staff authorised by the vehicle's constructor/outfitter.
- The fixing points are shown in Figure 7a, Figure 7b and Figure 7c.
- M5 Allen screws (UNI 5931-DIN 912-ISO 4762) should be used since they undergo deceleration testing at 10g as required by the EN 1789 standard. Otherwise, screws (with the same mechanical characteristics) having (non-tapered) heads, max. Ø10, and no more than 7 mm in height (including the washer, if fitted) can be used.
- See Figure 7d for the fixing centre distances.

Wiring connection/polarity

- The bracket is supplied complete with the ambulance battery connection cord (5) (already fitted), which powers the device and charges its battery when it is connected to the bracket.
- Connect the wires to the ambulance following the polarity shown in Figure 7e.
- If the cord is disconnected from the bracket for any reason, reconnect it following the polarity shown in Figure 7f.

N.B.: Take great care over the polarity (+ and -) because if it is reversed the device's electronic circuit board and battery will be damaged beyond repair.

PROBLEMS, CAUSES AND SOLUTIONS

| PROBLEMS | POSSIBLE CAUSES | SOLUTIONS |
|--|---|--|
| Excessive noise. | Damaged pump or blockages in the internal aspiration circuit. | Send to the assistance circuit. |
| The unit switches on but does not aspirate. | - Damaged pump. - Vacuum regulator fully open. Connection tubes disconnected and/or badly connected, broken connection tubes. Container not in a vertical position, full, or defective overflow valve. Possible blockage of the hydraulic circuit inside the unit. | - Send to the assistance circuit. - Check the position of the vacuum regulator. Check the connections and the integrity of the tubes. Position the container in a vertical position, check the overflow valve (blocked) and/or replace the silicon tubes. |
| The vacuum rate cannot be regulated. | Damage to the internal hydraulic system or blockage of the connection tubes to the aspiration unit. | Send to the assistance circuit. |
| When the appliance is switched on, the protection fuse always trips. | Pump probably damaged or in shortcircuit. | Send to the assistance circuit. |
| The vacuum gauge does not work. | Liquid penetrating the pneumatic circuit. | Send to the assistance circuit. |

Note: if you experience faults or malfunctioning problems different to those listed above, always and exclusively contact authorised assistance centres.

TECHNICAL CHARACTERISTICS

Professional portable medical device for secretion removal, suitable for fixed and transportable use. Complete with flame retardant ABS V0 casing and piston-cylinder electrical compressor lubricated for life.

| | |
|---|---|
| Risk class under Directive 93/42/EEC: | IIb |
| Suction class: | High vacuum / High flow |
| Protection level against entry of liquids and solids: | IP44 |
| Voltage: | 12VDC |
| Absorption: | 3A |
| Adjustable vacuum level: | from 0 to -0.85 bar (-85KPa) |
| Vacuum gauge precision class: | 2.5% under UNI EN 837 |
| Air flow: | 30 lt/min operating without restrictions +/- 10 % |
| Temporary use: | max. 45 minutes |
| Dimensions: | 38 (length) x 13.5 (width) x 26 (Height) cm |
| Weight: | 4.6 kg approx |
| Noise level: | 55dBA (1m) |
| Multi-voltage switching battery charger: | PRI: 100 ÷ 240V~ 50 – 60Hz SEC: 14VDC 4.28A |
| Internal battery: | 12VDC 4Ah Lead, hermetic |
| Battery charge duration: | 45 minutes at maximum suction power |
| Operating conditions: | |
| Temperature: | min. 0° C; max 40° C |
| Air humidity: | min. 10 %; max 95 % |
| Storage conditions: | |
| Temperature: | min. -10° C; max 50° C \nearrow |
| Air humidity: | min. 10 %; max 95 % RH |
| Operating-storage conditions: | min. 690 hPa; max 1060 hPa hPa |

SYMBOLS USED

| | |
|---|---|
|  Type BF appliance |  Alternate current |
|  It is compulsory to carefully read the instructions before using this device |  Do not use the unit when taking a bath or a shower |
| IP44 The device's casing is protected against solid particles having diameter of 1 mm or above, spraying water and access to hazardous parts with a wire |  0051 Conforms to Directive 93/42/EEC for medical devices |
|  Switch On |  The device contains a hermetic lead battery. It must be disposed of in accordance with current regulations on the disposal of toxic-harmful waste |
|  Switch off (or battery on charge) |  STERILE EO Sterilisation by ethyl oxide |
|  Direct current |  Disposable |

Electromagnetic Compatibility Compliance levels according to EN 60601-1-2:2015 standard

| |
|---|
| - ESD immunity: 15 kV air, 8 kV contact (EN 61000-4-2) |
| - Burst immunity: 2 kV/100 kHz (EN 61000-4-4) |
| - Surge immunity (EN 61000-4-5): 1 kV common mode /2 kV differential mode |
| - Magnetic field (EN 61000-4-8): 30 A/m |
| - Immunity to rf currents in the 150 kHz-80 MHz range (EN 61000-4-6) 3 V modulation 80% 1 kHz |
| - RF emissions, CISPR 11: Class B |
| - Harmonics emissions, EN 61000-3-2: Class A |

Rf field immunity (EN 61000-4-3):

| Field (V/m) | Frequency | Modulation |
|-------------|-----------------|--------------|
| 3 | 80MHz 2700MHz | 1kHz AM 80% |
| 27 | 380MHz 390MHz | 18Hz PM 50% |
| 28 | 430MHz 470MHz | 18Hz PM 50% |
| 9 | 704MHz 787MHz | 217Hz PM 50% |
| 28 | 800MHz 960MHz | 18Hz PM 50% |
| 28 | 1700MHz 1990MHz | 217Hz PM 50% |
| 28 | 2400MHz 2570MHz | 217Hz PM 50% |
| 9 | 5100MHz 5800MHz | 217Hz PM 50% |

Warnings:

Although compliant with the EN 60601-1-2 standard, the MINIASPEED BATTERY EVO PLUS medical device may interfere with other devices in the vicinity. The device must not be used in proximity to or stacked on top of other equipment. Install the device well away from other equipment that emits high frequencies (short waves, microwaves, electric scalpels, cell phones).

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are under control. The customer or user can help prevent electromagnetic interference by maintaining a minimum distances between mobile and portable RF communication equipment (transmitters) and the medical device as recommended below, according to the maximum output power of the radio communication equipment.

| Rated maximum output power of transmitter (W) | Separation distance (m) in relation to transmitter frequency | | |
|---|--|--|---|
| | from 150 kHz to 80 MHz $d = 1,2 \sqrt{P}$ | from 80 MHz to 800 MHz $d = 1,2 \sqrt{P}$ | from 800 MHz to 2,5 GHz $d = 2,3 \sqrt{P}$ |
| 0,01 | 0,12 | 0,12 | 0,23 |
| 0,1 | 0,38 | 0,38 | 0,73 |
| 1 | 1,2 | 1,2 | 2,3 |
| 10 | 3,8 | 3,8 | 7,3 |
| 100 | 12 | 12 | 23 |

For transmitters with rated maximum output power not listed above, the recommended separation distance d in metres (m) may be determined using the equation adopted for the transmitter frequency, where P is the maximum rated output power of the transmitter in Watts (W) stated by the transmitter manufacturer.

Notes:

- (1) At 80 MHz and 800 MHz the highest frequency range applies.
- (2) These guidelines might not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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DISPOSAL PROCEDURE (Dir. 2012/19/Ue-WEEE) The symbol on the bottom of the device indicates the separated collection of electric and electronic equipment. At the end of life of the device, do not dispose it as mixed solid municipal waste, but dispose it referring to a specific collection centre located in your area or returning it to the distributor, when buying a new device of the same type to be used with the same functions. This procedure of separated collection of electric and electronic devices is carried out forecasting a European environmental policy aiming at safeguarding, protecting and improving environment quality, as well as avoiding potential effects on human health due to the presence of hazardous substances in such equipment or to an improper use of the same or of parts of the same. **Caution!** The wrong disposal of electric and electronic equipment may involve sanctions.

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3A HEALTH CARE S.r.l.