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**EXTERNAL MALE CATHETER**

**CATHÉTER MASCULIN EXTERNE**

**CATÉTER EXTERNO MASCULINO**

**CATETER MASCULINO EXTERNO**

**URINALKONDOM**

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Mizra Medical Ltd  
21 Hasadna Street, Industrial Park, Holon, 5881411, Israel  
Made in Israel

RQMIS AREU SLU, Av Diagonal 409,  
Barcelona 08008, Spain

**Gima S.p.A.**  
Via Marconi, 1 - 20060 Gessate (MI) Italy  
gima@gimaitaly.com - export@gimaitaly.com  
[www.gimaitaly.com](http://www.gimaitaly.com)

Check Whether the package is in good condition before use, and confirm the package mark, production batch number or production date, and use within the validity period.

**Recommendation:**

- Do not apply the adhesive strip or reusable catheter fixing strap too tightly, to avoid circulation restriction.

**Storage conditions:**

- Before use, please store the device in the original packaging away from contaminated areas, dust, sunlight, excessive moisture and temperature and damaging chemicals. The product shall be stored in a ventilated, dry room away from heat source and free of corrosive gas. Always inspect product and conduct a user seal check before use as specified in the user instructions.

All serious accidents concerning the medical device supplied by us must be reported to the manufacturer and competent authority of the member state where your registered office is located.

**GIMA WARRANTY TERMS**

The Gima 12-month standard B2B warranty applies.

**ENGLISH**

**INSTRUCTIONS**

external male catheter is intended for collection of urine in males suffering from incontinence. MEC and related accessories are compliant with Regulation (EU) 2017/745  
Please read the instruction before using the product

**Description:**

The non-invasive non sterile, Male External Catheter, are intended for use to divert urine in male patients with urinary incontinence. This catheter is only for external and single use. This product contains Latex.

**Intended Use:**

This catheter is only for external and single use. Diverts urine from penis to extension tubing and urinary collection bags.

**Indications for use:**

Male external Catheters (MEC) is intended to be designed to safely conduct urine through funnel tube, keeping fluid pathway open.

**Applying the Catheter:**

1. Clean and dry the penis, trimming excess pubic hair.
2. Hold the catheter away from the head of the penis approx. 2 cm. away from the head of the penis.
3. Gently unroll the catheter along the length of the penis.
4. Secure the catheter can be done by using adhesive strip or the reusable catheter fixing strap to prevent slippage of the catheter.
5. Make sure that the adhesive strip or reusable catheter fixing strap is not applied too tightly to avoid constriction of circulation.
6. Attach the urine bag to the catheter.

**Removal of the catheter**

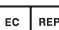
1. Disconnect the catheter from the urine bag.
2. Gently unroll the catheter away from the penis.

**Warning:**

- Male External Catheter is made of medical grade natural latex. Do not use the product if you are sensitive to latex. In case of suspected latex sensitivity stop using the catheter and refer to your doctor.
- Do not reuse the catheter.
- The Catheter must be Changed every 12 to 24 hours.

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