

EN - DISPOSABLE MOUTHPIECES: INSTRUCTIONS FOR USE

INTENDED USE: The disposable mouthpieces are intended for use together with pulmonary function testing equipment and allow the patient to take the exam without having direct contact with the equipment. The patient personally manipulates the mouthpiece by inserting it into the equipment. The mouthpieces should be applied only to intact, clean and healthy skin.

PRODUCT IDENTIFICATION: The identification data (reference, lot and date of production) is on the label attached to the box.

USAGE:

- Always follow the instructions of the spirometer.
- Insert the mouthpiece into the device sensor (turbine or pneumotacometer) without forcing. Never re-use the mouthpiece to avoid exposing the patient to cross-contamination.
- Do not open before use.
- Do not use if the packaging is damaged.
- The device is not suitable for use in neonatal paediatrics.
- Store in a cool and dry place.

DISPOSAL AND SAFETY: After use, the product must be disposed of in normal containers of non-toxic solid waste or in general hospital waste.




CONTRAINDICATIONS: Based on the experience and results of different analysis, there are no adverse effects as a result of proper use of the product.

STORAGE: Store in a cool and dry place.

WARNING:

- Packaging in individual heat-sealed envelope (single-wrap) guarantees the absence of bacteria, which is regularly monitored according to the manufacturer's quality system. If the medical device is reused, there is a risk of cross-infection for the patient who reuses it.
- If the device is reused, there are risks of cross-infection for the patient who re-uses it.

Manuale d'uso LUMED Rev. 4.0 01/09/2023

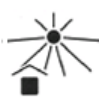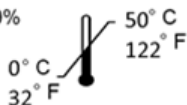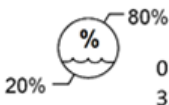

LUMED Srl
Via Staffora 18/9
20073 Opera ITALIA
www.lumed.com
tel +39 (0)2 57606750