



GIMA

PROFESSIONAL MEDICAL PRODUCTS

GIMA 39792

Fabbricante / Manufacturer / Fabricant / Fabricante:
Medistar Instruments Co.
New Miana Pura East, Ror
as Road, Sialkot - Pakistan
Shams@medistarinst.com
Made in Pakistan

REF MI-78951

EC REP Lorcan & Fyon B.V.
Europalaan, 40, 3526KS Utrecht, Netherlands

Importato da / Imported by / Importé par / Importado por:

Gima S.p.A.
Via Marconi, 1 - 20060 Gessate (MI) Italy
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com



ENGLISH

PEDICURE FILE PLASTIC


Detailed Description of Device Name:

Pedicure File Plastic, Double Sided, 23.5cm
Finish: White Plastic Color.

Instructions of Use

Pedicure File is made of high impact Polystyrene plastic. It is used for removing rough and excess hard skin (rubbing purpose).

Chemical Composition

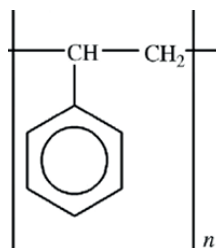

Indication of Use



- Put away from the reach of children.
- Device is used by the qualified professionals.

Cleaning/disinfection and storage conditions

Medistar's powder / paper coated instruments are not autoclavable, but they can be cleaned / disinfected with disinfectants like Endo Star (high level instruments disinfectant) or with any other better disinfectant/sterilizing solution. 20ml Endo Star in 1 liter tap water (2% dilution) requires 30 minutes immersion and 30ml Endo Star in 1 liter tap water (3% dilution) requires 15 minutes immersion.

In general, U.S. Environmental Protection Agency (EPA) registered disinfectants require 10 minute immersion; isopropyl and ethyl alcohol require 5 minutes.

Repeated processing has minimal effect on the instrument life. End of useful life for metal surgical instruments is normally determined by wear and damage due to the intended surgical use.

Store the instruments in an area protected from light sources and moisture, keep away from sunlight, keep dry, between a temperature from -10°C and 40°C.

REF	Product code	LOT	Lot number
	Keep in a cool, dry place		Keep away from sunlight
	Caution: read instructions (warnings) carefully		Consult instructions for use
	Temperature limit		Non-sterile
EC REP	Authorized representative in the European community		Date of manufacture

GIMA WARRANTY TERMS

The Gima 12-month standard B2B warranty applies.

