English Instructions for Use

Vent Filters and Sets.

Intended use: Filters, hydrophobic and non-hydrophobic for use within oxygen concentrators and other medical equipment such as ventilators. Filters are typically inserted into the gas sampling line to prevent a measuring equipment and it's system from becoming contaminated to reduce the risk of patient cross-contamination.

CAUTION: A CAUTION statement provides important information about a potentially hazardous situation which, if not avoided, may result in harm to the patient.

Any use of the medical device requires full understanding and strict observation of all portions of these Instructions for Use. The medical device may only be used for the purpose specified under "Intended Use". Observe all CAUTION statements throughout this manual and all statements on medical device labels. The manufacturer disclaims any liability for patient injury caused by noncompliance with these statements.

Following installation of the filter, check that all connections are secure.


The medical device IFU is not available for each individual filter. Only one copy of the Instructions for Use is included in the distribution package and should therefore be kept in an accessible location for users.


Medical device is only to be used by healthcare professionals.

User has the responsibility to report any serious incident that has occurred in relation to the medical device, to the manufacturer, EC REP and the competent authority of the Member State in which the user and/or the patient is established.

To avoid contamination and soiling, the product should remain packaged until ready to be used.

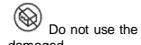Do not use the product if the packaging and/or product is damaged.

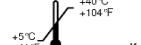If product was kept outside the stated storage temperature, product functionality cannot be guaranteed. Storage Temperature: +5°C to 40°C.

Store the product at room temperature, in dry place, free of dust and not exposed to direct sunlight.

Handle filter with care. Abusive handling can damage the filter.


For Your Safety and that of Your Patients, follow the Instruction for Use of the machine!

Before installation, check that all system components are free of obstructions and of foreign bodies. Otherwise, use is limited, or faulty operation is possible.


Inappropriate connection to other devices can result in a patient hazard.

Medical device products have been designed, tested and manufactured exclusively for disposable use and for a period of use not exceeding 24 hours following the first opening of the primary packaging.

Do not use after expiry date stated on primary packaging.

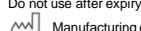Manufacturing date

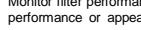Expiry date

Monitor filter performance during the use. If any changes to filter performance or appearance are observed, change the filter to new one.

Always keep spare filter ready to exchange with the one in use, if any of the below are observed:

- Increase in resistance;
- Filter is visibly soiled or damaged or if it rattles;
- Excessive condensation;
- Associated alarm.


The products must not be reused, cleaned or sterilised / not to be sterilised after use in order to re-use. Exceptions from this rule will be stated within the Instruction For Use for the particular products code if applicable. Reuse, cleaning or sterilization may lead to failure of the products, risk infection and injury to the patient.


This product is for single patient use only.

This product is latex and PHT free.

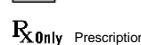Batch code symbol.

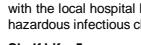Prescription only.

Following use, the products must be disposed off in accordance with the local hospital hygiene and waste disposal regulations for hazardous infectious clinical waste.

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