
Vasofix® Safety

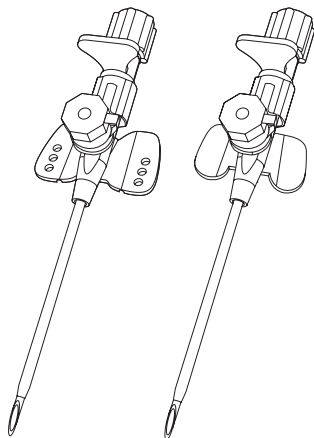B | BRAUN

LLDorder 5647 - G 221048

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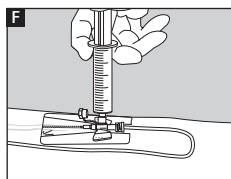
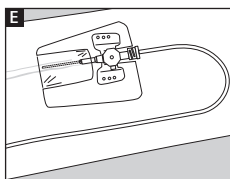
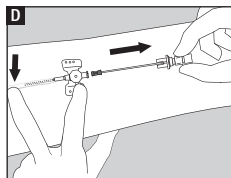
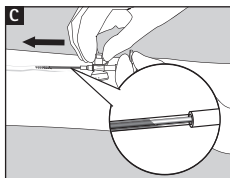
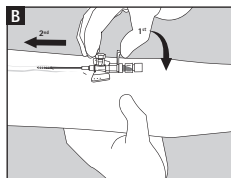
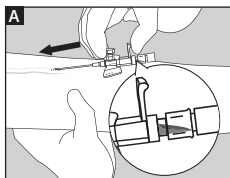
en Instructions for use

Description

The Vasofix® Safety IV Catheter with injection port is an indwelling peripheral catheter to obtain intravenous access for short-term use. It is designed to:

- Minimize inadvertent needlestick injuries with a passive safety shield
- Allow quick needle-free bolus injections via injection port
- Improve catheter stability and minimize catheter movement with fixation wings.

The Vasofix® Safety IV Catheter can be used with Luer access devices which are in compliance with ISO 80369-7.

This device is intended for single use only. It is provided sterile and non-pyrogenic.

Materials Used

- PUR, PP, PE, Silicone rubber, ABS, Chrome-Nickel steel.
- FEP, PP, PE, Silicone rubber, ABS, Chrome-Nickel steel.

Components are not made with natural rubber latex, PVC, or DEHP.

Indication

The Vasofix® Safety IV Catheter is inserted into a patient's peripheral venous system for short-term use to

- sample venous blood, or
- administer infusion solutions, blood products or pharmaceuticals (e.g. in or without carrier solutions) according to the SPC of the drug/ solution.

Patient population & Intended User

The catheter may be used for any patient population with consideration given to adequacy of vascular anatomy and appropriateness for the procedure being performed, the solution being infused and duration of therapy.

The catheter is intended to be used by qualified Healthcare Professionals who are trained on the use of peripheral IV catheters according to national regulations.

Contraindications

- The Vasofix® Safety IV Catheter should not be used in patients with known hypersensitivity to any of the materials employed.
- The Vasofix® Safety IV Catheter shall not be used for arterial access, subcutaneous access and high pressure power injection applications.

Precautions & Warnings

- Observe standard precautions on all patients. Aseptic technique, proper skin preparation, and continued protection of the insertion site are essential.
- When placing or maintaining any IV Catheter, standard precautions must be adhered in accordance with Centers for Disease Control and Prevention/ Occupational Safety and Health Administration (CDC/ OSHA) standards for blood-borne pathogens to avoid the risk of exposure to contaminated blood.

- Use only if packaging is intact. This device is sterile unless packaging has been opened or damaged.
- Do not re-use. Re-use of single-use devices creates a potential risk of patient or user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.
- In the case of an unsuccessful IV Catheter placement attempt, remove the needle first to activate safety mechanism, then remove catheter from patient and discard both.
- Never reinsert the needle inside the catheter once the needle has been partially or completely withdrawn as it may pierce and/or sever the catheter.
- Do not attempt to bypass the safety mechanism.
- In the unlikely event that the safety mechanism has not been triggered/activated, keep the tip of the needle away from the body and fingers at all times and immediately dispose the IV Catheter into a suitable sharp-proof container. Care must be taken to avoid needlestick injuries.
- Report needlestick injuries immediately and follow established institutional protocols.
- Use of only ISO Luer slip and Luer lock connectors ensures proper functioning.
- Luer Slip connections should not be left unattended due to potential of disconnection.
- Extreme care should be taken not to damage, pierce, cut or sever the catheter. Therefore, do not bend the catheter and/or needle during insertion, advancement, or removal of the needle.
- Do not use scissors or sharp instruments at or near the insertion site.
- If blood is present, flush device according to institutional protocol.
- Monitor tightness of all connections during the whole infusion procedure.
- The catheters may be used with access pressure of up to 3 bar.
- The catheters may be used to facilitate the placement of vascular access devices such as guidewires, indwelling central venous catheters (CVC), peripherally inserted central catheters (PICC) and midline catheters into the vascular system.

Residual risks/ Side effects

- General risks related to IV Catheters: infiltration, extravasation, leakage, air embolism, phlebitis, thrombophlebitis, thrombosis, catheter-related bloodstream infection, local infection, inflammation, catheter occlusion, catheter fracture, catheter embolism or clot formation.

Application

Use aseptic technique. Do not rotate catheter hub prior to insertion.

1. Gather all necessary items for insertion and stabilization procedure.
2. Select and disinfect per institutional protocol.
3. Apply tourniquet.
4. Remove protective guard in a straight outward motion.
5. Anchor the vessel with gentle skin traction and access vessel. Observe blood return in flashback chamber to confirm successful vessel entry (see figure A).
6. Lower and advance the entire catheter and needle unit slightly to ensure catheter tip is in vessel (see figure B).
7. Using the injection port, advance catheter off the needle (approximately 1/8" or 3mm). You can observe blood return between the needle and catheter to confirm catheter is in the vessel (see figure C). After confirmation, continue advancing catheter off the needle and into vessel.
8. Release tourniquet.
9. Before removing the needle, occlude vein distal to catheter tip to prevent blood spillage. At the same time, stabilize the catheter hub to prevent catheter dislodgement during needle removal. Withdraw the needle straight out with a controlled and continuous

motion (minimize rotation or bending of the needle). The metal safety shield will automatically attach to the needle tip as the needle tip exits the catheter hub (see figure D).

10. Dispose of the needle immediately into a sharps-proof container.
11. Immediately connect the infusion line or accessory device and cover the puncture site with a sterile and transparent dressing (see figure E) per institutional protocol.
12. For bolus injections open the injection port and attach syringe to injection port with a slight twisting movement (see figure F). The port opens and closes automatically. After injection immediately close protective cap of injection port.
13. After removal, dispose catheter according to local guidelines and/or institutional protocols.

Duration of Use

- Duration of use is dependent on the applied therapy according to the SPC of the drug/solutions and should follow the national guidelines and/or hospital protocols.
- The insertion site should be checked frequently and at regular intervals. The catheter should be removed in event of local or systemic signs of infection.

Disposal Disposal

Disposal according to local guidelines and/ or institutional protocols.

If, during the use of this product or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

en	Do not re-use
de	Nicht wiederverwenden
bg	Да не се използва повторно
cs	Nepoužívat opětovně
da	Må ikke genbruges
el	Να μην επαναχρησιμοποιείται
es	No reutilizable
et	Mitte korduvkasutada
fi	Ei saa käyttää uudelleen
fr	Ne pas réutiliser
hr	Nemojte ponovno upotrebljavati
hu	Ne használja fel újra
id	Jangan gunakan kembali
it	Non riutilizzare
kk	Қайтадан қолдануға тыйым салынған
lt	Negalima naudoti kartotinai
lv	Neizmantot atkārtoti
nl	Niet opnieuw gebruiken
no	Skal ikke gjenbrukes
pl	Nie używać ponownie
pt	Não reutilizar
ro	A nu se reutiliza
ru	Не использовать повторно
sk	Opätovne nepoužívať
sl	Samo za enkratno uporabo
sr	Ne koristiti ponovo
sv	Får ej återanvändas
tr	Yeniden kullanmayınız
uk	Повторно використовувати заборонено
vi	Không tái sử dụng
zh	不得二次使用

Consult instruction for use
Gebruiksaanwijzing beachten
Вижте инструкциите за употреба
Čtěte návod k použití
Læs brugsanvisningen
Συμβουλευτείτε τις οδηγίες χρήσης
Consulte las instrucciones de uso
Vt kasutusjuhendit
Katso käyttöohje
Consulter les instructions d'utilisation
Pogledajte upute za uporabu
Nézze meg a használati utasítást
Baca petunjuk penggunaan
Consultare le istruzioni per l'uso
Пайдалану нұсқаулығын қараңыз
Žr. naudojimo instrukcijas
Lasīt lietošanas instrukciju
Raadpleeg gebruikersinformatie
Se i bruksanvisningen
Przed użyciem należy zapoznać się z instrukcją
Consulte as instruções de utilização
Consultați instrucțiunile de utilizare
См. руководство по эксплуатации
Pozri návod na použitie
Glejte navodila za uporabo
Pogledajte uputstvo za upotrebu
Se bruksanvisningen
Kullanma kılavuzuna bakınız
Ознайомитись з інструкцією для застосування
Tham khảo hướng dẫn sử dụng
查阅使用说明

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B. Braun Medical AG
B. Braun Melsungen AG