- **I.V.STAND ON 5 WHEELS TROLLEY**
- **ASTA PORTAFLEBO SU CARRELLO A 5 ROTELLE**
- **POTENCE PORTE-SÉRUM SUR CHARIOT 5 ROULETTES**
- **SOPORTE PARA IV SOBRE BASE DE 5 RUEDAS**
- **SUPORTE DE SORO COM CARRINHO DE 5 RODAS**
- **INFUSIONSHALTER AUF 5 RÄDERN**
- **ΣΤΑΤΩ ΟΡΟΥ ΣΕ ΤΡΟΧΗΛΑΤΗ ΒΑΣΗ 5 ΤΡΟΧΩΝ**
- **STOJAK DO INFUZJI NA WÓZKU Z 5 KÓŁKAMI**
- **INFUZNÍ STOJAN NA 5 KOLECH VOZÍK**
- **DROPPSTÄLL PÅ EN FEMHJULS BAS**
- **STATIV PERFUZIE CU 5 ROȚI**
- **INFÚZNY STOJAN NA 5-KÓLESOVOM VOZÍKU**
- **INFUUSSTANDAARD OP 5 WIELEN, WAGEN**
- **I.V.-STATIV MED 5 HJUL**
- **IV-TELINEN 5-PYÖRÄISELLÄ VAUNULLA**
- **I.V. KOLICA S POSTOLJEM NA PET TOČKOVA**
- **INFÚZIÓS ÁLLVÁNY 5 KERÉKES ALAPON**
- **INTRAVENINIŲ INFUZIJŲ STOVAS SU 5 RATŲ VEŽIMĖLIU**
- **I.V. STATĪVS UZ 5 RITEŅU RATIŅIEM**
- **STOJALO ZA INFUZIJO NA VOZIČKU S 5 KOLESI**
- **KÄEPIDEME JA RIILIGA VIIE RATTAGA TILGUTI STATIIV**
- **I.V. САМОСТОЯЩА КОЛИЧКА С 5 КОЛЕЛА**

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GIMA 27815



 **Jiangsu Saikang Medical Equipment Co., Ltd**
No. 35 Lehong Road, Modern Agriculture Demonstration
Park, Zhangjiagang City, Jiangsu Province, China
Made in China

REF SKH041(13)

EC **REP** EVERBIZ GmbH LANDSBERGER STR. 155
80687 MÜNCHEN - GERMANY

Gima S.p.A.
Via Marconi, 1 - 20060 Gessate (MI) Italy
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

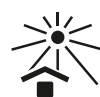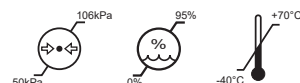INTENDED USE

Intravenous (IV) stands or poles are medical devices that are designed to hang bags containing intravenous fluids or medicines which need to be administered to the patient.

1. INSTALLING GUIDE

1. Install the base with casters.

- Connect the 5 feet of the base with casters to the plastic base connector with screws and fix.

2. Install the lower rod to the base.

- Connect the lower rod to the 5-feet base by screws and fix.

3. Install the plastic platform to the lower rod.


- Insert the plastic platform to the lower rod.

4. Install the upper rod to the plastic platform.

- Connect the upper rod to the plastic platform by screws and fix.

5. Install the 4 hooks on the top.


- Insert the 4 hooks on the top of the upper rod.

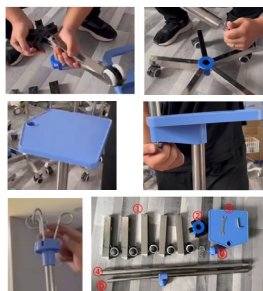2. USE METHOD

2.1 Use method of height adjustment

Adjust the knob, adjust the infusion stand to the appropriate height, check whether the screws are tight, ready to use.

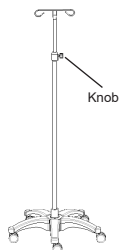3. MAINTENANCE

1. In order to use the IV pole safely, it is necessary to perform regular safety inspections on the IV pole. It is recommended to check it every six months to ensure that the connections are not loose and all operations are normal.
2. When the IV pole frame ages and it reaches a certain servicing time, which the metal part of the IV pole and plastic parts of the stainless steel basin frame can be recycle.
3. Avoid scratching the panel with sharp-angle appliances or knives during use, and clean it frequently to keep it clean and dry.
4. If the panel is accidentally stained with stains, it is recommended to clean it. Do not use alkaline or corrosive chemicals to clean the IV pole, which will cause the stainless steel surface to rust.
5. Regularly check the casters to avoid hitting or overloading, which will cause damage to the casters.
6. It is forbidden to force the casters to brake when braking, which will cause damage to the brake system.
7. The wear of the tire surface of the wheel can be detected by visual inspection. Some thread spun yarn and other accumulations may be wrapped around the wheel, remove the bolts and nuts of the wheel, clean up the debris, and check whether the bearing of the wheel is damaged. If the parts are not damaged, can be reassembled and used.
8. Judging whether the caster is replaced according to the time used, wear and tear, you can replace the wheels and make sure that the axles bolt and nuts are well put together and connected. If possible, use a lock washers or locknuts to help you with this process.
9. For the casters equipped with brakes, it is necessary to regularly check whether the brakes are working properly. Check the brakes once a day or before each use. For a IV pole equipped with multiple brake casters, only one brake caster can be locked at a time which will allow you to try and push the IV pole also check whether one of the brake performance of each caster is good. If the function of the brake fails due to wear or damage to the wheels, immediately replace the wheels and then retest the brake system.
10. If the mechanism of the brake system of the caster is damaged, and the brake needs to be repaired or replaced, please contact our after-sales service department or authorized dealer provide by our company. Every time the brakes are replaced, the brake performance of the casters must be retested.

4. PACKAGING, TRANSPORT, STORAGE

1. The packaging of stainless steel I.V. pole is carried out according to the contract or product standard.
2. Avoid rushing out, violent vibration, and prevent from sunlight and rain during the transportation of the stainless steel I.V. pole.
3. The IV stand should be stored:
 - a) Ambient temperature: $-40^{\circ} \sim +70^{\circ}$.
 - b) Relative humidity: $\leq 95\%$.
 - c) Atmospheric pressure: $500 \text{ hPa} \sim 1060 \text{ hPa}$.

5. MAX LOAD

Hooks: max load for each hook: 3kg

Pump Tray: max load 10kg

Handle: max load 10kg

Overall Structure: max load 32kg

6. AFTER-SALES SERVICE

1. Please keep the files that comes with the machine and invoices of this product properly, and you need to present these files when the company performs warranty and maintenance for the product.
2. If there is any problem in the process of using, please contact our company in time, so that our company can provide you with accurate and fast technical support and maintenance services in a timely manner.
3. From the date of sale, if the product is broken or damaged due to the correct installation and use according to the specification, the product will enjoy one-year free warranty and lifelong maintenance service with the "certificate" or invoice.
4. Do not disassemble the internal device of this product by yourself to avoid unnecessary damage. If you find that the use is affected due to quality problems, please contact our after-sales service department or our authorized dealer.
5. From the date of purchase, if it is indeed damaged or does not work normally due to quality problems within one year, the company will provide to the users a free repairing of the product.
6. Lifetime service from manufacturer: JIANGSU SAIKANG MEDICAL EQUIPMENT CO LTD

GIMA WARRANTY TERMS

The Gima 12-month standard B2B warranty applies

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