

READ ALL INSTRUCTIONS PRIOR TO USE. CONSULT PHYSICIAN BEFORE USE.

CONTRAINDICATIONS: Not intended for persons over 242 lbs (110 kg) - Not intended for icy or slippery conditions. - Not intended for persons with latex allergies.

INSTRUCTION GUIDE: To unfold your cane hold on to the handle and let the foldable shaft segments unfold towards the floor. To prevent falls, you should always check to make sure that the individual shaft segments (Fig. 1) are securely locked into place.

WARNING DO NOT disassemble the body of the cane. The cane contains an elastic cord inside. Uncontrolled release the elastic cord inside can cause severe injury to unprotected body parts, particularly eyes, if the elastic cord is cut, frayed or damaged discontinue use. **DO NOT** stretch the cord more than is needed to disassemble or assemble the segments.

TO ADJUST THE HEIGHT OF YOUR CANE: Loosen the twist lock by rotating it clockwise several turns. Depress the locking button and pull up on the handle to lengthen, or push down on the handle to shorten. Remember to retighten the twist lock after adjusting the cane to your desired height.

WARNING To prevent falls, always ensure individual shaft segments are securely locked into place and locking button is fully engaged.

Always make sure the pivoting foot of the stick is securely screwed to the cane to avoid falls.

TO COLLAPSE THE CANE: Your cane is easily foldable into a convenient travel size by separating the 4 individual shaft segments from each other. Hold each set of shaft segments near the joint and pull them apart. The individual shaft segments will neatly fold together



TO USE THE LIGHTS ON YOUR CANE: The button on the top of the cane handle activates the lights. Pressing the button once will activate the light in front of the cane, pressing it again will turn on also the red light at the base of the handle, pressing it a third time will turn them off. The front light can be rotated to direct the light where you need

SOS SIGNAL: To activate the SOS signal, press the button on the side of the handle, the light at the base will flash and the stick will emit an alarm sound.

BATTERIES REPLACEMENT: Open the battery compartment at the back of the handle. Remove any dead batteries and properly dispose of them. Using 2 AAA batteries, insert them with the correct polarity (positive pole to the outside where the spring is located). Close the battery compartment cover. Regularly inspect cane for signs of wear, including any fraying to the elastic cord or wear to the rubber feet. If you see signs of wear, discontinue use immediately.

WARNING Do not mix old and new batteries. Do not mix alkaline, standard (carbon zinc) or rechargeable (nickel - cadmium) batteries.

GIMA WARRANTY TERMS: The Gima 12-month standard B2B warranty applies.

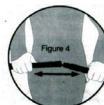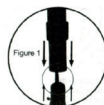

Gima 43050 - 43052

REF GB-806(1FOOT) GB-806(SKID PROOF)

Fabricato da/Manufactured by/Fabrique par/Fabricado por/Hergestellt von/
Fabricado por/Κατασκευάστηκε από/Произведено от/Vyrobce/Fremstillet af/
Tootja/Vaalistaja/Proizvodat/Gyártó/Ražots/Gamintojas/Produsert av/
Gefabriceerd door/Wyprodukowane przez/Fabricat de/Vyrába/Proizvajalec/
Tilbered av/تصنيع بواسطة
Ningbo Gaobo Sports Co., Ltd
No.1, Xiajin, Dali, Shenzhen Town, Ninghai County
Ningbo 315514, China
Made in China

Importato da/Imported by/importé par/importado por/importiert von/
Importado por/Εισαγυγί από/Довезено от/Dovezeno od/Importat de/
af/Importija/ Tuoja:/Uvezeno od strane/Importálta:/importēja/importavo/
Importert av/Geimporteerd door/importowane przez/importat de/Dovážal/
Uvozi/importerad av/توريد بواسطة
Gima S.p.A.
Via Marconi, 1 - 20060 Gessate (MI) Italy
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

