


ENGLISH

The Contec Elizabethan Collar is made of durable plastic, providing protection for your pet during recovery. Its sturdy design helps prevent scratching, biting, or chewing, ensuring a safe healing process.

How to use

1.

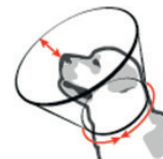


Measure the circumference of the pet's neck and select the appropriate Elizabeth ring based on the size chart. Ensure that after wearing, the pet's nose does not exceed the edge of the circle to effectively prevent licking and biting wounds.

2.

Gently trap the collar around the pet's neck, pointing it in the correct direction (with the wide end facing forward), ensuring that it does not affect the breathing or diet of pets.

3. F


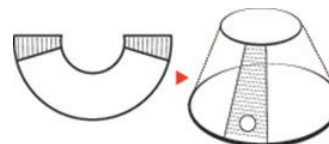


after adjusting the collar to the appropriate tightness, the Velcro is firmly adhered. After wearing, ensure that the pet can move freely, but cannot easily remove or break free.

4.

When worn for the first time, pets may have uncomfortable reactions. Give the pet some time to adapt, while observing its movements to ensure that wearing the loop does not affect normal activities.

Clean PVC Elizabethan collar: Wipe with mild soapy water, rinse with damp cloth, dry in shade. No soaking, harsh chemicals, or sun exposure.

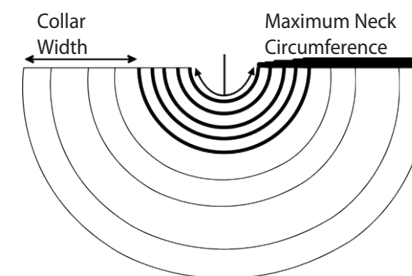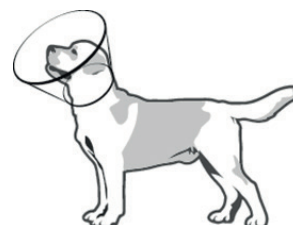

Measure with a ruler at home

Specification	Suitable for neck	depth	Recommended weight
Extra Large	49-56cm	28cm	Within 40kg
No. 1	43-50cm	22cm	Within 30 kg
No. 2	41-48cm	20.5cm	Within 15 kg
No.3	36-42cm	18cm	Within 10 kg
No. 4	31-36cm	13.5cm	Within 7.5 kg
No. 5	26-29cm	13cm	Within 5 kg
No. 6	19-23cm	11cm	Within 4 kg
No.7	17-20cm	9cm	Within 3 kg

How to clean

Clean PVC Elizabethan collar: Wipe with mild soapy water, rinse with damp cloth, dry in shade. No soaking, harsh chemicals, or sun exposure.

REF	VET0007
	VET0006
	VET0005
	VET0004
	VET0003
	VET0002
	VET0001



Contec Medical Systems Co.,Ltd.
Address: No.112 Qinhuang West Street,
Economic & Technical Development Zone, Qinhuangdao, Hebei Province
PEOPLE'S REPUBLIC OF CHINA
Tel: +86-335-8015430
Fax: +86-335-8015588
E-mail: cms@contecmed.com.cn
Website: <http://www.contecmed.com>
Made in China