



PROFESSIONAL MEDICAL PRODUCTS

## DISPOSABLE STANDARD SMMS GOWN

## ENGLISH

**Product Name** Disposable Standard SMMS Gown

**Product Code** FGNS-SG-SM40-M, FGNS-SG-SM40-L,

FGNS-SG-SM40-XL, FGNS-SG-SM40-XXL

**Product Size** Unit: cm

GIMA CODE	21040	21041	21042	21043
FULL GUARD CODE	FGNS-SG-SM40-M	FGNS-SG-SM40-L	FGNS-SG-SM40-XL	FGNS-SG-SM40-XXL
Sizes	120x140 (M)	120x150 (L)	130x150 (XL)	140x160 (XXL)
Length (Shoulder to hem)	120	120	130	140
Width (monolayer)	140	150	150	160
Belt length	65	65	65	69
Sleeve Length	58	59	59	61
Deviation	±2	±2	±2	±2

Product Model Gown Type.M/L/XL/XXL

**Structure and Composition** The product uses 40g SMMS non-woven fabric, It is composed of front,back,sleeve,ties and cuff. the technology process is by cutting, ultrasonic welding, ties, velcro, cuff, folding,sealing.

**Intended use** The products are mainly used for general isolation in outpatient clinics, wards, inspection rooms, cleaning room, laboratory etc. It is intended to be used by healthcare professionals in hospitals and clinics to maintain a good level of hygiene and contribute to infection control.

**Contraindication** The products can be used only once and should be used immediately after opening the package. If the package is damaged,the product is wet or moldy. Don't use them.

**How to use the product**

Open the package and take out the product, gently shake the gown, lift the sides of the collar, insert the hands into the sleeves, extend the arms forward and tie the belts.

**Storage** It should be stored in a clean room with relative humidity not exceeding 80%, no corrosive gas, stored in a cool, dry and well ventilated.

**Warning**

1. This product is disposal product, don't reuse.
2. This product is single use, never share your product with another.
3. It is forbidden to use when the packaging or the products are damaged.

**Notice: Any serious incident that has occurred in relation to the device should be reported to the manufacturer and competent authority of the Member State in which the user or patient is established.**













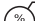



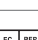



**Validity** 3 years.

**Date of production** See date on packaging

**Date of expiration** See date on packaging


**GIMA WARRANTY TERMS**


The Gima 12-month standard B2B warranty applies.

 REF	Product code		Keep away from sunlight		Keep in a cool, dry place
 LOT	Lot number		Medical Device compliant with Regulation (EU) 2017/745		Caution: read instructions (warnings) carefully
	Manufacturer		Date of manufacture		Medical Device
	Non-sterile		Expiration date		Don't use if package is damaged
	Humidity limit		Consult instructions for use		Disposable device, do not re-use
	Temperature limit		Authorized representative in the European community		Fragile, handle with care
	This side up		Unique device identifier		

GIMA 21040  
GIMA 21041  
GIMA 21042  
GIMA 21043

 Xianning Full Guard Medical Products Co., Ltd  
Yongan East Road Wenquan Economic  
Development Zone, Xianning City, Hubei  
Province, China  
Made in China

 REF FGNS-SG-SM40-M  
FGNS-SG-SM40-L  
FGNS-SG-SM40-XL  
FGNS-SG-SM40-XXL

 EC REP Humiss Beratung  
GmbH Gneisenaustraße 8.40477,  
Düsseldorf, Deutschland

Importato da / Imported by:

Gima S.p.A.

Via Marconi, 1 - 20060 Gessate (MI) Italy

gima@gimaitaly.com - export@gimaitaly.com

[www.gimaitaly.com](http://www.gimaitaly.com)

