

CHIAPPA ARREDAMENTI S.r.l.

Via Monte Rosa, 15
20014 Nerviano (MI) – ITALY
T: +39-0331-587504

E: web@chiappaarredamenti.com
W: www.chiappaarredamenti.com

Mod. 120765 SLB RD
EXTENDABLE BARIATRIC BED - ELECTRIC FOUR SECTIONS
BACKREST, LEGREST, HEIGHT, TR & RTR ADJUSTMENT



<mattress platform from 200cm x 90cm to 200cm x 120cm>

Technical Specification

Bed Upper Frame Dimensions	
1. Fully Retracted	230 cm x 106 cm (side rails included)
2. Fully Extended	246 cm x 136 cm (side rails included)
Mattress Platform Adjustments	
1. Fully Retracted	200 cm x 90 cm
2. Intermediate Position	200 cm x 105 cm
3. Fully Extended	200 cm x 120 cm
4. Foot Extension	from 200 cm to 220 cm
Highest Position	75 cm (79 cm with integrated scale)
Lowest Position	39 cm (43 cm with integrated scale)
Bed Weight	175 kg (195 kg with integrated scale)
Electrical Shock Protection	Class I
Liquid Ingress Protection (electric set)	IP 66
Maximum User Weight	450 kg
Safe Working Load	500kg

Product Features

Four section electrically profiling bed frame with heavy steel construction for years of service
Trendelenburg & Reverse Trendelenburg 12°
Cardiac Chair
Mattress Platform Positioning <ul style="list-style-type: none"> • Backrest 0° - 52° (64° with accessory Backrest Auto-Regression: mod. 500937) • Kneebreak 0° - 32°
Floor clearance 16,5 cm
Manual Handset with Locking Key + Supervisor
Electric CPR function
Armchair position with excellent comfort and ergonomic side rails configuration
Carer controlled lock out function for all controls
Battery backup
x3 standard twin castors Ø 150 mm and x1 Ø 150 mm Anti Static and steering twin castor
NEW Central Braking System with a three position pedal at foot side for a greater control: steering / swivelling / braked
x4 NEW multi position and removable solid side rails included colour co-ordinated HANDLES to bed frame (complying with 60601-2-52)
x2 NEW EXPANDABLE & easily removable head/foot boards with HPL (light blue, blue, yellow, ivory or green colour available) (complying with 60601-2-52)
x4 corner bumpers integrated into the frame
x2 sockets for I.V. Pole at head side
ANTIMICROBIAL White Epoxy Powder Coated



<backrest, legrest, height 39-75cm (43-79cm with integrated scale), cardiac chair>



<Trendelenburg >



<Reverse Trendelenburg>



<integrated foot extension from 200cm to 220cm >



<blanket tray >



steering



swivelling



braked

< NEW Central Braking System with a three position pedal on foot side for a greater control >

Accessories Not Included

Heavy Duty Lifting Pole with integrated IV hooks (mod. 500933)



I.V. Pole Fixed Height (mod. 500699)




I.V. Pole Additional Two Sockets at Foot Side (mod. 500698)



Additional Antistatic Castor (mod. 120864)

Backrest Auto-Regression (mod. 500937)

Backrest X-ray Cassette Holder (mod. 500953)

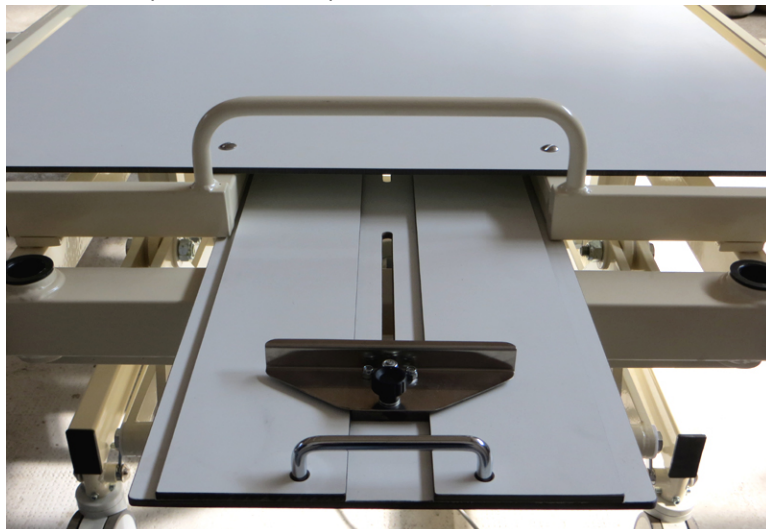

Bariatric Mattress (mod. 500453)

BARIATRIC BED AVAILABLE ALSO with **INTEGRATED SCALE**

Mod. 120765 SLB RDS



The **EXTENDABLE BARIATRIC BED** provided with **digital in-bed Scale** allows the staff to weigh patients as they lay there on the bed without using a standing scale. The staff can quickly gather data from many patients, without wasted motion or patient discomfort!
Now you can set up an efficient system for monitoring the weight.


Load Capacity: 500kg

Division: 200g

Accuracy Class III according to Nawi

(Non Automatic Weighing Instrument)

European Directive 90/384/EEC.

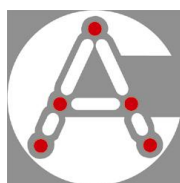

EC DIRECTIVES 93/42 EEC – MEDICAL DEVICES CLASS I
and further modifications 2007/47/EC
CEI EN 60601-2-52

Italian Ministry of Health RDM REGISTRATION # 106290/R
CND REGISTRATION #: z12030799
GMDN REGISTRATION #: 35563

RDM **106290-R**

00000001



::: made in italy :::