



# GIMA

PROFESSIONAL MEDICAL PRODUCTS

**CINGHIE MONITOR FETALE SINGOLO PAZIENTE**  
**SINGLE PATIENT FOETAL TRANSDUCER BELTS**  
**SANGLES POUR TRANSDUCTEUR FOETAL À USAGE UNIQUE**  
**CORREAS PACIENTE INDIVIDUAL PARA Sonda FETAL**  
**CORREIA FETAL TRANSDUTORA DESCARTÁVEL**  
**EINWEG-GÜRTEL FÜR FETALEN SCHALLKOPF**  
**ΖΩΝΗ ΜΟΡΦΟΤΡΟΠΕΑ ΕΜΒΡΥΪΚΟΥ ΕΛΕΓΧΟΥ ΜΙΑΣ**  
**ΧΡΗΣΗΣ**  
**EENMALIGE FOETALE TRANSDUCER RIEM**  
**PASEK PRZETWORNIKA PŁODU JEDNORAZOWEGO**  
**UZYTUKU**  
**CUREA ABDOMINALĂ PENTRU MONITORIZARE FE-**  
**TALĂ, DE UNICĂ FOLOSINȚĂ**  
**FOSTERGIVAREBÄLTE FÖR ENGÅNGSBRUK**

**REF B01 (GIMA 29544)**


Fabbricante / Manufacturer / Fabricant / Fabricante / Fabricante / Hersteller /  
Κατασκευαστής / Fabrikant / Producent / Producător / Tillverkare:  
Dongguan Belcare Medical Technology Co., Ltd.  
405, Building 5, #8 of Qingfeng South Road, Tangxia Town,  
Dongguan, Guangdong, China 523710.  
Made in China

Importato da/Imported by/ Importé par/ Importiert von/ Importado  
por/ Importado por/ Εισάγεται από/ Geïmporteerd door /  
Importowane przez / Importat de / Importerad av:

**Gima S.p.A.**

Via Marconi, 1 - 20060 Gessate (MI) Italy  
gima@gimaitaly.com - export@gimaitaly.com

[www.gimaitaly.com](http://www.gimaitaly.com)

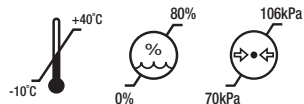

### Intended use:

- This product works according to adjust buckle's place on the belt to fix the fetal transducer.

### Product introduction:

- Product name: Disposable fetal transducer belt.
- Product structure: It's made up of elastic belt and buckle.

### Product range:

- It used with fetal transducer when pregnant in fetal monitoring, to fix the fetal transducer.

### Use condition:

- Temperature: 0°C~40°C
- Relative humidity: 15%~95%
- Atmospheric pressure: 70kpa~106kpa

### Package and storage:

- Disposable fetal belt should be single packaged, packaged product should be stored in a non-corrosive gas, good ventilation room.

### Store environment as below:

- Temperature: -10°C~+40°C
- Relative humidity: ≤80%
- Atmospheric pressure: 70kpa~106kpa

### Installation/use:

- Attach one end of the belt to the transducer by engaging one of the holes in the belt over the button on the top of the transducer.
- Keeping the transducer in position, tension the other end of the belt and engage the belt over the button, ensuring sufficient tension to keep the transducer in firm contact with the abdomen.
- Avoid over-tightening as this will cause unnecessary discomfort to the patient.

### Product clean:

- Wash belt in warm water, do not rub the edge fetal belt to avoid removal, do not exposure.

### Validity period:

- Validity period: 2 years.

### Replace and abandon:

- Replace: Under these situations such as fracture, elastic damage, expiry date, please change new fetal belt to use.
- Abandon: Abandoned belt can be processed as normal waste.

### Contraindications, precautions, warnings and cautionary statements:

- Contraindications: The belts must not be fixed at the injury site.
- Spandex allergies to stop using.
- Cautionary statements: This product is a single-use products to avoid cross infection. Please read the instructions before use.

### GIMA WARRANTY TERMS

The Gima 12-month standard B2B warranty applies

