

BODY FAT ANALYZER WITH CARRYING CASE

• 27322 GIMA BODY FAT ANALYZER

Ideal for dieticians, gymnasias and sports centres, it is perfect also for homecare. It will help manage nutrition, fitness and weight control and can indicate if you are at the right weight and also whether or not you are nutritionally healthy.

Advanced technology in your hands.

GIMA Body Fat Analyzer measures the flow of electrical signals (total body electrical resistance or impedance) as they pass through FAT, LEAN (muscle, bones, vital organs) or WATER. When the amount of fat, lean or water changes, so does the Body Fat signal, giving you a highly reliable and accurate measure of the amount of each of these components that make up your total weight. Specifically, 4 electrodes are placed over metacarpus and metatarsus where a 50 kHz current is introduced. The detected signal level will be an index of total body water and fat free mass.

Only 2 minutes for an essential and easy test

Testing and results take less than two minutes and have the same accuracy as other more expensive devices.

A lot of meaningful data

FAT: actual weight of fat on the body

LEAN: actual weight of lean, mostly muscle and vital organs

WATER: water content in the body

% FAT: percentage of fat weight in total body weight

% LEAN: percentage of lean weight in total body weight

BMR: basal metabolic rate

Target Weight: weight goal under optimal condition

Target % Fat: normal percent for sex and age

Portable, for performing test everywhere

Small size and light weight for portability.

Body Fat Analyzer helps to get and stay fit.

It helps to maximize personal level of fitness if you want to lose weight, build muscle, tone and trim, or simply be sure you are maintaining your current level of fitness:

- enables you to check effectiveness of weight loss program, telling exactly where your weight loss is coming from (fat-lean-water).
- it helps monitoring improvements in muscle mass as a result of exercise and/or diet.

- establishes a baseline for your body composition, so that nutritional and fitness improvements can be regularly monitored. Body fat is sold with a carrying case, kit of 2 sensor cable, sensor pads (box of 20), 2x9V batteries, transformer and manual (English, Italian).

• 27349 GIMA BODY FAT ANALYZER + SOFTWARE

Same as 27322 but provided of software to transfer and save data into PC. Compatible with Windows 98, XP, Vista.

• 27321 SPARE SENSOR CABLES (kit of two) units sold after 06/02

• 27323 SENSOR PAD (box of 20) - replacement

Each sensor can be reused 5-6 times.

• 32205 9V ALKALINE BATTERY

• 27324 CALIBRATION KIT

Frequent calibration is necessary for accurate results.



27322
CE
89/336



(PC not included)

27349

TECHNICAL SPECIFICATIONS

Size: 66 x 176 x 208 mm
 Test frequency: approx. 50 kHz
 Test current: less than 1mA
 Weight: 1,500 g (unit, battery and probe).
 Power source: two 9 Volt DC batteries or 230V (by transformer)
 Norms: 89/336/EEC/ (it is not a medical device)

BODY FAT MONITOR



• 27348 BODY FAT MONITOR

Portable and convenient body fat monitor for private use.

Main features:

- body fat analyzer function
- body fat ratio function
- BMI function
- calorie calculation function
- auto power off (when no operation is detected for approximate 2 minutes)
- memory function stores up to 9 personals data
- 1 personal data for guest (will be initial after power off)

TECHNICAL SPECIFICATIONS

Size: 240 x 160 x 43 mm
 Weight: 305 g
 Power supply: 2xAAA batteries - not included

PROFESSIONAL SKINFOLD CALIPERS: MECHANICAL AND DIGITAL



- 27320 PLICOMETRO - Professional skinfold caliper - mechanical
- 27346 PLICOMETRO - Professional skinfold caliper - digital with PC cable for data transfer

This high precision instrument measures the thickness of a fold of skin with its underlying layer of fat giving physicians the situation of the subcutaneous tissues of the patient. Possibility of re-calibration. Produced to meet International Standards, the Gima "Skinfold caliper" has springs which exert a pressure of 10 g/mm² on the skinfold and an accurate scale which measures thickness in mm. With possibility of re-calibration of the pressure it solves the major problem of traditional skinfold caliper that is reduction of the tension of the spring in the time.

Easy to use and read

With double reading scale (code 27320) or digital display (code 27346), it is supplied with picture showing advised positions for measurements. Tables will help to obtain correct % of bodyfat from measured data and make a correct evaluation of the patient.

Stylish and user-friendly

A special nylon finish makes the device feel skin, avoiding trouble of contact with metal.

Made in metal and sold in an attractive wooden box.

Made in Italy, provided of medical device CE certification.


User manual in English, French, Italian, Spanish.

Data management on PC (only 27346)

Compatible with Windows 2000 and XP, does not work with Windows 98.



SCALE WITH BF/BW

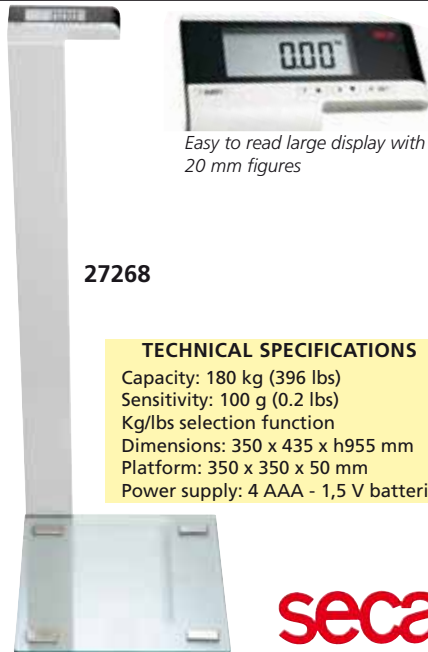

TECHNICAL SPECIFICATIONS

- Capacity: 180 kg (396 lbs)
- Sensitivity: 100 g (0.2 lbs)
- Kg/lbs selection function
- Dimensions: 350 x 435 x h955 mm
- Platform: 350 x 350 x 50 mm
- Power supply: 4 AAA - 1,5 V batteries

• 27268 SECA SUPRA PLUS 720 DIGITAL SCALE WITH BF/BW

Column scale with body fat (BF) and body water (BW) analysis and automatic BMI calculation. The hip-high keyboard makes bending over a thing of the past. The ITO-coated glass platform, which is joined to the base by four stainless steel elements, ensures accurate measurements and great stability.

- Analysis of BMI, body fat and water. For better control
- Big display on a high column. For easy reading
- Keyboard below display. For easy data entry.
- Functions: auto-HOLD, step-off, BMI, BF/BW-measurement, user memory, automatic user recognition, tip-on automatic switch-on, automatic switch-off.

FAT-1 SKINFOLD CALIPER


• 27344 FAT-1 SKINFOLD CALIPER

Economic but precise and accurate skinfold caliper for professional use, made in Italy. Measures the thickness of skin folds over the body, based on the characteristics developed by Tanner Tanner and Whitehouse:

- pressure on the fold: approx 10 g/mm²
- contact surface on the fold: approx 6x11 mm

This device allows a correct evaluation of the nutritional status and the sectoral distribution of adipose tissue. Raises a fold of skin with your fingers and measures its thickness by the forceps: on the top edge a scale in millimeters and an arrow measures the skin fold (triceps, subscapular muscles, etc...) with an accuracy of about 1/2 mm. The coloured bands placed near the scale allow a simultaneous comparison of the measured value with referential percentages related to men and women. This information is included in the manual provided with the unit, which helps an easy interpretation of the values. Advised in dietetics, cardiology, sports medicine, aesthetic medicine, gynecology, pediatrics, ecc...

Multilanguage manual: English, Italian, French, German, Spanish and Greek.

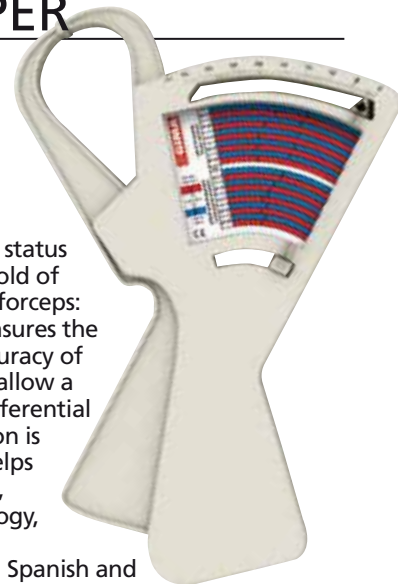

BMI TAPE



• 27341 BMI - Tape Measure

Pocket retractable tape measure 150 cm. BMI (Body Mass Index) calculator on front panel.

• 27342 BMI - Tape Measure

Similar to 27341 but made in China. Available in 40 days only for orders of 1,000 pcs.



BODY TAPE

• 27343 BODY TAPE MEASURE

Track your body measurements and progress.

Measures any body part: arm, thigh, calf, chest, waist, hips....

Measures in inches and centimeters; total measurement length: 60 inches / 152 cm.



Stylish design, easy measurement. Sturdy vinyl tape. Push-button retraction and locking feature ensure easy measurement and superior accuracy.



HEIGHT MEASURE



LCD digit height: 9 mm

Size: 345x92x28 mm

Power supply: 2x1.5V AA batteries

SOEHNLE
PROFESSIONAL
27333

• **27335 TAPE HEIGHT MEASURE**

Clear scale 0-200 cm. Easy wall fixation.

• **27333 ELECTRONIC HEIGHT ROD - Soehnle**

Measuring range: 50-300 cm (20-118 inch), division: 1 cm.

PODOSCOPE



27353

Size: 46 x 55 x h 33 cm

• **27353 PODOSCOPE**

The podoscope is a device for the plantar analysis of the feet. It has a lacquered wooden structure and a crystal top with an inclined mirror below it for a real view of the sole of the feet, which highlights the points of greatest and least load. Double lateral illumination with high level effect lamps. Max Load 135 kg. Strong reliable and professional product, made in Italy.

POSTURAL ANALYZER



27352

Size: 80x72x h 225 cm

• **27352 POSTURAL ANALYZER**

The postural analyser consists of a bilaminated a footboard with adjustable feet. The two lateral aluminium rods are equipped with millimetric scales along with measuring sliders and strings to allow correct body alignment (Barré's vertical). Adjustable upper mirror.

SCOLIOMETER

• **27351 SCOLIOMETER**

Plastic instrument for the evaluation of the scoliosis (deviation of the spinal column on the lateral, left and right planes). The instrument detects the anomaly and its gravity. It also allows to compare reports during the time and define the evolution of the measured pathological process, thus reducing the need of frequent radio graphical tests. The small sphere oscillating within an air space indicates the spine inclination degree, compared to the horizontal plane.



27351


PEDIATRIC MEASURES

seca

Although rugged, the surface is easy to care for and can be cleaned with all commercially available disinfectants

CE 0068

27326

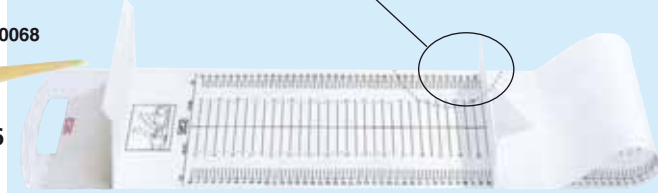

27329

• **27329 SECA 210 - BABY MEASURING MAT**

Mobile measuring mat for babies and toddlers. Kind to the skin, washable, foldable. The space saving solution for easy and precise measurement of baby length. Range 10-99 cm graduation: 5 mm.

• **27326 BABY MEASURE "CALIBRO" Weight 300 g.**

Made of oxidized aluminium, metric graduation 0-900 mm, division 5 mm. Certified by Italian "weights and measure institute". Weight: 300g.

CE 0068

27327



• **27327 CRANIOMETRO - HEAD MEASURE**

Made of special, soft non stretchable plastic. Measures diameter of child head up to 60 cm.

• **27337 ORCHIDOMETER - wood**

Used by pediatricians for screening of children's testicles



27337

A string of 12 numbered wooden beads of increasing size from about 1 to 25 mm used to determine size of testicles.

GONIOMETERS

Used by orthopaedic surgeons to measure the range of movement of joints: elbows, knees and fingers.

27340

• **27340 GONIOMETER - 205x45 mm**

With circular scale used to check the size of hip joints and the angle of the femoral neck on x-rays. Personalisation with your logo for minimum 2,000 pcs.



27339

• **27339 GONIOMETER WITH PAIN SCALE RULER - 203x45x7 mm**



27338

• **27338 GONIOMETER WITH ARMS (tolerance 1°)**

To measure analytical motility of various body segments in degrees (from 0 to 360 clockwise or anti-clockwise direction). Consists of two perspex arms (3mm thick) that will remain open at any given position due to the friction at the joint between them. Two graded scales: one increasing scale and one decreasing scale.

BONE MEASURE


• **27325 BONE MEASURE**

Slide-knot caliper in anticorrosive to measure knee, shoulder, wrist, elbow...

Measuring

capacity: 150 mm. Burner length: 100 mm.

Tested by Italian "weights and measures institute".

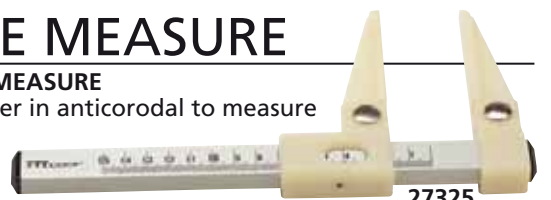

27325

CE 0086

SECA ELECTRONIC BABY SCALE



The weighing tray can be removed quickly and easily at the touch of a button

seca

27283

27282

- 27282 SECA 354 - electronic
- 27283 SECA 384 - electronic - class III



Two scales in one: with a comfortable tray for weighting in a reclining position, easily detachable to be converted into a floor scale. A further advantage is its breast-milk-intake function which determines how much the baby has consumed. Other functions: tare, auto-hold, kg/lbs switch over automatic switch off.

TECHNICAL SPECIFICATIONS

Capacity: 20 kg	Sensitivity: 10 g < 10 kg > 20g
Dimensions: 552 x 156 x 332 mm	Weight: 2.3 kg (27282), 2.8 kg (27283)
Batteries: operated	Made in Asia

SECA MECHANICAL BABY SCALE



CE

seca

27299

- 27299 SECA 725 - mechanical
- 27294 SECA 745 - mechanical - class III



The sturdy classic: with its softly curved weighing tray, the stable all-steel casing and the precise, smooth running weights, this model is the indestructible original baby scale.

TECHNICAL SPECIFICATIONS

Capacity: 16 kg	Sensitivity: 5 g (27299) 10 g (27294)
Dimensions: 550 x 190 x 305 mm	Weight: 6 kg
Made in Europe	

HOME CARE BABY SCALE



CE

90/384/EEC
89/336/EEC

SOEHNLE
PROFESSIONAL

27309

- 27309 SOEHNLE 8310 DIGITAL BABY SCALE

Well-designed, easy to handle baby scale made out of plastic

- Detachable baby tray scale for children up to 20 kg.
- Quattro-Tronic technology with 4 load cells for fast readings.
- Fully-electronic scale: no wear and tear parts, no maintenance.
- Easy to read LCD-display with 20 mm figures.
- One button control for on-tare-off functions.
- "Weight-lock" function for precise and reliable readings.
- Tare function to calibrate for baby blankets/textiles.

- 27304 TRANSPORT CASE
- 27311 TRAY FOR 8310 - spare

TECHNICAL SPECIFICATIONS

Capacity: 20 kg	Dimensions: 530x240x100 mm
Battery operated 9V - not included	Sensitivity: 10 g
Warranty: 2 years	Weight: 4 kg
	Made in Germany

HOSPITAL DIGITAL BABY SCALE



Mains or battery operation

CE III M

27263+27264


- 27263 SOEHNLE 7725 HOSPITAL SCALE - digital

A two-range scale with high precision and reliability suitable for hospital use (approved class III, MGP class I). The damping and hold function guarantee a stable weight indication. With high stability for safe weighing of our babies. Easy to clean plastic and ergonomic forms without chamfers and components allow a simple and hygienic cleaning.

- High-contrast LCD display with 20 mm high figures
- Large reading angle
- Automatic zeroing at switch-on
- Tare function key (laid on pad can be tared)
- Hold function

- 27264 BABY HEIGHT ROD

Digital and mechanical rod, 35 to 80 cm, division 1 mm.

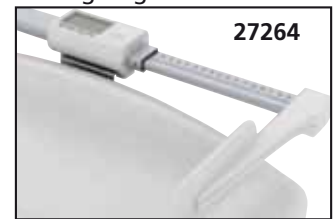

27264

GIMA ELECTRONIC BABY SCALE

- 27313 GIMA ELECTRONIC BABY SCALE

Handy, reliable and well-designed electronic baby scale. Provided with large LCD display 65x26 mm and kg/lb switch. Supplied with height tape meter and washable sheet.



GIMA

27313

Multilingual user manual and gift box

TECHNICAL SPECIFICATIONS

Capacity: 20 kg / 44 lbs	Sensitivity: 10 g
Dimensions: 595 x 340 x 90 mm	Weight: 4 kg
Battery operated: 4x1.5V AA	Auto-off function



FAMILY BABY SCALE

- 27301 FAMILY BABY SCALE

Baby scale suitable for domiciliary use.

Capacity 20 kg/44 lbs
Scale division: 50 g (0-10 kg)
100 g (10-20 kg)

Plastic tray
Dialplate in kg and lbs
Manual zero adjustment



CE 0044

27301

SOEHNLE PROFESSIONAL DIGITAL SCALES

- 27275 SOEHNLE 7830
- 27276 SOEHNLE 7831 with heightmeter
- 27277 SOEHNLE 7730

The new range of Soehnle professional personal scales with stand is a special highlight with its gently curved glass indicator. While stepping on the scale, the display is activated automatically and the weight is shown with redbacklighted numbers.

It does not require any operation and is perfectly suitable for simple weighing in medical institutions, wellness areas and pharmacies.

Scale with integrated height rod Soehnle 7831 (27276) with stand and height rod switches automatically from kg to cm if the height rod is raised.

Thus, patients can be simultaneously weighed and measured.

Features:


- Indication of weight, height and BMI (27276 only)
- Integrated height rod from 120 cm to 202 cm (27276 only)
- No keys, simple operations
- Glass scale indicator (glass thickness 8 mm)
- LED display: 5 digits, digit height: 15 mm
- RS232 interface for PC/EDP included

Made in Germany.

- 27278 SOEHNLE 7708 multifunction scale

With 3 handrails for handicapped and adiposity patients and drive-on-ramp for wheelchairs. Same features as above scales.

SOEHNLE
PROFESSIONAL

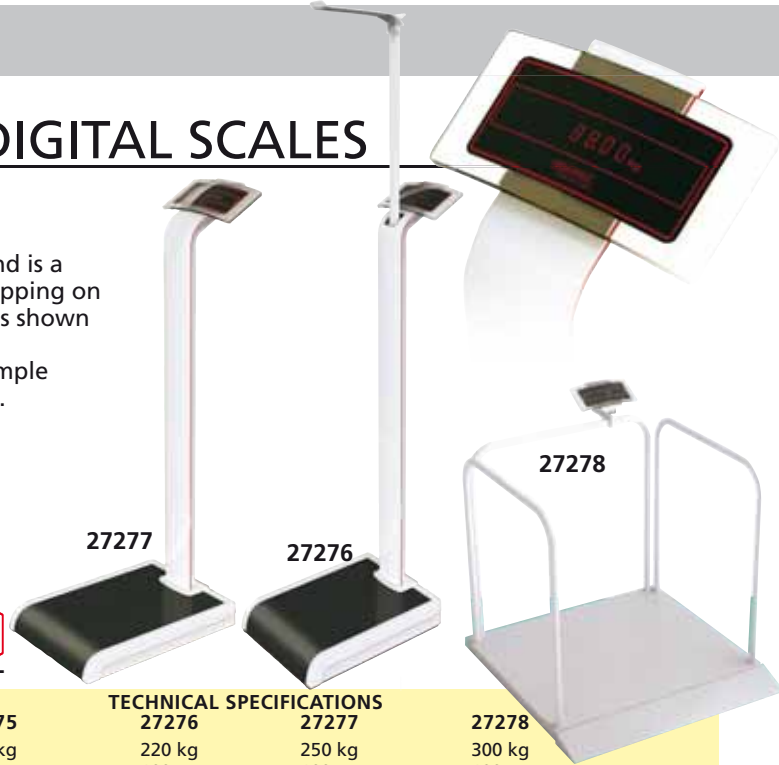

	27275	27276	27277	27278
Capacity:	220 kg	220 kg	250 kg	300 kg
Sensitivity:	100 g	100 g	100 g	100 g
Dimension mm:	369x609x893	369x609x1279	369x609x893	840x880x55 mm (platform)
Weight:	20.8 kg	20.8 kg		
Power Supply:	230 V - 50-60 Hz (available on request 110 V)			

DIGITAL SCALES

7 digit LCD backlit display with 20 mm high figures

Easy to read with a 180° rotary display

Power supply or rechargeable batteries

Dual kg/lb weight unit

Large platform 890 x 890 mm

Platform 310x275 mm

UP TO 330 KG 27279

All models with 5 languages manual (English, French, German, Spanish, Italian)

SOEHNLE
PROFESSIONAL

SOEHNLE DIGITAL SCALES - Made in Germany
All scales are provided with dual kg/lb readings.

- 27290 SOEHNLE 7800 - 200 kg
Digital scale with 850 mm stand, easy to read display with 20 mm figures. Easy one button handling.
- 27269 HEIGHT METER - optional for scale 27290, 27291, 27295
Maximum height 210 cm
- 27279 SOEHNLE 7703 scale with holding device, for handicapped and adiposity patients.
Function key either for BMI (Body Mass Index), hold, and dialysis function.

TECHNICAL SPECIFICATIONS

	27290	27279
Capacity:	200 kg (440 lb)	330 kg (792 lb)
Sensitivity:	100 g	100g<200 kg>200 g
Dimension mm:	350 x 450 x 860	890 x 800 x 1,484
Weight:	10 kg	35 kg
Power Supply:	230V - 50/60 Hz	
Battery life:	200 hours	40 hours
	with four 1.5 V mignon batteries	

Keys to symbols:

= Metrological Approval for use in medical field
 = Calibration according to class III sensitivity required by 90/384/EEC

DIGITAL CHAIR SCALE

SOEHNLE
PROFESSIONAL



- 27273 CHAIR SCALE SOEHNLE 7802

Designed for weighing sick or handicapped patients. Fully electronic scale: no wear and tear parts, no maintenance. Easy to read display with 20 mm high figures. One button control for on-tare-off functions. Tiltable arms for easy access for bed-ridden and frail patients. Powered by either mains or by batteries (4x1.5V mignon batteries) Battery operating life up to 150 hours.

TECHNICAL SPECIFICATIONS

Capacity:	200 kg (440 lb)
Sensitivity:	100 g
Total weight:	27 kg
Dimensions:	580 x 1,100 x 930 mm
Operating voltage:	230 V / 50-60 Hz, 10 mA, 2.2 VA
Made in Germany	

VETERINARY SCALE
SEE PAGE 244

BED DIGITAL SCALE

• 27270 SOEHNLE 7711 BED SCALE BLUETOOTH
Modern bed scale with bluetooth™ link, is perfectly suited for use in hospital. Its four radio-controlled interlocking weighing elements make it easy to install as no complicated cables are necessary. Features:
Approved to Class III, MPG Class 1
Wireless transfer of measurement data
Easy-Lift-System - innovative lift technology
Low power requirement of foot sensors
BMI, hold and dialysis function.

TECHNICAL SPECIFICATIONS

Capacity:	500 kg / 1100 lb
Sensitivity:	0-300 kg: 100 g/0.2 lb 300-500 kg: 200 g/0.4 lb
Power supply:	230V/50Hz
Made in Germany	

SOEHNLE
PROFESSIONAL



WHEELCHAIR SCALE



Capacity: 200 kg
Sensitivity 100 g
Size: 800x800x40 mm

SOEHNLE
PROFESSIONAL

• 27274 SOEHNLE 7805 WHEELCHAIR SCALE
Transportable by means of 2 rollers and handle. Extremely lightweight platform weighing only 14 kg, with aluminium frame. Weighing platform made of waterproof bonded laminated wood. Anti-slip surface.

= Approved scales according to 93/42 EEC

DIGITAL SCALE



- **27288 PEGASO DIGITAL SCALE with height meter (70x190 cm)**
Professional digital scale supplied with integrated height meter. This scale can be operated by mains or batteries.
- **27281 POWER SUPPLIER**

TECHNICAL SPECIFICATIONS

Capacity: 200 kg Sensitivity: 100 g
 Dimensions: 535 x 275 x 940 mm
 Platform: 372 x 275 mm Weight: 15 kg
 Power supply: 110-240 V - 50/60 Hz
 Battery accumulator: 6V 5AH
 Recharging time: 4-5 hours
 Operating time when fully charged: 36 hours

- **27310 ASTRA SCALE - 200 kg mechanical with height meter (75-200 cm)**
Traditional scale with robust lever system, powder coated platform and column, provided with wheels for easy transport. User manual in English, Spanish, French, Italian, German, Greek, Portuguese, Arabic.

TECHNICAL SPECIFICATIONS

Capacity: 200 kg
 Sensitivity: 100 g
 Weight: 15 kg
 Dimensions: 275 x 1,485 x 530 mm
 Platform: 375 x 275 mm

ASTRA SCALE

Reading on both sides of measuring rod



CE 0044

Platform with antiskid cover

27310

SECA MEDICAL SCALES - CLASS III

CE 0109
 III M



Dial with BMI



27284

- **27284 SECA 756 - mechanical with height meter - (60-200 cm)**
Tradition combined with modern design. With a very sensitive precision mechanism to show weight through the aid of a well visible red pointer on the circular dial. User manual in English, Spanish, French, German, Italian.

TECHNICAL SPECIFICATIONS

Capacity: 160 kg
 Sensitivity: 500 g
 Weight: 18.5 kg
 Dimensions: 628x360x1,000 mm
 Platform: 360x350 mm with wheels

CE 0109
 III M



Supplied with height meter

27296

- **27296 SECA 711 - mechanical with height meter - class III - (60-200 cm)**
- **27297 SECA 700 - mechanical with height meter - (60-200 cm)**
Mechanical column scale with weight adjustment at eye level. Weight can be measured in a conventional way: by aligning two chrome plated indicators. User manual in English, Spanish, French, German, Italian.
- **27298 SECA220 HEIGHT METER - spare**

TECHNICAL SPECIFICATIONS

Capacity: 220 kg
 Sensitivity: 100 g (27296) - 50 g (27297)
 Weight: 21.7 kg
 Dimensions: 520x1,556x520 mm
 Platform: 350x360 mm

CE 0109
 III M
 Made in Europe



27292

The non slip surface on the platform guarantees a safe and secure foothold



All it takes is the press of a key to find out BMI



27292 + 27298

- **27292 SECA 799 - electronic with BMI - class III**
- **27298 SECA 220 HEIGHT METER**
Versatile, low maintenance and site independent with economic battery operation, ideal for daily use in hospital and doctors practices. Column scale is easy to move around with its transport castors. BMI function permits a reliable evaluation of the nutritional conditions of the patient. Equipped with several functions: TARE, HOLD, weighing range switch-over automatic switch off. User manual in English, Spanish, French, German, Italian.

TECHNICAL SPECIFICATIONS

Capacity: 200 kg
 Sensitivity: 100 g < 150 kg > 200 g
 Weight: 6.3 kg
 Dimensions: 294x831x417 mm
 Platform: 272x280 mm

Keys to symbols:

III = Metrological Approval for use in medical field
 CE III = Calibration according to class III sensitivity required by 90/384/EEC

CE 0109 = Approved scales according to 93/42 EEC

GIMA PROFESSIONAL SCALES



27239

TECHNICAL SPECIFICATIONS

Capacity: 160 kg Dial Ø: 180 mm
Sensitivity: 1,000 g
Weight: 3 kg
Dimension: 430 x 380 x 70 mm



- **27239 GIMA PROFESSIONAL SCALE**
Mechanical scale with classic design: a large and clear dial, a robust housing made of steel. Supplied in multilingual gift box. Made in Europe.



- **27240 DIGITAL PROFESSIONAL SCALE**
Digital scale with extra large, easy to read LCD display (65 mm). Strong and durable steel construction comfortable platform. Lb/kg/st switchable by user. Supplied in multilingual gift box. Made in Europe.



27240

TECHNICAL SPECIFICATIONS

Capacity: 200 kg
Sensitivity: 100 g
Power: lithium battery (included)
Dimension: 430 x 380 x 70 mm



GIMA BATHROOM SCALES



27230

- **27230 MECHANICAL SCALE**
Analogue bathroom scale with stable construction and plastic cover.

TECHNICAL SPECIFICATIONS

Capacity: 125 kg (275 lb)
Sensitivity: 1,000 g (2 lb)
Warranty: 2 years



LCD display
76x36 mm

27232

- **27232 DIGITAL SCALE**
Electronic bathroom scale with plastic platform. 4 sensor high precision system with automatic switch off. Lb/kg/st switchable by user.

TECHNICAL SPECIFICATIONS

Capacity: max 150 kg (330 lbs) - min. 3 kg (6.6 lb)
Sensitivity: 100 g = (0.2 lb)
Dimension: 302 x 315 x 40 mm
Power: battery: 2xAAA - not included
Warranty: 2 years



LCD display
60x26 mm

27233

- **27233 DIGITAL GLASS SCALE - white**
- **27234 DIGITAL GLASS SCALE - orange**
Electronic bathroom scale with 6 mm thick tempered glass platform. 4 sensor high precision system with automatic switch off. Lb/kg/st switchable by user.

TECHNICAL SPECIFICATIONS

Capacity: max 150 kg (330 lb) - min. 3 kg (6.6 lb)
Sensitivity: 100 g (0.2 lb)
Dimension: 285 x 285 x 32 mm
Power: 1x CR2032 battery - not included
Warranty: 2 years



DIGITAL BATHROOM FASHION SCALES



ULTRA SLIM


27235

Foldable
LCD display
80x44 mm

- **27235 SLIM DIGITAL SCALE**
Ultra flat digital scale in an elegant silver colour glass design. Only 18 mm high for excellent stability. Automatic on/off function weighs immediately upon stepping on the scale. Shuts down when not in use.

TECHNICAL SPECIFICATIONS

Capacity: 150 kg (330 lbs) - min. 3 kg (6.6 lb)
Sensitivity: 100 g (0.2 lb)
Platform: tempered glass 300 x 300 x 6 mm
Power: 2x CR2032 battery - not included
Warranty: 2 years

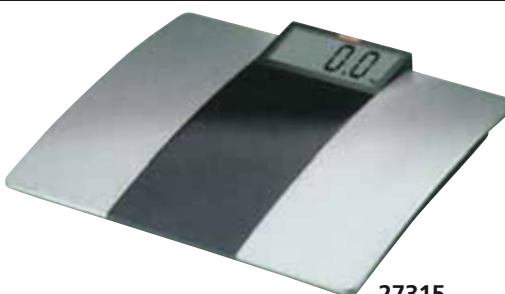

27315



- **27315 VERONA DIGITAL SCALE**
Digital bathroom scale with easy-clean plastic platform cover and bright LED digits. Automatic Start.

TECHNICAL SPECIFICATIONS

Capacity: 150 kg (330 lbs)
Sensitivity: 100 g (0.2 lb)
Power: 2x3V lithium battery - not included
Warranty: 3 years



LCD display
78x38 mm

27236

- **27236 GLASS SCALE**
Electronic bathroom scale with 6 mm thick safety glass platform. 4 sensor high precision system. Automatic switch off. Lb/kg/st selectable by user.

TECHNICAL SPECIFICATIONS

Capacity: max 150 kg (330 lbs) - min. 3 kg
Sensitivity: 100 g (0.2 lb)
Platform: tempered glass 300 x 300 x 6 mm
Power: 2x CR2032 batteries - not included
Warranty: 2 years



PROFESSIONAL SCALES

ANALOGIC AND DIGITAL 180 KG

• **27242 MEDISCALE SOEHNLE 7840**
Mechanical-digital scale with high precision and elegant design. The weight is displayed on a LCD display and on the glass scale.



Large display with 25 mm high digits



27242

TECHNICAL SPECIFICATIONS

Capacity: 180 kg LCD display and glass scale
Sensitivity: 100 g Dimension: 550 x 325 x 120 mm
Weight: 4 kg Battery: 4x1.5V AA - not included

SOEHNLE PROFESSIONAL



27242

• **27285 SECA 762 SCALE**
Traditional model with large easy-to-read dial. Pointer professional scale in strong metal structure white painted. Consistent accuracy on any surface. Spacious platform design with no slip rubber mat, stable for safety. Tare function.



27285

TECHNICAL SPECIFICATIONS

Capacity: 150 kg
Sensitivity: 500 g ± 2 g
Weight: 4 kg

seca

BODY COMPOSITION ANALYZERS - SCALES



16.5 mm

ULTRA SLIM

27238

TECHNICAL SPECIFICATIONS

Capacity: max 150 kg (330 lb) min 3 kg (6.6 lb)
Sensitivity: 100 g (0.2 lb)
LCD Display: 80 x 44 mm
Tempered glass: 300 x 300 x 6 mm
Battery: 2x CR2032 - not included

• **27238 WELLNESS SCALE**
Measuring body weight, fat, water, muscle content (in percentage) and daily caloric consumption. Overload and low battery indication, weight in kg/lb/st 10 memories.

• **27237 BODY FAT SCALE**
Measuring body weight, fat, water, muscle content, bone mass weight and daily caloric consumption. Overload and low voltage indication. 4 memories with individual keys for easy operation. Display all data at the same time.



27237

TECHNICAL SPECIFICATIONS

Capacity: max 150 kg (330 lb) min 3 kg (6.6 lb)
Sensitivity: 100 g (0.2 lb)
LCD Display: 93 x 45 mm
Battery: 4xAAA - not included

PRECISION SCALE



27246

• **27246 PRECISION SCALE**
With tempered glass Ø 160 x 4 mm tare function, kg/lb-oz

TECHNICAL SPECIFICATIONS

Capacity: 5 kg LCD display: 35x15 mm
Sensitivity: 1 g Size: 190 x 160 x h 50 mm
Battery: 1xCR2032 - not included

SLIM SCALE



ULTRA SLIM

• **27247 SLIM PRECISION SCALE**
With safety tempered ultra slim glass 248x180x3 mm. Touching key, overload indication, tare function, back light.

TECHNICAL SPECIFICATIONS

Capacity: max 5 kg (11 lb) min 2 g (0,207 lb)
Sensitivity: 1 g (0,05 oz)
LCD display: 65 x 35 mm
Size: 248 x 180 x 16 mm
Battery: 3xAAA 1.5V - not included

DIET SCALE

990 DIFFERENT FOODS


• **27249 DIET SCALE**
Measuring weight and nutrition value of foods. Just input food code, you can easily know calories, fat, cholesterol, fiber, protein, sodium, carbohydrates, of the food. Analysis of 990 kind of foods, ingredients and liquids with a high precision strain gauge sensor system. Ultraslim design: 17 mm and large LCD display.



TECHNICAL SPECIFICATIONS

Capacity: max 5 kg (11 lb)
Sensitivity: 1 g (0,05 oz)
Lithium battery: not included

THERMOFOCUS® PROFESSIONAL "NO CONTACT" THERMOMETERS

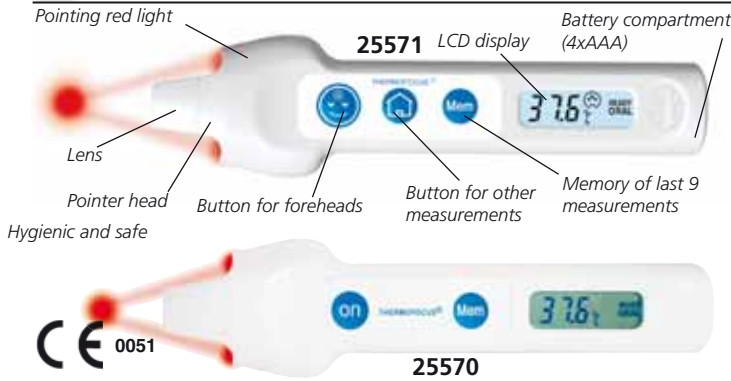

THERMOFOCUS THERMOMETERS

The most advanced thermometers in the world! ThermoFocus read infrared radiations emitted by the surface of the skin and calculates accurately the body temperature.

- **without touching:** never in contact to avoid any risk of infection
- **non invasive:** do not require patient collaboration, advised for children
- **accurate:** a sophisticated microprocessor warrants high precision
- **economic:** no need of expensive probes
- **fast:** takes the temperature in one second
- **1°C-55°C:** for measurement of wide range of objects
- **body temperature scan (25571/2):** can be done quickly
- **25570 THERMOFOCUS - mod. 0700 - basic - blister**
2 buttons - for the private user who needs an easy, fast and precise thermometer. Simply press the ON button, hold it perpendicular to the center of the forehead and release the button as soon as the two red lights converge in a single point. It can also be used to take the temperature of objects such as a baby's feeding bottle. Not recommended for doctors and hospitals.
- **25571 THERMOFOCUS - mod. 1500 - 6 in 1**
3 buttons - for the family or the private doctor who require maximum accuracy and flexibility of the temperature taking. Same as 25570 with an additional button ("home") to allow the user to read the ambient temperature and the one of other areas of the skin, as well as the temperature of objects and liquids within a range from 1 to 55°C (33.8 to 131°F). Provided with the exclusive Manual Quick Calibration System (MQCS) technology to quickly stabilize the device if moved between rooms with the different temperatures.
- **25572 THERMOFOCUS - mod. 1500 H1N1 - emergency version**
Specially made for emergency situations can take thousands of measurements every hour. Same as 25571 but with compulsory MQCS function. It is supplied with a special lanyard to be worn around the neck.
- **32200 "AAA" BATTERIES**

TECHNICAL SPECIFICATIONS

Forehead measuring range:	34.0/42.5°C (93.2/108.5°F)
Measuring range (apart from forehead):	1.0/55.0°C (33.8/131°F)
Ambient temperature working range	16/40° (60.8/104.0°F)*
Precision:	0.1°C (0.2°F)
Accuracy level according to ASTM E 1965-98 :	
from 34 to 35.9°C (93.2/96.6°F)=	+/-0.3°C (+/-0.5°F)
from 36 to 39°C (96.8/102.2°F)=	+/-0.2°C (+/-0.4°F)
from 39.1 to 42.5°C (102.4/108.5°F)	+/-0.3°C (+/-0.5°F)

Distance of operation: 3 cm (1.2 inch) fixed through optical signal
 * The thermometer can work also in environments with a temperature between 5 and 16°C (41/60.8°F), but we do not guarantee the accuracy.
 Battery: 4xAAA (1.5V) included Weight: 99 g Size: 165x40x22 mm
 Made in Italy

THERMOFOCUS 5 IN 1



EAR THERMOMETER

- **25580 IR EAR THERMOMETER**
- **25581 PROBE COVERS - box of 40 pcs.**
- **32212 SPARE BATTERY**
- Fully meet "ASTM" standard and CE EN 12470-5 (only few brands passed it) in measurement accuracy (35,5 °C to 42 °C ± 0.2 °C, others ± 0.3 °C).
- Low customer-return-rate, under 0.3%
- Unique and handy design
- One-second measurement, fast and easy
- 9 memory for temperature comparison and evaluation, fever alarm, easy reading LCD display
- Self-diagnostic and infrared measurement

TECHNICAL SPECIFICATIONS 25580-25585

Measurement Range:	- forehead 34.0°C~42.2°C (93.2°F~108.0°)
	- scan mode: -22/80°C (-7.6/-76°F) - 25585 only
Response time:	1 second (25580) 5-30 sec. (25585)
Time between Measurement:	< 5 sec
Display Model	ORAL °C/°F
Battery Type:	CR2032 (3212)
Dimension (LxWxH):	14 x 3.8 x 3 cm - 70g
Norms:	IEC60601-1, IEC60601-1-2 (EMC), ASTM E1965-98




25580

FOREHEAD THERMOMETER

- **25585 INFRARED MULTI-FUNCTION FOREHEAD THERMOMETER**
- **32212 SPARE BATTERY**
- Multifunctional innovative infrared thermometer with 5 applications: forehead temperature, room temperature, scan mode for surface temperatures (bath water or milk), real time clock, stopwatch to monitor the heart beat
- New 3 LED zone indicator to get fever status friendly (see picture)
- Instant measurement
- Large LCD to read easily
- Detailed 9 memories with time/date & site
- Fever alarm with beeper
- Selectable °C/°F mode
- Small size to carry conveniently when traveling
- Accurate: ± 0,2/0,3 °C
- Self-diagnosis warning when device is abnormal.
- Packed in deluxe carton box.



25585

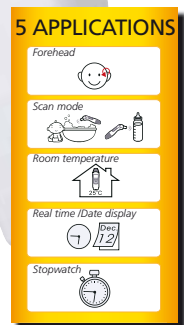

3 LED Fever Indicator



43

25593

25594

25597

LIQUID CRYSTAL THERMOMETERS

Non toxic, mercury free, safe, unbreakable, reliable, long lasting, easy to read, liquid crystal thermometers. 10 languages manual. Made in Italy.

REUSABLE FOREHEAD THERMOMETERS

- **25593 FEVERTEST - display of 25 blister of 2 pcs.** It indicates the normal temperature or fever in 15 seconds.
- **25596 FEVERTEST BULK - box of 100 - hospital use**
- **25597 FEVERCONTROL - display of 25 blister** It indicates temperature in 15 seconds with graduation of 0.5°C.

BODY TEMPERATURE MONITOR

- **25594 48h FEVER MONITOR - display of 25 blister of 3 pcs.** Innovative "continual" fever indicator. For 48h monitoring.



25597

DIGITAL THERMOMETERS

30-60 SECONDS CE 0476



- 25559 GIMA DIGITAL THERMOMETER - °C - box - 30 seconds
Sold in a counter display of 36 pcs. (English - Italian)
- 25560 GIMA DIGITAL THERMOMETER - °C - hang box - 60 seconds
- 25561 GIMA DIGITAL THERMOMETER - °F - hang box - 60 seconds

30 SECONDS CE 0476



- 25562 W-P DIGITAL THERMOMETER °C - hang box water-proof

30 SECONDS CE 0476



- 25563 FLEXI DIGITAL THERMOMETER °C - hang box flexible tip, water-proof

30 SECONDS CE 0476



- 25565 BACKLIGHT DIGITAL THERMOMETER °C - hang box flexible tip, backlight display



- 25609 GALAXY DIGITAL THERMOMETER °C - flexible tip, water-proof, with large LCD display (21x7.4 mm) and hang box

BATH THERMOMETER



25567

- 25567 BATH THERMOMETER °C - Star Measure water temperature in 90 sec. Measurement range 0° C to 70° C. Supplied in hang blister.
- 25598 BATHERMO °C Display 25 blister. Bath thermometer to check water temperature in few seconds. 10 language manual.



BABY THERMOMETERS

30 SECONDS 25566



- 25566 FROG DIGITAL THERMOMETER °C - hang box - flexible tip
- 25568 NIPPLE-LIKE DIGITAL PACIFIER THERMOMETER °C Individually packed in carton box, memorizes last temperature. Auto power-off. Display range 32.0°C-42.0°C, with accuracy ±0.1°C. 1.5V button battery



25560



25565



25609

GIMA DIGITAL THERMOMETERS

- Gima digital thermometers, available in °C or °F, deliver accurate and quick read-out. Safe and easy to use with clear and large display, have following features:
- display range 32.0°C-42.0°C (90.0-107.6°F)
 - accuracy ±0.1°C (±0.2°F)
 - memory: last reading
 - beeper function
 - auto shut-off
 - plastic case
 - battery LR/41 - 1.5 V
 - multilanguage manual and box
 - CE and FDA approved
 - 25549 BATTERY LR/41 - box of 20 25559


Counter display of 36 pieces



BASAL THERMOMETER

- 25608 DIGITAL BASAL THERMOMETER °C - hang box - waterproof With 4 digits, 2 decimals, for fertility control

4 DIGITS, 2 DECIMALS

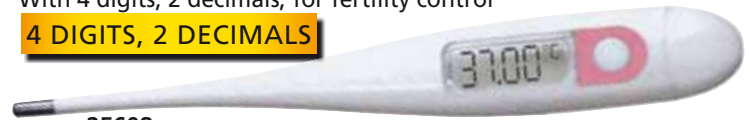

25608

MERCURY FREE THERMOMETERS

MERCURY FREE

CE 0476



25614

25613

25586

- 25613 ECOLOGICAL OVAL THERMOMETER °C - hang box
- 25614 ECOLOGICAL OVAL THERMOMETER - box of 12
- 25586 ECOLOGICAL OVAL THERMOMETER - with shake down aid - hang box. Available january 2012. Mercury free clinical thermometer with gallium alloy instead of mercury. Safe, precise and functional. Glass structure. Temporarily not for sale in Germany.