

# MULTICARE "IN" - SELF CONTROL OF RISK FACTOR


- 23965 MULTICARE "IN" KIT - Italian manual
- 23966 MULTICARE "IN" KIT - English manual
- 23967 MULTICARE "IN" KIT - CD multilanguage manual

\*Includes 10 languages: GB, IT, FR, ES, GR, PT, PL, SL, HU, TR

A new approach for the periodical self-monitor of risk factors - cholesterol, triglycerides and glucose. The only device on sale with reflecto-metrical and electro-chemical bio sensor suction. Multicare "IN" is supplied in carton box with a nylon bag, 10 glucose strips, lancet device, 10 lancets, rapid guide and user manual. Multicare "IN" has 7 language box: Italian, English, French, Spanish, Greek, Portuguese, Polish.


**The first and only multi-parameter device with double technology**  
**BIOSENSOR SUCTION STRIP - glucose**

- Electro-chemical technology for a more practical use.
- Requires very little blood samples (1 µl).
- Short time for the results (5 sec.)



**REFLECTOMETRIC SUCTION STRIP - cholesterol and tryglicerides**

- multiCare in is a unique technology.
- It allows a more functional self-monitoring of triglycerides and cholesterol.

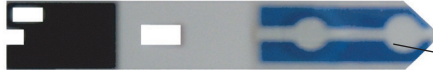

GIMA code	SPARE STRIPS FOR MULTICARE "IN" - 7 LANGUAGE BOX
23969	Glucose strips biosensor - tube of 25 strips + 1 chip
23970	Glucose strips biosensor - tube of 50 strips + 1 chip
23971	Cholesterol strips - tube of 5 strips + 1 chip
23972	Cholesterol strips - tube of 25 strips + 1 chip
23973	Triglycerides strips - tube of 5 strips + 1 chip
23974	Triglycerides strips - tube of 25 strips + 1 chip
CONTROL SOLUTION FOR MULTICARE "IN"	
23978	Glucose control solution
23979	Cholesterol control solution
23980	Triglycerides control solution
LANCETS FOR MULTICARE AND MULTICARE "IN"	
23940	Lancets 28G - sterile - bag of 25 pcs
23916	Lancets 28G - sterile - box of 100

- 23982 SOFTWARE FOR MULTICARE "IN" - English, French, Italian, Spanish
- Operating system necessary to installation: Windows XP, Vista, 16 Mb Ram, 10 Mb remaining space on hard disk



CE 0344

**multiCare in**

**3 IN 1: CHOLESTEROL, TRIGLYCERIDES, GLUCOSE BIOSENSOR**

7 Languages box:  
 English, Italian, French, Spanish, Greek, Polish, Portuguese.



23965

**TECHNICAL SPECIFICATIONS**

	Measuring range	Measuring time	Blood sample
Glucose electro-chemical:	10-600mg/dl	5"	1 µl
Cholesterol	130-400 mg/dl	30"	10 µl
Triglycerides	50-500 mg dl	30"	10 µl
Memory:	500 test results with date and time		
Batteries:	2 lithium batteries Cr 2032 3 Volts (3212)		
Size:	97x49x20.5 mm		
Weight:	65 g		

# CHOLESTEROL, TRIGLYCERIDES, GLUCOSE POCKET MONITORS

- 23948 MULTICARE KIT - English/Spanish, manual
  - 23950 MULTICARE KIT + 10 glucose strips - English/Italian manual
- Unique diagnostic device for self-check of glucose, cholesterol and triglycerides. Multi-parameter technology allows to measure, with the same reflectometer, different parameters by simply changing type of strips: a chip in each pack of strips is specific for each analysis. Regular measuring of these three parameters is important to **prevent cardiovascular pathologies** if we consider that intolerance to glucose, hyper-cholesterolemia and hyper-triglyceridemia are three of the main cardiovascular risk factors, along with age, obesity, lifestyle, smoking and arterial hypertension. Multicare is supplied in 3 language carton box (English, Italian, Spanish) with a nylon bag, lancet device, 10 lancets, user manual (English, Italian, Spanish as above). Available on request French and Portuguese manual.


• 32212 LITHIUM BATTERY CR2032

GIMA code	STRIPS AND CONTROL SOLUTION FOR MULTICARE
23955	Glucose strips - tube of 25 strips + 1 chip
23956	Glucose strips - tube of 50 strips + 1 chip
23957	Cholesterol strips - tube of 5 strips + 1 chip
23958	Cholesterol strips - tube of 25 strips + 1 chip
23959	Triglycerides strips - tube of 5 strips + 1 chip
23960	Triglycerides strips - tube of 25 strips + 1 chip
23961	Glucose control solution
23962	Cholesterol control solution
23963	Triglycerides control solution

**Practical**  
 Device activation by inserting the strip or by pushing the button

**Versatile**  
 Blood sample can be dropped both on the inserted strip and on the external strip.

**Innovative**  
 Reduced blood volume 3 microlitres of blood are needed for the glucose test.

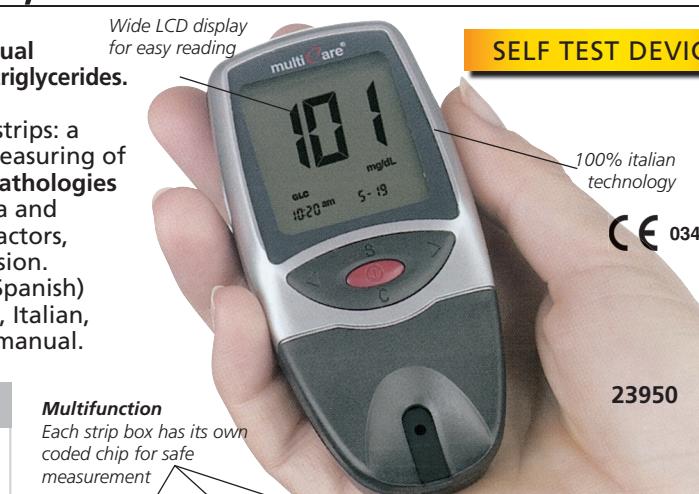

Wide LCD display for easy reading

**SELF TEST DEVICE**

100% italian technology

CE 0344

23950

**Multifunction**  
 Each strip box has its own coded chip for safe measurement



**TECHNICAL SPECIFICATIONS**

Measuring range glucose:	25-500mg/dl (1.3-28 mmol/l)
triglycerides	50-500 mg/dL (0.5-5.6 mmol/l)
cholesterol	130-400 mg/dL (3.3-10.2 mmol/l)
Memory:	255 test results with date and time
Size:	48x104x20.5 mm
Weight:	65 g
Batteries:	2 batteries