


FAMILY AND PINGOO ULTRASONIC NEBULIZERS



GIMA code	ULTRASONIC NEBULIZERS
28106	Family ultrasonic nebulizer
28117	Pingoo ultrasonic nebulizer
28113	Re-chargeable battery Ni Cd 12 VDC - Family
28115	Re-chargeable battery Ni Cd 12 VDC - Pingoo
28114	Spare Kit: includes adult and paediatric facemask, nasal prong, mouthpiece and tube
28121	Medication cups Family (box of 15 pcs.)
28122	Medication cups Pingoo (box of 15 pcs.)

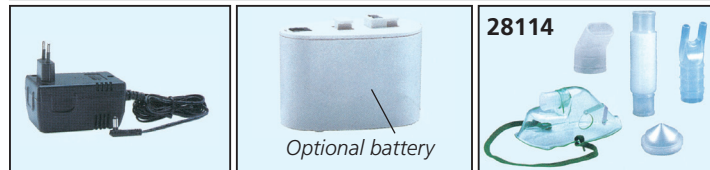

ULTRASONIC NEBULIZERS

Thanks to the innovative electronic technologies these easy-to-use ultrasonic nebulizers provide a quiet, fast and portable operation with optimal particle size.

Portable: small and lightweight, it allows patients to perform treatments anytime and anywhere.

Optimal particle size: particle size determines the extent to which the medication can be inhaled by the patient into the deep lung tissue.

Approx 90% of particles generated by Family or Pingoo have the optimal 1-5 μ (microns) size.

Fast treatment: with its high-efficiency design, it allows a quick treatment time, most within 6 minutes, with limited power consumption.

Low noise level: treatments are not only fast but are virtually silent due to its advanced ultrasonic technology.

8 languages user manual: English, French, German, Italian, Spanish, Portuguese, Greek and Arabic.



TECHNICAL SPECIFICATIONS

Operating voltage: 220-240 V 50/60 Hz Power: 10 W
 Ultrasonic frequency: 2.5 MHz
 Nebulization time: 0.2-0.7 ml/min depending on medication density
 Medication cups capacity: max 8 ml Family, max 4.5 ml Pingoo
 Automatic switch-off
 Particle size: 5 μ (microns) average
 Three programmable flow levels
 Noise level: <40 dB

STANDARD ACCESSORIES

Adult facemask (only Family)
 Pediatric facemask
 Nasal prong, Mouthpiece
 Connection tube

Medication cups
 Transformer (220-240V-12V DC)
 Carrying bag
 User manual



POCKET ULTRASONIC NEB



- 28080 POCKET ULTRASONIC NEBULIZER
 - 28081 MEDICINE CUPS - box of 50
 - 28082 RESPIRATORY ACCESSORY KIT (adult and pediatric mask, mouthpiece, flex tube)
- Small, handy, light weight, portable nebulizer, that can operate with standard AA batteries (conventional or rechargeable) instead of expensive dedicated rechargeable batteries or AC adaptor 100-240 V. Two adjustable nebulization rate. Provides optimal treatment for all kinds of medication with different density. Supplied in a nylon carrying bag with AC adapter adult and child mask, mouthpiece 10 medicine cups, and transformer.

STANDARD AA BATTERY OPERATED

TECHNICAL SPECIFICATIONS

Operating voltage: DC 6 V (4xAA alkaline batteries)
 AC adaptor: input: 100-240 V - 50/60 Hz AC/Output: 12 V / 700 mAh
 Rechargeable battery: 4xAAA NiMH rechargeable batteries / 2000 mAh (90 minutes)
 Ultrasonic Frequency: 2.5 MHz
 Medication cups capacity: $\leq$ 6 ml
 Nebulization rate: 2 speed adjustable
 - 1st level - 0.3 ml/min - 2nd level - 0.5 ml/min
 Particle size: 3-5 μ (microns)
 Treatment time: 5-6 min (2-3 ml)
 Sound level: $\leq$ 50 db



HOMECARE ULTRASONIC NEB

- 28085 HOMECARE ULTRASONIC NEBULIZER
- 28086 MEDICINE CUPS - box of 50
- 28082 RESPIRATORY ACCESSORY KIT (adult and pediatric mask, mouthpiece, flex tube)

Desktop ultrasonic nebulizer for homecare. Supplied with flexible nose, adult and child mask, mouthpiece and 10 medicine cups.

- **3-speed nebulization rate control** - provides optimal treatment for all kinds of medication with different density.
 - **Specialized design** - with swivel air outlet to suit all kinds of treatment.



TECHNICAL SPECIFICATIONS

Operating voltage: 230 V-50 Hz
 Power consumption: $\leq$ 30 VA
 Medication cups capacity: $\leq$ 10 ml
 Nebulization rate: 3 speed adjustable
 Particle size: 3-5 μ (microns)
 Ultrasonic frequency: 2.5 MHz $\pm$ 100Khz
 Size: 220x110x205 mm - 1.3 Kg



28085