


GIMA STRETCHERS AND SPINAL BOARDS


CE



34101

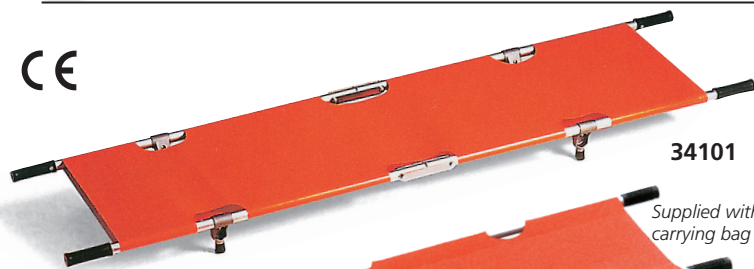

34107

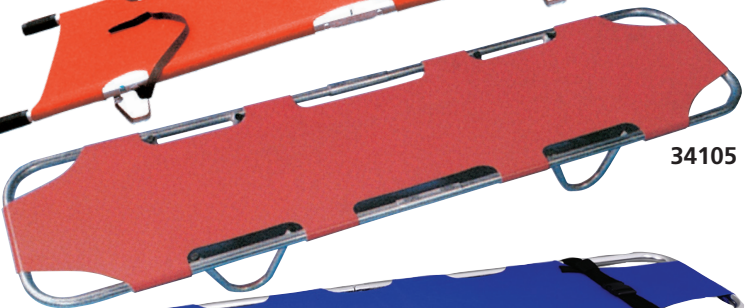


34105



34076



34077

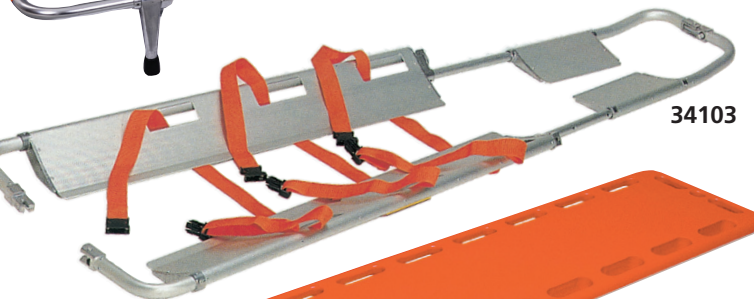

34103



34083

34097+
3401134023+
34024

• **34101 GIMA STRETCHER 2 - foldable in 2**
Emergency stretcher is sturdy, but lightweight in construction, with aluminium frame and three solid iron cross-bars for heavy use. Reinforced nylon canvas, fluid repellent and non-tear.

• **34107 GIMA STRETCHER 4 - foldable in 4**
Made of strong alloy and water-tight material. Light weight, easy to carry, foldable in 4. Length foldable in 2 parts, width foldable in 2 parts.

• **34105 STACKABLE STRETCHER**
Used in small rooms where it is necessary to store a large number of stretchers. PVC does not absorb blood or body fluids. Strong lightweight aluminium construction with vinyl coated polyester

• **34076 BLUE STRETCHER - foldable in 2**
Light-weight, small sized, easy to carry. High-strength aluminium poles and vinyl coated waterproof nylon canvas.

• **34077 WHEEL STRETCHER - foldable in 2**
Wheels, outriggers and adjustable backrest. High strength aluminium poles and vinyl coated, waterproof nylon canvas. Light weight, small sized, easy to carry.

• **34103 SCOOPING STRETCHER**
A wonderfully practical scooping stretcher allows to put the stretcher under the patient without moving him/her much. The shape of the stretcher has been designed to immobilize the patient properly during transport.

- effective hooking and unhooking device for stretcher uncoupling
- adjustable stretcher length with automatic locking device
- bag with 3 quick unhooking safety belts for patient immobilization
- once folded can be easily stored inside the ambulance
- completely built in light and sturdy aluminium


• **34083 SPINAL BOARD**
X-ray translucent, washable, buoyant and impervious to all fluids. Made in linear high density polyethylene with plastic pins. Compatible with all head-immobilization devices and strapping system. CT and MRI transparent. Colour orange.

• **34097 SPINAL BOARD WITH SAFETY PINS**
Strong and lightweight high density PE spinal board, with two fibre reinforced plastic stays. Provided with metal pins for better fixation of immobilisation belts. Excellent in water rescues. X-ray translucent, except 2 plastic stays.

• **34023 BABY GO PEDIATRIC SPINAL BOARD**
Baby go is a professional board, designed for children, which maintains the spine in neutral position with the ideal alignment for air flow. Board has four different profiles which guarantee the correct position of the child according to his body size.

• **34024 PEDIATRIC HEAD IMMOBILIZER**
Ensure optimal immobilization of the head of children perfectly fitting their pronounced occipital prominence.

• **34025 PEDIATRIC SPIDER STRAP SYSTEM**

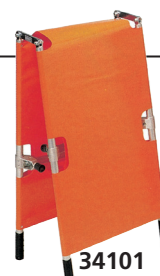

34101




34107




34076




34077

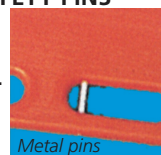

Metal pins

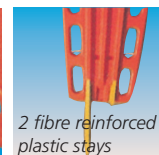

2 fibre reinforced plastic stays

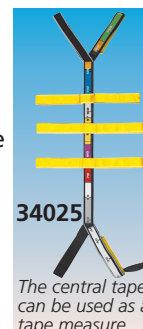

34025

The central tape can be used as a tape measure

TECHNICAL SPECIFICATIONS

	34101	34107	34105	34076	34077	34103	34083	34097	34023
	GIMA 2	GIMA 4	STACKABLE	BLUE	WHEEL	SCOOPING	SPINAL	SPINAL WITH PINS	PAED SPINAL
Open size cm:	203x50x14	221x53x16	195x54x3	188x45x7	188x53x21	min 171x43 max 212x43	184x44.5x6.5	183x46x5	119.5x32x4.5
Folded size cm:	96x50x17	104x19x12	96x54x3	94x48x10	92x53x12	120x43	-	-	-
Weight:	6 kg	5 kg	4 kg	4 kg	8 kg	8 kg	7.6 kg	7 kg	3 kg
Load:	120 kg	150 kg	200 kg	150 kg	150 kg	150 kg	150 kg	200 kg	30 kg