

BASKET STRETCHER

- **34094 TWIN SHELL BASKET STRETCHER - orange**
- **34098 BASKET STRETCHER - orange**

Robust, resistant and flexible, it is ideal for a wide range of rescue operations in mines, at altitude and in water. A revolutionary production technique ensures uniform thickness of the entire length of the stretcher. The shell is in high-density polyethylene and is joined to an aluminium framework for better strength. The handles for transport, integrated in the structure, are placed along the entire perimeter. The grommets for the fixing of the spring catches of the harnessing are in stainless steel. Inside the shell there is a removable mattress, applied with Velcro®, made in closed-cell EPDM, which makes it impermeable to both blood and water. High density polyethylene guarantees shatterproofness and sanitization. It can contain a spine board. Comes with nylon rope along entire perimeter, three 50 mm nylon belts with quick-release buckles, water resistant and closed-cell mattress, footrest in pre-formed polyethylene with position regulator.

Twin shell basket stretcher

It can be transversally divided in half to improve the storage and transport. To ensure a perfect hold between the two shells, have been inserted four nylon cones and a security system protected by a protuberance placed in the mould itself.

• 34099 LIFTING BRIDLES

Adjustable, can be used for horizontal and vertical patient lifting. Manufactured in durable woven nylon with sound steel attachment hardware and brass locking mechanism.

• 34100 FLOTATION SYSTEM

• 34096 CARRYING BAG for 34094




34098


TECHNICAL SPECIFICATIONS

	34098	34094
Size (LxWxH):	215x65x19 cm	118/215x64x19 cm
Load capacity:	280 kg	280 kg
Weight:	13.5 kg	17 kg

34098+34100

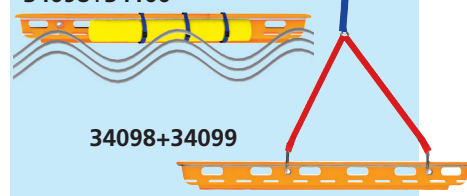

34098+34099

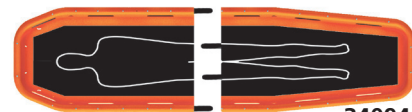

34094

WHEELCHAIR STRETCHERS

- **34060 WHEELCHAIR STRETCHER - 2 wheels**
- **34074 WHEELCHAIR STRETCHER - 4 wheels**

Specially designed for confined areas such as narrow hallways or multiple landing stairways; anywhere access is seriously limited. Made in light aluminium with waterproof orange sheet. Two rear stationary wheels Ø 12 cm (34060) or 4 wheels, 2 with brake (34074) two back handles and four front handles, 2 adjustable 30 cm, facilitate transport on most floor types. 2 quick lock safety belts immobilize patient during transport. Locking safety latch maintains chair position during use.

TECHNICAL SPECIFICATIONS

Folded size:	91x52x19.5 cm	Open size:	91x52x71 cm
Wheelchair size:	- back height: 91 cm - width: 52 cm		
	- seat height: 48 cm - width: 52 cm		
Loading:	150 kg	Weight:	8.5 kg

Foldable
down handleFoldable
up handle

CE



34060



34060

Adjustable
handle (+30 cm)

34074

FOLDING STAIR CHAIR AND EVACUATION CHAIR

- **34058 EVACUATION CHAIR - orange/chrome**

Designed and built for emergency patient removal from confined spaces. Folded it can be stored in most ambulances. Manufactured from chromed steel with brakes on the front wheels. Seat cover made of vinyl coated nylon, easily washable. Folds in one movement, equipped with two quick-release patient restraints.

- **34059 FOLDING CHAIR - black/yellow (orange/chrome on request)**

A cost effective choice for moving patients through difficult situations. Four non-slip grip handles make lifting the patient safer and easier. Two large rear wheels provide easy manoeuvring. Tough, black, vinyl coated seat cover easy to clean, resists stain and provides bacteriostatic protection. Supplied with a set of 3 plastic buckle belts. Made in Italy

- **34061 SKID EVACUATION CHAIR - with braking slide**

Skid-Ok a new ultra compact evacuation chair with innovative aluminium frame and tubular trellis swing fork. New concept materials and components make the chair lightweight and easy to manoeuvre to offer top performance in any scenario.

Reclining armrests with rubber finish improve comfort. The innovative self braking slide system offers the operator total control during the stair descent. Manufactured in welded aluminium with polished finish.

Adjustable seat strengthened with Spentex®.

Skid-OK is supplied with 2x50 mm polypropylene belts. Made in Italy.

TECHNICAL SPECIFICATIONS

	34059 FOLDING	34058 EVACUATION	34061 SKID
Height (open/closed):	92 cm	98/78 cm	154/104 cm
Width:	48 cm	41 cm	52 cm
Depth with retracted handles:	59 cm	-	
Depth with extended handles:	125 cm	open 80 cm	
Depth when folded handles:	19 cm	closed 31 cm	
Weight:	10 kg	12 kg	10 kg
Load:	165 kg	116 kg	150 kg

CE



34058

Height
open/closed
154/104 cm

34061



34059