

ANTI-DECUBITUS AIR MATTRESSES TO PREVENT SORES

- **28482 MASTER PUMP SYSTEM** - with pressure adjustment
- **28483 KIT MASTER (AIR MATTRESS 28492 + MASTER PUMP 28482)**

Made of stiff and strong ABS, it may be placed on the ground, as it has rubber feet, or hung at the edge of the bed.
A control handle makes pressure regulation easier: air amount is adjustable to suit patient weight from 40 to 120 kg.
The mattress should inflate totally in about 20 minutes.
Air output 2.5 l/min. Master pump can be used with bubble or interchangeable section mattress.
Flame retardant PVC.

TECHNICAL SPECIFICATIONS FOR MASTER PUMP - RINT PLUS PUMP

Size: 330x150x120 mm	Suitable for mattress: bubble or section
Weight: 1.95 Kg	Regulations: IEC-601-1 class II Type B
Operating voltage: 230 V - 50 Hz	Made in Italy
on request 110 V - 60 Hz	User manual in English Spanish, French, Italian, German
Power Consumption: 8 W - CLASS II	
Working pressure: 60-135 mm/Hg	



Low noise, long lasting pump made with high quality materials

28482



28492

**MADE IN ITALY
0.4 mm PVC THICKNESS**



28493

Patient Weight	Pressure
50 kg - 110 lbs	40 mm/Hg
50-80 kg - 110-176 lbs	70 mm/Hg
80-150 kg - 176-330 lbs	100 mm/Hg

- **28492 AIR MATTRESS WITH BUBBLES**

This air mattress is divided in bubble sections which inflate and deflate alternatively latitudinally through the alternative compressor. The continual movement given by the pressure pumped inside is useful to prevent minor bedsores problems.

- **28493 AIR MATTRESS - with interchangeable sections**

It has 20 sections which can be inflated or deflated alternatively. Easy and practical to use, it is ideal for patients with minor bedsores problems. Elements are jointed in-between and fixed to a PVC sheet.

- **28494 INTERCHANGEABLE CELL - spare for 28493**


Press studs which attach each section to the supporting sheet allow an easy exchange in case of damage.

- **28495 RINT PLUS AIR MATTRESS - interchangeable cells**
Rint Plus Mattress is made of a strong and soft PVC, realized with 17 big interchangeable bars. The integrated CPR function (Cardio Pulmonary Reanimation) allows the immediate emptying of the 17 bars. It is supplied with bi-elastic, liquid proof and permeable to water-steam cover. Made in Italy.

- **28498 RINT PLUS PUMP SYSTEM - with pressure adjustment**
Same as Master pump but with special valve for use with Rint Plus Mattress.

TECHNICAL SPECIFICATIONS MATTRESSES

	28492	28493	28495
Total size:	190x85 cm	190x85 cm	200x90x10.5 cm
Inflatable size:	160x85 cm	160x85 cm	170x90 cm
Weight:	3 Kg	5.2 Kg	5.5 Kg
PVC thickness:	0.4 mm	0.4 mm	0.4 mm
Number of elements:	-	20	17
Suitable for persons up to:	80/90 kg	100 kg	120 kg

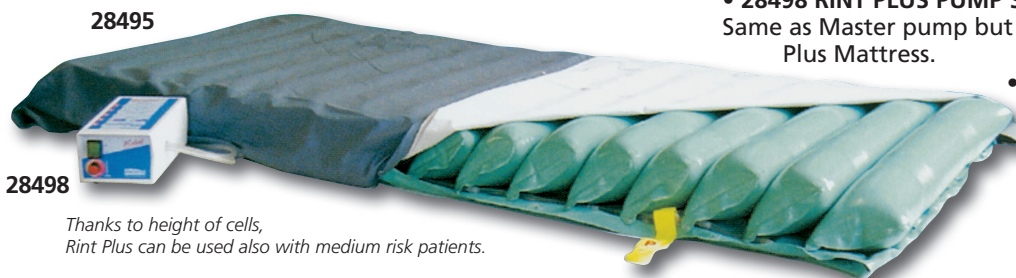

28498

Thanks to height of cells, Rint Plus can be used also with medium risk patients.

- **28496 INTERCHANGEABLE CELL - spare for 28495**
- **28497 COVER - spare for 28495**

ANTI-DECUBITUS PVC WATER MATTRESS

- **28499 WATER MATTRESS**

It is made in PVC with a stabilising inflatable external structure to prevent excessive water "wave motion" and of 3 independent units. The water mattress distributes the total pressure on the whole supporting surface, making this less than the average arterial pressure, avoiding the exposure of the skin to pressure risk. Easy to use, it is placed over the usual mattress and it is washable. It is sold with a manual pump, 2 disinfectants and an inflation tube.
Size: 190x85 cm. Weight 4 kg. PVC thickness 0,4 mm



28499

Each water section measures 60x70 cm

PVC rigid edge inflatable with air