

PILL CRUSHING DEVICE

- **25872 SILENT KNIGHT HOSPITAL PILL CRUSHING DEVICE**
 - **25873 PILL CRUSHER POUCHES - 1,000 pouches (20 packs of 50)**
- Professional, silent crusher for hospital to effortlessly and easily reduce pills in fine powder in few-seconds. Pouches avoid any cross contamination and eliminate loss of medication. Silent Knight is supplied with 50 pouches.



Step 1: place pills into pill crusher pouch.
Step 2: press handle down to crush pills.
Step 3: press the handle down to activate grinding movement. Press handle down several times to grind pills into fine powder.
Step 4: remove the pouch.




ELECTRIC PILL GRINDER

- **25850 SEVERO ELECTRIC PILL GRINDER**
- Professional silent pill grinder for fast and effective crushing of medication in hospitals, nursing homes and home care organizations. It avoids exposure to dust or cross contamination. The device works through autonomous grinding cycles via combined compression and rotation, through which it is possible to grind multiple standard pills and coated pills per cycle. Once the medicine is put into the 30 ml special cup, it takes only 8 seconds (for regular pills) to obtain it as powder. Pulverized pills are ideal for people who have trouble swallowing pills or to be administered via feeding tubes. Supplied with stainless steel grinding head, 30 cups, 30 cellophane grinding sheets, AC/DC adapter and GB, FR, IT, PT manual. Available DE, NL on www.gimaitaly.com.

MADE IN NETHERLANDS

**FROM PILL TO POWDER
IN FEW SECONDS**

TECHNICAL SPECIFICATIONS	
Cup volume:	30 ml
Grinding cycle time:	8 seconds for regular pills
Sound level:	52dB
Battery:	rechargeable lithium-ion, 12 V 3 Ah
Battery last:	500 crushing's (three days of continuous working)
Charging time:	2 hours
Size - weight:	80x70x217 mm - 1.5 kg

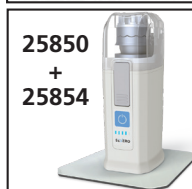

- **25851 SEVERO MEDICINE CUPS - box of 2100**
 - **25852 SEVERO CUP LIDS - box of 2000**
 - **25853 SEVERO CELLOPHANE GRINDING SHEETS - box of 2100**
 - **25854 SEVERO MAGNETIC STABILITY PLATE**
- Powder coated metal plate and magnetic sticker to be inserted underneath to stabilize Severo on medication trolleys.

PILL DISPENSER - WIRELESS

- **25871 WIRELESS PILL DISPENSER**
- Automatic wireless pill dispenser to manage and record pills' dosage data with Healbox App and cloud server. It lets patients and caregivers to program dosages up to a maximum of 28 consecutive days. It also allows to give right dose to right patients at right time. The device has audio and visual alerts up to 9 alarms per day, which sound for 30 minutes or until the dispenser is flipped. It automatically rotates when it's time for the next dose. 28 compartments with 24 aspirin sized pills load capacity. Lid can be locked with metal key. For hospital, community service or personal use. Supplied with 1 medication tray, 1 set of dosage rings, 1 metal key, 4xAA batteries and multilingual instructions (GB, FR, IT, ES, PT)

TECHNICAL SPECIFICATIONS			
Compartments:	28 slots	Connection:	10 meters limited range
Capacity:	24 Aspirin-sized pills	Battery:	4xAA batteries
Material:	ABS	Size:	225x h 68 mm
LED display size:	60x30 mm	Weight:	600 g



GRASPING TOOLS

- **28960 GRASPING TOOL - 3 WAYS - 66 cm**
- Aluminium body, PP handle. Three grasping ways: steel hook, magnetic tip, rubber head. Multilingual box and instructions: GB, FR, IT, ES.

- **28963 GRASPING TOOL - foldable**
- Extended reach grasping tool helps to eliminate bending or climbing for hard-to-reach items. Length: 80 cm, foldable in 2. Locking tab holds the item in place. Multilingual box and instructions: GB, FR, IT, ES.