

ANTI-DECUBITUS CUSHIONS


MADE IN FRANCE



28530



28534

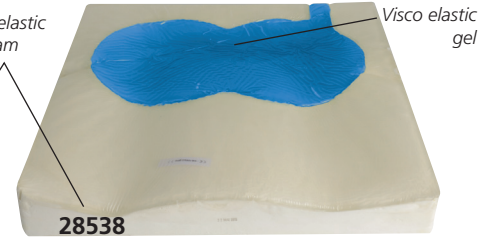

28538



PREVENTION FOR LOW TO MEDIUM
DECUBITUS RISK



PREVENTION FOR HIGH
DECUBITUS RISK



PREVENTION FOR HIGH TO VERY HIGH
DECUBITUS RISK

GIMA line of seat cushions is designed to prevent decubitus risk for patients on wheelchair or permanently on chairs. All cushions in Gel or Foam are supplied with a waterproof PU cover with zipper and handle.

• 28530 BULGEL CUSHION 41x41x3 cm

• 28531 BULGEL CUSHION 43x43x3 cm

This anti-decubitus cushion is made of a flame-resistant pocket including a PE foam core covered by an aqueous gel.

• 28534 VISCO-MOUSS CUSHION 41x41x7* cm
• 28535 VISCO-MOUSS CUSHION 43x43x7* cm

Pre-shaped moulded cushion made of viscoelastic foam with memory effect. Memory foam is pressure sensitive and moulds quickly to body shape, returning to its original shape once the pressure is removed. Anatomical shape with raised lateral edges and front-back slope improves positioning, stability and comfort when seated. Foam density: 86 kg/m³

*Minimum 7 cm, maximum 8 cm

• 28538 GEL AIR 2D CUSHION 41x41x7.5 cm
• 28539 GEL AIR 2D CUSHION 43x43x7.5 cm

Moulded viscoelastic foam base (i.e. Visco Mouss) with a polyurethane gel insert. The viscoelastic fluid gel (blue pad) ensures maximum pressure relief of vulnerable bony areas while the regular gel and the foam increase the contact surface and maximize pressure relief. The cushion's shape ensures a good posture and prevents from sliding forward. Foam density: 86 kg/m³

AIR CUSHIONS

• 28544 GIMA AIR CUSHION - 40x40x6 cm
• 28546 GIMA AIR CUSHION - 45x45x6 cm
with polyester cover, pump and repair kit
- single compartment
- air cells, 64 or 81, prevent pressure sores
- maximizes pressure distribution and blood flow

- conforms to the body's shape
- adjusts to body changes over time.

Max load: 150 kg

Multilingual manual: GB, FR, IT, ES, PT, GR, PL, DE, SE, RO.



Valve to adjust air
and softness of cushion
Black rope for carrying the cushion
PREVENTION FOR HIGH TO VERY HIGH
DECUBITUS RISK


HEEL PROTECTOR

• 28548 HEEL PROTECTOR - 50x22x16 cm

It relieves the heel pressure and prevents decubitus ulcer by increasing the contact surface.

Anatomical shape visco elastic foam including a bactericide and fungicide agent for better hygiene and disinfection.

Density: 90 kg/m³

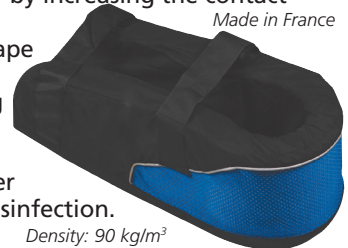

Made in France

ANTI-DECUBITUS CUSHIONS AND NECK PILLOWS

MADE IN ITALY

Products to prevent low decubitus risk, made of silicone fibre, coated with 100% preshrunk cotton cover. Main property of

silicone fibre is to maintain original shape also after strong pressure. This allows to maintain comfort and proper weight

distribution in the time. Breathable cotton cover provides maximum cool comfort. Washable at 30-40°C.



28549

Size: 43x45 cm



28550

Size: 45x45 cm



Size: 45x45 cm



Size: 45x45 cm


• 28549 RING SHAPED CUSHION

• 28550 SHAPED CUSHION WITH HOLE

• 28551 3-SECTION CUSHION


• 28552 CUSHION WITH CONCAVE CENTRE

All cushions and pillow are supplied with multilingual instructions:
GB, FR, IT, ES, PT, DE

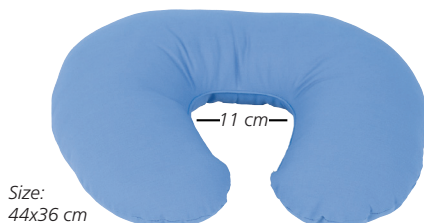

Size:
44x36 cm

• 28555 TRAVEL NECK PILLOW

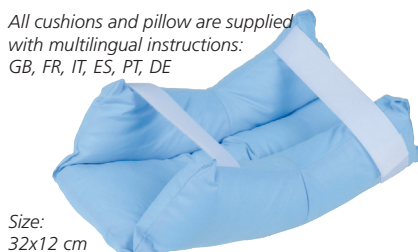

Size:
32x12 cm

• 28556 HEEL CUSHION



Size:
32x15 cm

• 28557 HEEL CUSHION PLUS