

NEONATAL CRADLE

• 43500 NEONATAL CRADLE with trolley

Transparent unbreakable plastic cradle with round edges and ventilation holes. Epoxy powder painted mainframe made of curved and welded tubular steel. Trolley shape allows horizontal and Trendelenburg inclinations and bedside positioning. 4 synthetic rubber castors Ø 75 mm, two with brake. Weight: 15 kg. Max load: 10 kg. Manual: GB, FR, IT, ES. Delivered disassembled.

• 43501 CRADLE - spare

• 27685 MATTRESS 63x37x5.5 cm - fireproof
With removable fire retardant cover.

• 27686 COVER - water proof



Dimensions:
83.5x52.5xh 90 cm

MADE IN ITALY



27400

SWADDLING TABLE

• 27400 SWADDLING TABLE - beige

Baby changing table with Itallumag tubular frame with 2 shelves. Padded worktop 115x60 cm in soft washable fireproof (Class 1 IM) simulated leather and lower shelf 100x40 cm, in laminated chipboard panel.

Max load: 50 kg

Size: 115x63xh 85 cm

Weight: 28 kg

Delivered disassembled.

MADE IN ITALY

BABY CRADLES

• 43502 BABY CRADLE - light blue

• 43503 BABY CRADLE - pink


Cradles with stylish and ergonomic design to provide a comfortable experience for babies. Transparent baby basin made of plastic. Wear-resistant, antiskid, shock absorption double side castors Ø 76 mm (3") can freely rotate 360°. Easy to operate handle to control gas spring. 12° reverse Trendelenburg. Adjustable height by hydraulic pump. Size: 875x525xh 650 to 900 mm. Packing size: 930x560x850 mm. Supplied with baby basin.

• 43504 MATTRESS+LINER - optional

Antimicrobial, fire-retardant, polyurethane mattress 40x76x2 cm. Zip L Tekma liner. Waterproof and antimicrobial.



43502

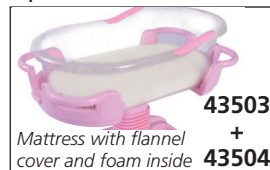

Mattress with flannel cover and foam inside
43503 + 43504



Adjustable height



+12° Reverse Trendelenburg

PEDIATRIC BED

• 27399 PAEDIATRIC BED 1-4 years - ivory

Epoxy powder painted mainframe made of steel tube Ø 30 mm, ready to be hooked into the mattress platform. Mesh mattress platform made of oval tube 40x20 mm, fixed to the mainframe with removable bolts. Side rails 62 cm h. sliding on chromed steel guides, with safety device. Legs with 4 castors Ø 80 mm, two with brake. Weight 35 kg. External measurement 140x70xh 130 cm. Inner size 130x65xh 63 cm.

• 27687 MATTRESS 125x62x10 cm - fireproof

• 27688 COVER - water proof



27399

MADE IN ITALY

CHILDREN BEDS - 2 CRANKS

• 43506 CHILDREN BED - 2 cranks - light blue

• 43507 CHILDREN BED - 2 cranks - pink

Beds for children with frame made of sheet metal bending and welding. Plastic parts made of ABS material and epoxy painting. Bed ends are made of PP blow moulding to protect bed and patient.

• 43508 MATTRESS - 1770x720x100 mm - optional

Mattress made of polyurethane foam, density 30 Kg/m³.

Radiotransparent, antimicrobial, latex and CFC free.

Odourless and fire retardant. Supplied with 100% polyester liner, washable.

Guardrail

Aluminium alloy liftable guardrail with 8 cm distance among protective tubes. 3 gears to adjust height, anti-pinch hands to protect children.

43507

Hand control

Manual cranks to control back and leg liftings

Shoes' holder

Utility and shoes holder

TECHNICAL SPECIFICATIONS

Size:	1880x900xh 1260 mm
Bed height:	585 mm
Mattress platform:	1770x720 mm
Back-rest adjustment:	0-75° (± 10°)
Knee-rest adjustment:	0-35° (± 10°)

IV Pole
Stainless steel IV pole, bearing 15 kg



Castors Ø 12.7 cm (5") castors with brakes

Platform

Concave vent hole design of the panel allows strong air permeability and can effectively prevent bedsores.

STANDARD ACCESSORIES

1 set full length side rail	4 buffer air springs
1 stainless steel IV pole bearing 15 kg	Shoe holder
4 IV pole provision (2 on each head)	
2 drainage hooks	