


MANUAL PATIENT BEDS WITH 3 OR 4 JOINTS

Guardrail

Aluminium alloy side rail, 6 staff, collapsed type with safe lock

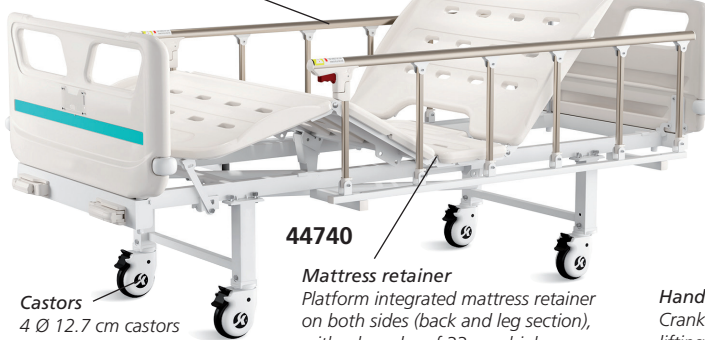

44740

Mattress retainer

Platform integrated mattress retainer on both sides (back and leg section), with a boarder of 33 mm high

Castors
4 Ø 12.7 cm castors

Hand control
Crank to control lifting

- **44740 3 JOINTS PATIENT BED - 2 cranks - castors**
- **44742 4 JOINTS PATIENT BED - 3 cranks - variable height - castors**

Manual patient beds with 3 or 4 joints adjustable by means of 2 or 3 cranks. PP head and foot board equipped with bumpers wheel Ø 65 mm, preventing collision. PP bed platform, divided in 4 or 5 sections with integrated mattress retainer. Plastic parts and epoxy painting.

Plastic parts and epoxy painting.

TECHNICAL SPECIFICATIONS

Size:	2130x995 mm
Height:	500 mm (44740), adjustable 420-730 mm (44742)
Mattress platform:	1925x900 mm
Max. patient load:	200 kg
Movement:	2 cranks (44740), 3 cranks (44742)
Back-rest adjustment:	0-70° (± 5°)
Knee-rest adjustment:	0-38° (± 5°)

STANDARD ACCESSORIES

2 foldable side rails
IV pole 4 hooks
4 IV pole provisions
2 drainage hooks

IV Pole

Stainless steel IV pole, bearing 15 kg




44742

Wall bumper wheel
Equipped with 4 bumpers to protect the bed and wall



44742 Height adjustment

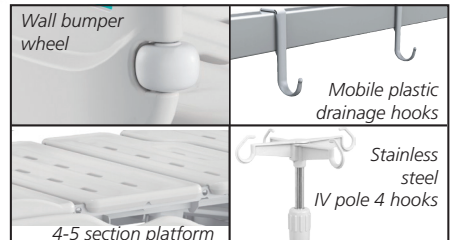

Wall bumper wheel

Mobile plastic drainage hooks

4-5 section platform

Stainless steel IV pole 4 hooks

ELECTRIC PATIENT BEDS

Mattress retainer


Platform integrated mattress retainer on both sides (back and leg section) with a boarder of 33 mm high

Guardrail

Aluminium alloy foldable side rail, 6 staff, collapsed type with safe lock

IV Pole

Stainless steel IV stand, bearing 15 kg

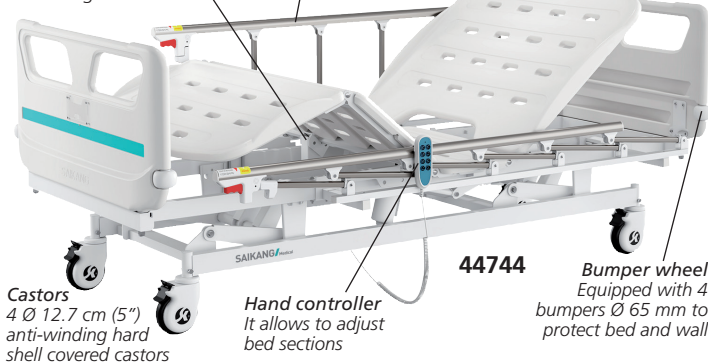

44744

Castors
4 Ø 12.7 cm (5") anti-winding hard shell covered castors

Hand controller
It allows to adjust bed sections

Bumper wheel
Equipped with 4 bumpers Ø 65 mm to protect bed and wall

- **44744 ELECTRIC PATIENT BED - adjustable height**
 - **44746 ELECTRIC PATIENT BED with TR and RTR - adjustable height**
- Electric patient beds with integrated PP bed platform and electric motor to facilitate different movements of bed sections. Hand controller allows nurses to control standard functions. Both models have an American motor with UL/ROHS/EN standard. Plastic parts and epoxy painting. PP head and foot board equipped with four bumpers wheels Ø 65 mm. It is provided with auto-regression 10+9 cm (44746 only) to minimize the pressure and shear forces in patients' pelvic area and to reduce development of pressure ulcers.

STANDARD ACCESSORIES

- 1 hand controller
- 1 set foldable side rails (44744 only)
- 1 set PP side rails extra height (44746 only)
- 1 IV stand with 4 hooks
- 4 IV pole provisions
- 2 drainage hooks
- 4 angle indicators (44746 only)

Angle indicator


Integrated angle indicator to clearly show back-rest and Trendelenburg position

Guardrail

Tuck away side rail with integrated buffer. Handle allows patient to get up and off the bed easily

Hand controller

Controller with flexible spiral cable for easy operation. Nurse can control all standard functions

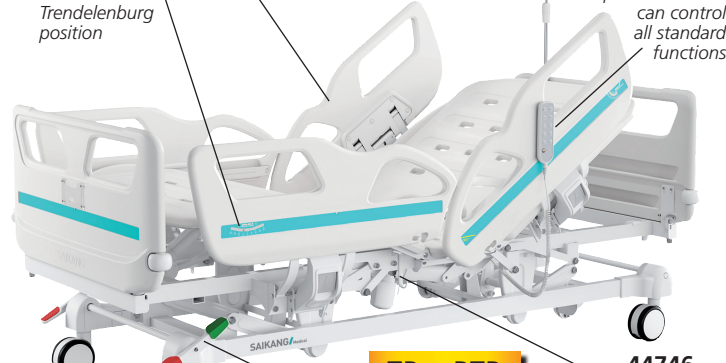


TR + RTR

44746

Castors
4 Ø 12.7 cm (5") double side castors with central locking. Anti-winding hard shell cover

Brake system
Brake spring retains quality even after 10.000 uses

Drainage hooks
Plastic drainage hooks integrated inside the siderails

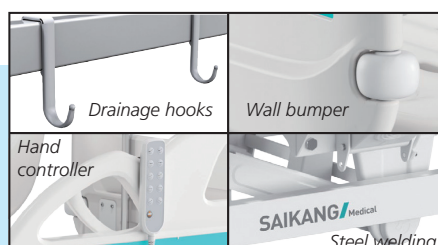

Drainage hooks

Wall bumper

Hand controller

SAIKANG Medical
Steel welding

TECHNICAL SPECIFICATIONS

	44744	44746
Size:	208.5x99.5 cm	208.5x102.5 cm
Height adjustment:	from 42 to 73 cm	from 45.5 to 74.5 cm
Mattress platform:	192.5x90 cm	192.5x90 cm
Safe working load:	200 kg	250 kg
Back-rest adjustment:	0-75° (± 5°)	0-75° (± 5°)
Knee-rest adjustment:	0-46° (± 5°)	0-32° (± 5°)
Trendelenburg:	-	0-12°
Reverse Trendelenburg:	-	0-12°
Angle indicator:	-	0-90°
Auto-regression:	-	10+9 cm
Position line:	-	30°
Motor: life:	20.000 movements	20.000 movements
noise:	under 50 dB	under 50 dB
Voltage:	100-240 V - 50/60 Hz	100-240 V - 50/60 Hz