

PATIENT TROLLEYS - STRETCHERS



Removable top
on the trolley

Adjustable backrest
up to a 60°

MADE IN ITALY

Multilingual instructions:
GB, FR, IT, ES

27801



• 27800 PATIENT TROLLEY

• 27801 PATIENT TROLLEY - removable top

Made of chrome plated steel tube. Comes complete with mattress made of washable self extinguishing plastic material. Mounted on 4 castors Ø 100 mm two of which with brakes. Size: 180x54xh 80 cm. Load: 110 kg.

GIMA code	PATIENT TROLLEY ACCESSORIES - not included
27810	Side rails (couple)
27811	I.V. Stand
27812	Steel basket only for 27800
27813	Oxygen bottle holder

PROFESSIONAL WARD STRETCHER

MADE IN ITALY

• 44760 WARD STRETCHER

Main features:

- **Platform:** all-welded mainframe made of semi-oval steel tube 50x25x2 mm.
 - **Undercarriage:** made of steel tube Ø 30x1.5 mm and oval steel tube 40x20x1.5 mm
 - **Push/Pull handle:** made of steel tube Ø 30x1.5 mm.
 - **Surface finish:** epoxy gloss finish standard grey 7040
 - **Mattress platform section:** two section mattress platform fitted with grids made of 5 mm Ø steel rods and square mesh 50x50 mm.
 - **Backrest:** square cross section 25x25x1.5 mm, adjustable by means of a rack device (5 positions).
 - **Castors** - four castors Ø 200 mm made of grey synthetic rubber two of which with independent brake.
 - **Bumpers:** Four bumpers at the corners made of anti-acid and anti-shock grey resin.
- Supplied disassembled.
External size: 202x71xh 90 cm
Wired net height: 78 cm
Load: 135 kg.
Manual: GB, FR, IT, ES.



44760
+
accessories
(not included)

GIMA code	WARD STRETCHER ACCESSORIES - not included
44761	Side rails (couple) with automatic locking device
44762	I.V. Stand with two hooks
44763	Basket - grey painted steel
44764	Oxygen bottle holder* - chrome plated steel
44765	Fireproof simulated leather mattress - class I - 182x64x4 cm

*For oxygen bottle up to 41 cm, Ø 15,5 cm



TRENDELENBURG PATIENT TROLLEY - STRETCHER



Code 27810 side rails
(not included)

Bumpers at
all 4 corners

• 27850 TRENDELENBURG TROLLEY

Made of steel-alloy, oven painted ral 9010, with chrome plated handles.


Adjustable height from 50 to 70 cm, with oleodynamic pump. Trendelenburg up to 9°.

Anti-Trendelenburg up to 14° (with gas pump).

Plastic bumpers at four corners. Maximum load: 120 kg. 4 wheels of Ø 200 mm (2 with brake).

Size: 65x184 (internal) x h 50-70 cm, 72x191 cm (external).

GIMA code	ACCESSORIES
27810	Side rails (couple) 100xh 28 cm
27853	I. V. Stand
27854	Oxygen bottle holder
27855	Mattress - black sky 3 cm thick in flame retardant material - 184x64x3 cm

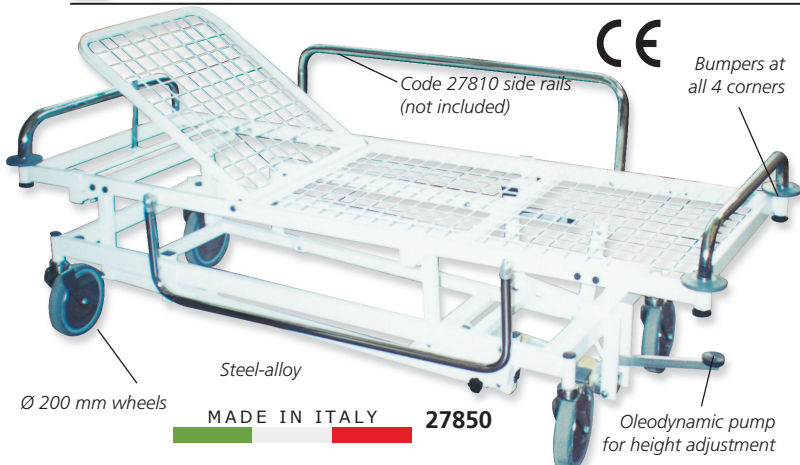

MADE IN ITALY

27850

Oleodynamic pump
for height adjustment