

PROFESSIONAL PATIENT TROLLEY - ADJUSTABLE HEAD SECTION

- **27804 PROFESSIONAL PATIENT TROLLEY** with adjustable head section, side rails and oxygen cylinder holder
- **27805 PROFESSIONAL PATIENT TROLLEY - RADIOTRANSSPARENT FRAME** with adjustable head section, side rails and oxygen cylinder holder

Innovative, stylish comfortable trolley easy to use.
Trolley made of sturdy tubular frame for optimal stability.

Supplied with 4 wheels Ø 200 mm (2 with brake), bumpers, push handle.

Detachable frame with adjustable head section (27804) or fixed frame, 12 mm thickness, with adjustable head section (27805).

Size: 196x76xh 75 cm. Net size: 193x68 cm

Maximum load: 200 kg.

Weight: 70 kg.

- **27809 MATTRESS - 190x68xh 8 cm - blue**

Made of 45% polyester, 55% polyurethane.


Fire retardant, Class 1.

Antistatic and antibacterial.

Density: 30 kg/m³ ±5%

- **27808 IV STAND with 2 plastic hooks.**

RADIOTRANSSPARENT FRAME

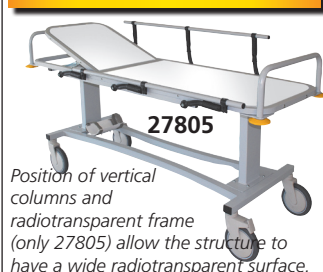

Position of vertical columns and radiotransparent frame (only 27805) allow the structure to have a wide radiotransparent surface.



MADE IN ITALY

ADJUSTABLE HEIGHT PATIENT TROLLEY WITH TR AND RTR

- **27806 PATIENT TROLLEY - with adjustable height, TR and RTR**

Stylish sturdy and comfortable trolley with innovative design.

Adjustable height by means of hydraulic pump with pedal placed on both sides. Lateral sides, easy and fast to fold.

The upper wire mesh is removable to guarantee a better hygiene. Equipped with push handles, 4 bumpers at the corners and a universal hold-cylinder.

Reclinable cushion by two racks easily actuated by the grip handle. Wheels 200 mm diam. with total and simultaneous central brake; one of 4 wheels is equipped with directional brake;

activated by means of 4 lateral pedals.

The joints of the structure rotate on self-lubricating sintered threaded inserts. The trolley can reach the Trendelenburg (-10°) and anti-Trendelenburg (+10°) position very quickly thanks to a handle placed on the upper lateral side: it activates a gas spring (600N) that allows millimetric adjustability. The entire structure, in rectangular and round tubular is epoxy painted in grey and white colour.

Size: 199x76xh 52/87 cm (without mattress).

Load: 130 kg. Weight: 90 kg

- **27809 MATTRESS - 190x68xh 8 cm - blue**

Made of 45% polyester, 55% polyurethane. Fire retardant, Class 1. Density: 30 kg/m³ ±5%. Antistatic and antibacterial.

- **27808 IV STAND with two plastic hooks**



MADE IN ITALY



HYDRAULIC HEIGHT ADJUSTABLE PATIENT TROLLEY - TR and RTR

- **44790 HEIGHT ADJUSTABLE PATIENT TROLLEY - TR and RTR**

Two-section variable-height stretcher with adjustable backrest by gas devices up to 65°. Top made of radiolucent Hpl, steel frame and long leg section made of painted steel.

Possibility of insertion of cylinder basket and four IV poles, two in the head and two in the foot area.

Side bumpers on the frame for protection of the stretcher.

Two collapsible handles on the head side and one on the foot side, all in painted steel, to provide a safe and comfortable grip.

Adjustable height 530 to 890 mm by a large pedal located in the foot area. All steel parts are epoxy powder coated.

Manual Trendelenburg and Reverse Trendelenburg (±18°) by lever located at the foot side of the stretcher.

Sturdy folding steel side rails with compass movement.

Four castors Ø 200 mm with simultaneous brake controlled by centralized pedal, directional ones at the head.

- **44793 WATERPROOF, FIREPROOF, RADIOTRANSSPARENT MATTRESS 6 cm**

Mattress with waterproof, washable, atoxic cover.

Mattress allows cardio massage.

- **44794 I.V STAND 4 HOOKS**

- **44795 OXYGEN BOTTLE HOLDER**

- **44796 STORAGE BASKET**

MADE IN ITALY



RADIOTRANSSPARENT FRAME

Overall size:

Table size:

Max. backrest inclination:

Safe working load:

Norms:

TECHNICAL SPECIFICATIONS

2000x825xh 530-890 mm

1900x630 mm

65°

230 kg

2017/745 EU, Medical Device Class I