


PATIENT TRANSPORTATION TROLLEY - WITH TR AND RTR

- **44792 HYDRAULIC HEIGHT ADJUSTABLE PATIENT TROLLEY ON COLUMNS WITH TR AND RTR - side rails included**
Patient trolley with 1 joint, 3 sections and HPL radiolucent top. Head section continuously adjustable by gas devices up to 80°; back section moved by gas assist and foot section moved by rastomat. Epoxy powder-coated steel frame with multifunctional supports at the corners allowing insertion of up to two accessories per corner such as drip poles, monitor holders or others. Base with protective cover and compartment for gas bottles and objects. Equipped with four 200 mm Ø swivel castors and one with directional lock. Centralised, simultaneous braking system by means of a pedal located on the foot side. Folding handles at foot end. Pair of folding side rails with compass movement made of extruded aluminium alloy with automatic lock and double lowered position. Optional height-adjustable IV pole with 4 hooks. Large foot pedal on both sides of the stretcher to adjust height from 550 to 850 mm. Pedals on both sides for lowering and reaching the Trendelenburg and Reverse Trendelenburg positions (- 18° / + 18°). Control pedal with three positions: 1 free, 1 directional and 1 total lock. Suitable for transport even over long distances.
- **44793 WATERPROOF, FIREPROOF, RADIOTRANSSPARENT MATTRESS 6 cm - optional**
- **44794 I.V STAND 4 HOOKS - optional**
- **44795 OXYGEN BOTTLE HOLDER - optional**
- **44797 STORAGE BASKET - optional**

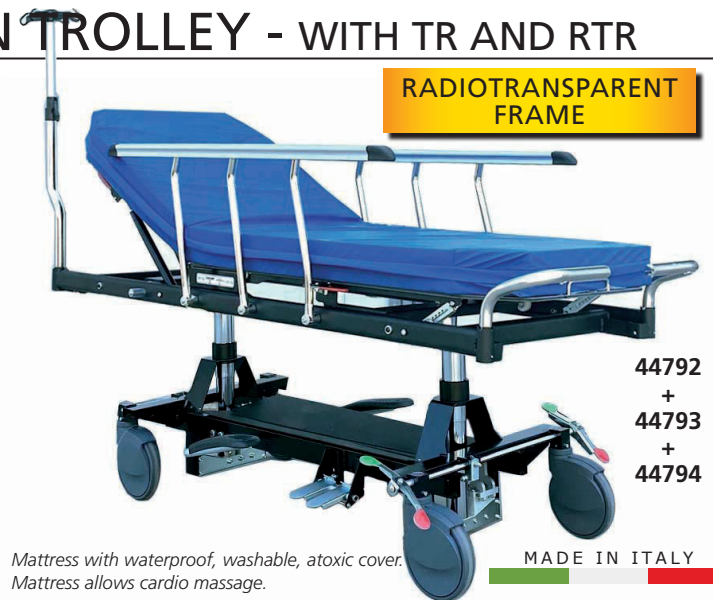

44792
+
44793
+
44794

Mattress with waterproof, washable, atoxic cover.
Mattress allows cardio massage.


MADE IN ITALY

TECHNICAL SPECIFICATIONS

Overall size:	2060x820xh 550-850 mm	
Size - weight:	1950x650 mm - 105 kg	
Max. backrest tilt:	80°	
Trendelenburg / RTR:	- 18° / + 18°	
Safe working load:	320 kg	Max. patient load: 260 kg
Norms:	MDR 2017/745; Medical Device Class I	

HYDRAULIC HEIGHT ADJUSTABLE PATIENT TROLLEY - WITH TR AND RTR

- **44782 HYDRAULIC HEIGHT ADJUSTABLE PATIENT TROLLEY with TR and RTR**
Metal patient trolley equipped with foldable aluminium alloy guardrails and two hydraulic pumps which allow height adjustment, Trendelenburg and reverse Trendelenburg. Four double castors Ø 20 cm, lockable by pedals on both sides. Waterproof PU mattress. Four detachable bumper wheels Ø 74 mm to protect bed and patients during transportation. Stainless steel IV pole with plastic hooks. X-ray translucent bed platform. Central wheel to facilitate direction control. Multilingual manual: GB, FR, IT, ES, PT, DE, GR, RO.

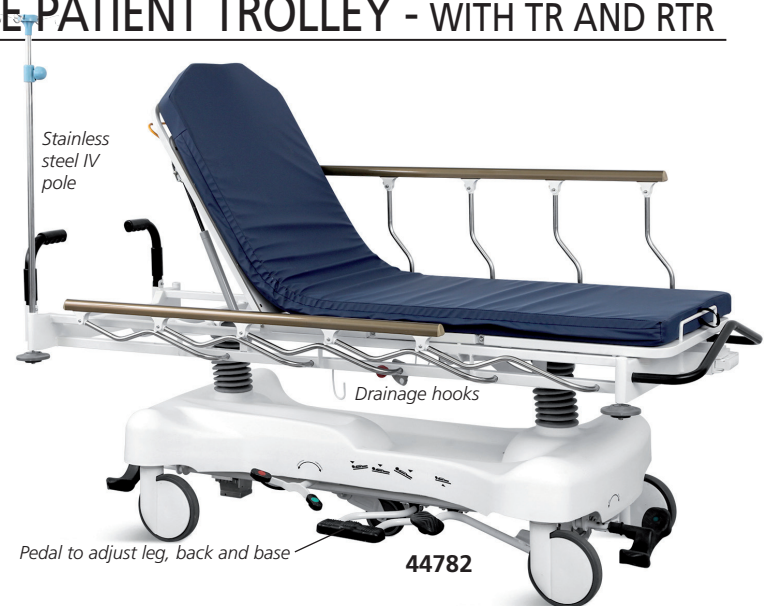

44782

TECHNICAL SPECIFICATIONS


Size:	2100x850 mm	Knee-rest adjustment:	0-45°
Height adjustment:	630-930 mm	Trendelenburg:	0-15°
Back-plate adjustment:	0-75° (± 10°)	Reverse Trendelenburg:	0-15°

STANDARD ACCESSORIES

1 set emergency push handle	4 IV pole provisions	1 set IV pole
1 set folding side rails	1 mattress	

HEIGHT ADJUSTABLE PATIENT TROLLEY

- **44780 HEIGHT ADJUSTABLE PATIENT TROLLEY**
Height adjustable patient trolley made of metal. Both bed frame and chassis have 30x50 mm square tube to make stretcher stronger and more stable. Bed platform and guardrail are made of high quality PP material. Back lift makes patient more comfortable. Equipped with strips to prevent patients from possibly sliding off the stretcher and 2 cm thickness mattress. Directional wheel and central locking pedals. Multilingual manual: GB, FR, IT, ES, PT, DE, GR, RO.

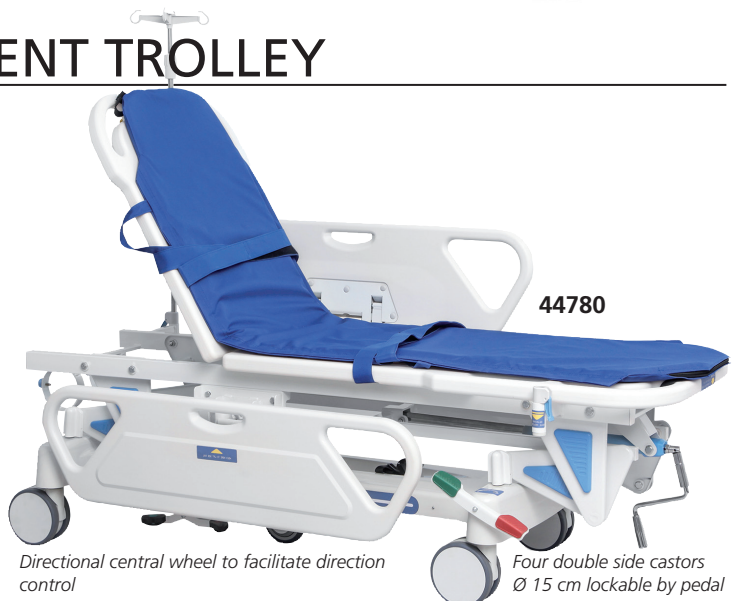

44780

TECHNICAL SPECIFICATIONS

Size:	1900x640 mm
Height adjustment:	550-850 mm
Back-plate adjustment:	0-70° (± 10°)

STANDARD ACCESSORIES

1 set IV pole	1 transfer mattress
2 IV pole provisions	2 strips
1 oxygen bottle holder for 2 l bottle	1 set ABS folding side rails