

MULTIFUNCTIONAL TROLLEYS WITH CENTRAL BRAKING SYSTEM

• 45711 MULTIFUNCTIONAL TROLLEY 65x88xh 100 cm

• 45713 DRESSING TROLLEY 65x88xh 100 cm with dispenser

Solid and robust construction with steel sheet painted with antimicrobial epoxy powders and technopolymers for maximum lightness, maneuverability and ease of use.

Upper part: large support surface in PST with ergonomic and integrated push handle. Push handle has a central mechanism for locking and unlocking wheels' brakes.

Front part: retractable worktop support with "push and pull" opening.

Four modular drawers with soft self-closure, ball bearings and central locking with key. Drawers are made of painted steel.

Left side: recessed compartment and retractable worktop made of technopolymer plastic; 10 L waste bin with knee opening and technopolymer catheter holder.

Right side: recessed compartment with three tilting pockets made of technopolymer material.

Supplied with four Ø 100 mm twinned castors with automatic locking system.



Left side




Central braking system



MADE IN ITALY

45711

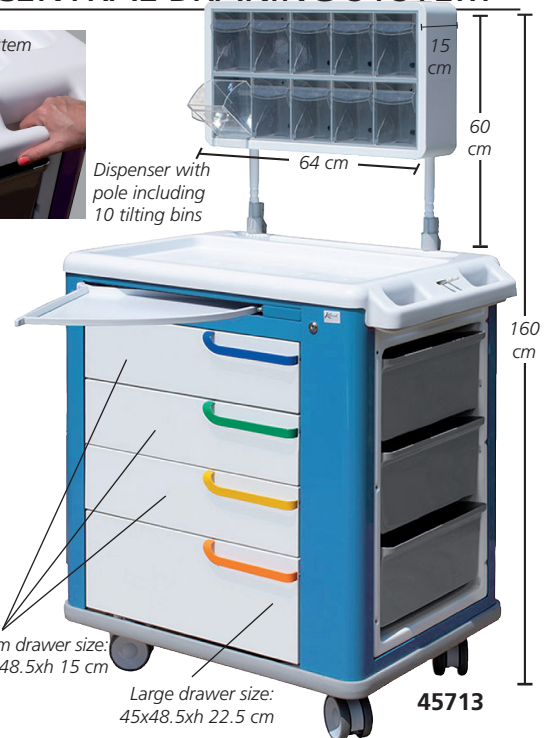

Medium drawer size:
45x48.5xh 15 cm

Large drawer size:
45x48.5xh 22.5 cm

45713

KS MULTIFUNCTIONAL TROLLEYS

MULTIFUNCTIONAL TROLLEYS: BASIC, DRESSING, EMERGENCY

Professional trolleys in 2 sizes (small and standard) and 3 models (Basic, Dressing and Emergency).

- antimicrobial epoxy painted body with seal lock (Emergency) or key lock (Basic and Dressing)
- top surface with 4-sided rim
- top surface with built-in compartment

MADE IN ITALY

- ball-bearing full-extension sliding drawer with self-closing system
- built-in handle
- 4 swivelling castors Ø 125 mm - 2 antistatic (only 45688/9) and with brake
- lateral waste container 15 l with lid and catheter holder (Dressing and Basic only)
- size: small 760x714xh 1,022 mm
standard 911x704xh 1,075 mm

Lateral side with waste container 15 l with lid and catheter holder



45680



Dispenser with pole including 10 tilting bins

45685

Back with large oxygen tank compartment and stainless steel holding rail with hooks

Lateral side with waste containers 15 l with lid

45689



45681

2 small drawers: 600xh 75 mm
2 medium drawers: 600xh 151 mm
1 large drawer: 600xh 228 mm

Writing surface

Lateral container unit with 4 supply trays



45684

3 medium drawers:
450xh 151 mm
1 large drawer:
450xh 228 mm

Writing surface

Lateral container unit with 4 supply trays



45689

2 small drawers: 600xh 75 mm
2 medium drawers: 600xh 151 mm
1 large drawer: 600xh 228 mm

- 45680 BASIC TROLLEY - small
- 45681 BASIC TROLLEY - standard

- 45684 DRESSING TROLLEY - small
- 45685 DRESSING TROLLEY - standard
Dressing trolley is supplied with upper dispenser (10 tilting bins).

- 45688 EMERGENCY TROLLEY - small
- 45689 EMERGENCY TROLLEY - standard
Supplied with rotating defibrillator/monitor tray, height adjustable I.V. pole with 2 hooks. Medical Device Class I.