


1 OR 3 ENGINES GYNAECOLOGICAL CHAIRS AND STOOLS

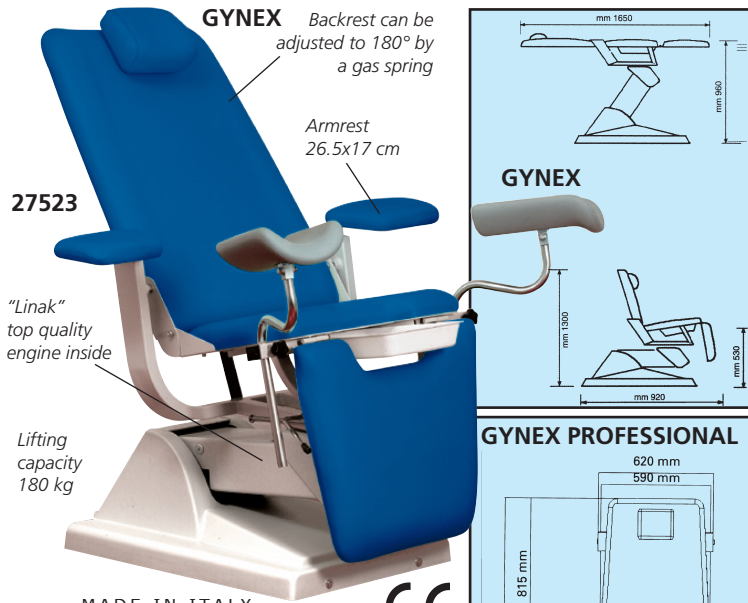

GYNEX - GYNAECOLOGICAL CHAIR
GYNEX: a chair offering a unique combination of functional, aesthetic and working characteristic.

Thanks to its self-balancing seat and backrest, GYNEX can be adjusted to any build for the greatest patient comfort. The chair is provided with an electromechanical lifting device, having a power supply at 24 V.

The electric actuators enable the movement up and down from 57 cm off the floor ("0" position) up to a maximum of 96 cm. The backrest is supplied with adjustable headrest and paper roll holder. The legrest can be independently adjusted and it has an automatic-locking device.

The chair can be flattened to a bed position for examination activity, the armrest can be removed and independently oriented.

The electric control of the chair for upright position is operated by means of a foot pedal.

TECHNICAL SPECIFICATIONS

	GYNEX	GYNEX PROFESSIONAL
Power supply tension:	230V - 50 Hz (110V - 60 Hz on request)	24 V
Motor tension:	24 V	24 V
Lifting motor push:	6,000 N	6,000 N
Backrest push:	800 N (gas spring)	2,500 N (motor)
Legrest motor push:	-	2,500 N
Lifting capacity:	180 kg	180 kg
Maximum power:	90 VA	90 VA
Fuses protection:	1.6 A	1.6 A

GYNEX CHAIRS AND STOOLS

	Light blue Parigi	Green Toronto	Blue Chicago	Orange Rotterdam	Metal Apricot	Metal Sea blue	Colour on request
Gynex professional	27514	27512**	27513	27512**	27515	27512**	27512**
Gynex	27521	27522	27523	27524	27525	27526	27520**
Stool with ring	27531	27532	27533	27534	27535	27536	27530**
Stool without ring	27527	27528**					

**Available in 60 days, specify colour with order

COLOUR ON REQUEST: CLASSIC COLOURS

White Helsinki	Beige Roma	Beige Tunisi	Brown Perugia	Yellow Parma	Orange Praga	Red Mosca	Pink Malaga	Bordeaux Bangkok
Green New Delhi	Green David	Green Melbourne	Green Ireland	Green Honolulu	Green Caracas	Violet Nizza	Violet Ankara	Light blue Palm beach
Blue Berna	Blue Baltico	Blue Anversa	Grey Oslo	Grey Amburgo	Grey Londra	Black Odessa		

COLOUR ON REQUEST: METALLIC COLOURS

Metal Ice	Metal Silver	Metal Bronze	Metal Gold	Metal Lemon	Metal Red	Metal Sky Blue	Metal Blue	Metal Dust Grey
--------------	-----------------	-----------------	---------------	----------------	--------------	-------------------	---------------	--------------------

GYNEX PROFESSIONAL with 3 engines



GYNEX PROFESSIONAL - GYNAECOLOGICAL CHAIR

Gynex professional chair offers functionality, aesthetics, ergonomics. It can satisfy all patient and user needs for a professional gynaecologic visit. Gynex professional is provided with versatile accessories: back paper-roll dispenser, metal basin, leg-supports. 3 engines allow, in addition to lifting, electric backrest and legrest movements. The combined movement of backrest and seat by means of self-balancing system (operated by a silent low tension actuator of 2,500 N) grants maximum comfort to the patient. The padded and pelted armrests allow a relaxing protecting sitting position, armrests are individually revolving and removable. Backrest is provided with removable height-adjustable headrest. The wide and soft backrest has independent adjustment movement through silent low tension actuator. This chair is moved by a silent low tension engine which allows a vertical excursion of 40 cm from zero position. The electrical drives are controlled by an ergonomic hand switch control, connected to the chair by a cable. Gynex professional can be easily transformed into an examination bed (horizontal position).



STANDARD ACCESSORIES

	GYNEX	GYNEX PROFESSIONAL
Polyurethane leg-rest	•	•
Padded leg-rest (colour coordinated)	•	•
Paper-roll holder	•	•
Basin	• (plastic)	• (metal)
Remote control	•	•

STOOL - with or without ring

For professional use, completely built in stainless steel, 360° pivoting self-braking castors, rubber-tyred.

Height can be comfortably adjusted by means of a gas pump; back-rest is 180° revolving and depth-adjustable.

Seat is balanced density padded and Peltex-lined.

Minimum/maximum height: 530/660 mm

