
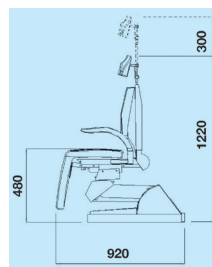


ENT CHAIRS

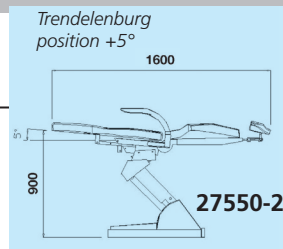

The seat is adjustable
(42 cm vertical excursion)

Width
650 mm

Backrest synchronized to the
leg-rest movement
can be adjusted at 90°
and inclined until 190°

Revolving and
removable
arm-rests in
light grey
expanded
polyurethane


Optional footrest
on request

Seat base in heavy metal, painted in
epoxydic powders, with support of pantograph
arm - 85° left and right rotation - lever locking device

- 27545 NEW OTO/PV ENT CHAIR with headrest - colour on request*
 - 27548 NEW OTO/PV ENT CHAIR with headrest - green Melbourne
 - 27549 NEW OTO/PV ENT CHAIR with headrest - light blue Parigi
- Patient chair with Peltex-lined seat and backrest.
Electric up and down movements by silent engine controlled through foot switch.
Back-rest can be adjusted up to -5° +90° by a gas spring.
Chair rotates 360° and is provided with a rotation locking device placed at the bottom.
Power: 230 V - 50 Hz (110 V - 60 Hz available on request).

- 27550 OTOPEX ENT CHAIR with headrest - colour on request*
 - 27551 OTOPEX ENT CHAIR with headrest - blue Chicago
 - 27552 OTOPEX ENT CHAIR with headrest - green Toronto
- ENT chair designed also for ambulatory operations.
Moved by very silent electrical engines operated by joystick foot-control. Also provided with button to re-set to previously selected working position. Upholstery made of simulated leather, high density undeformable foam padding (50 kg/cm³).
Seat and leg-rest are joined in a single piece.
Power: 230 V - 50 Hz (110 V - 60 Hz available on request).

ENT CABINETS

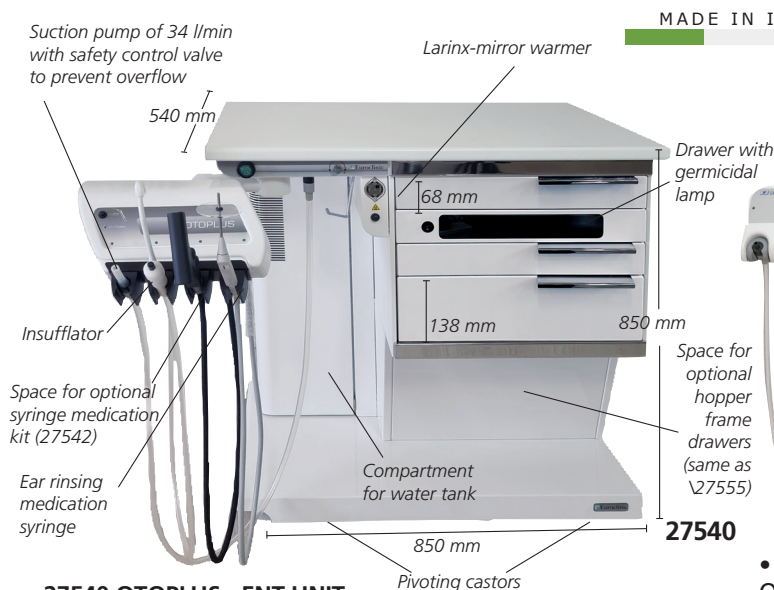

• 27540 OTOPLUS - ENT UNIT

The E.N.T. unit has rubber pivoting castors, a working surface made of scratch-proof rolled sections, 2 small and one large drawers provided of removable internal partitions, one drawer with inside germicidal lamps to keep instrument sterilization. Otoplus is provided with a silent suction pump of 34 l/min, with power supply complete with dim adjustment device for Clar/Storz lamp, larynx-mirror warmer, twin-cylinder dry compressor for tubal insufflation, 2-litres distilled water tank with acoustic signal for reserve level.
The articulated arm with instrument holding tray is complete with terminal suction cannula and terminal insufflation cannula, ear rinsing syringe with automatic water temperature control, (a microprocessor keeps temperature at a constant level of 37 °C).
Power 230 V - 50 Hz (Available 110 V - 60 Hz).

• 27539 OTOPLUS COMPLETE

Same as code 27540 with vasister drawers (code 27541) and endoscope support (code 27543) included.



• 27554 OTOPLUS DC ENT UNIT

Otoplus unit with column for diagnostic equipment (light source, camera or printer) and a support pole for monitor.

• 27555 OTOPLUS DC ENT UNIT - complete

Same as code 27554 with hopper frame drawers (code 27541) and endoscope support (code 27543) included.

• 27541 HOPPER FRAME DRAWERS - 2 independent drawers

• 27543 ENDOSCOPE SUPPORT - 2 transparent plastic tubes

1 tube 50 cm long and 1 tube 25 cm long both closed at the bottom.

• 27542 SYRINGE MEDICATION KIT including 1 syringe, 2 atomizers (1 for liquid and 1 for powder) and 3 ampoules in dark glass. To be ordered with ENT unit, should be installed at factory.


• 28239 FILTER - spare for 27539, 27540, 27554-5

• 27544 LED LIGHT SOURCE

Makes cold white light (5500 °K, 200 lumen) for endoscopy use.

It can be connected to endoscopes without need of F.O. cable. Handle connect directly to rigid endoscopes (Storz connection).

LCD display:
12.8x6.4 cm

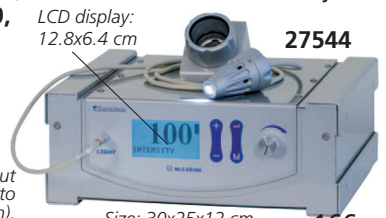

Size: 30x25x12 cm