

2 SECTION EXAMINATION COUCHES - CHROMED AND STAINLESS STEEL

- 44461 CHROME PLATED EXAMINATION COUCH
- 44463 STAINLESS STEEL EXAMINATION COUCH

Structure made of chrome plated steel (44461) or stainless steel (44463) tube Ø 25x25 mm.
 Backrest adjustable by rack.
 Padding and covering made of washable, flame retardant, plastic material. Rubber feet.
 Supplied with a paper roll holder, which can be placed either in the front or in the back part of the examination couch.
 Fireproof padding and covering, Class 1 IM.
 Dimension: 177x60xh 75 cm.
 Weight: 25kg
 Delivered disassembled.
In conformity with MDR.





KING HEIGHT ADJUSTABLE EXAMINATION COUCH

- 27613 KING HEIGHT ADJUSTABLE COUCH - blue - 9 - white frame
- 27612 KING HEIGHT ADJUSTABLE COUCH - any colour* - white frame
- 44477 KING HEIGHT ADJUSTABLE COUCH - with hole - mid light blue - black frame
- 44480 KING HEIGHT ADJUSTABLE COUCH - with hole - any colour* - black frame

*available in 30-40 days, see below available colours

Height adjustable (48-95 cm) examination couch by electric actuator (Linak engine 6.000 N) controlled by foot switch with head section easily adjustable (to +60° and -10°) by a gas pump. Strong and stylish leg structure in oval tubing 50x25x2 mm, central frame in rectangular tubing 50x20x1.5 mm painted with RAL 9003 white or black epoxy powder. Minimum thickness 70/80 micron, glossy finish. High density polyurethane foam platform with wooden base covered by leatherette upholstery. Both materials are fireproof Class 1 IM. Leatherette washable with the most common detergents.
 Size: 190x60xh 48/95 cm. Weight: 70 kg
 Maximum load: 150 kg.
 Operating power: 100-240 V, 50/60 Hz.
 Couch roll holder included.

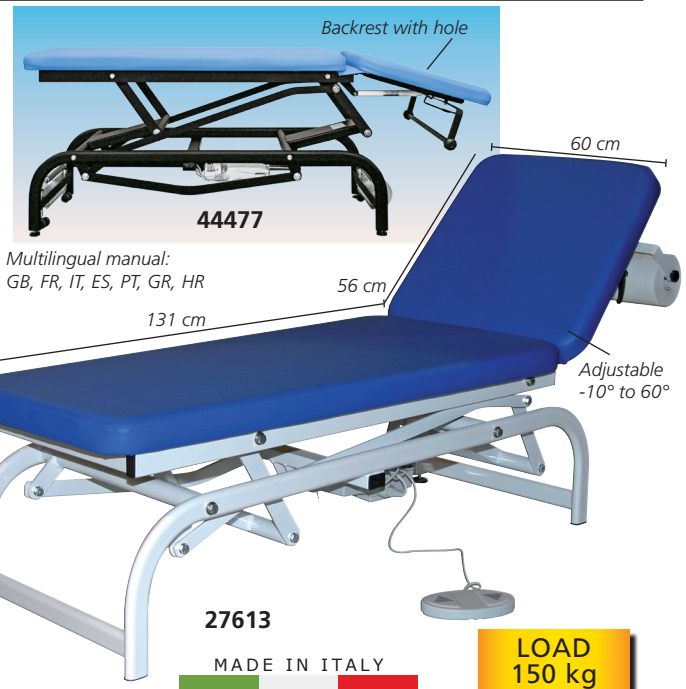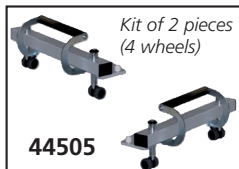

- 44505 RETRACTABLE/REMOVABLE WHEELS KIT - optional





4-SECTION VARIABLE HEIGHT TREATMENT COUCH

- 44531 4-SECTION VARIABLE HEIGHT TREATMENT COUCH - blue
- 44535 4-SECTION VARIABLE HEIGHT TREATMENT COUCH - any colour*

*available in 30-40 days, see below available colours

Couch made of high quality tubular that makes it very sturdy. The cushions can be reclined by a gas spring activated by an ergonomic lever under the top frame. Head section movable to +70° and to -35°. The two electrical columns allow to adjust height and reach TR/RTR position very quickly by pedals, push-button or two perimeter bars. High density, rubber padded, fire retardant and Class 1 IM top, covered in flame retardant Class 1 IM coated fabric (synthetic leather). Washable and disinfected with common detergents. Overall size: 195x68xh 48/89 cm. Weight: 80 kg. Complies with: EU 93/42 CEE. Maximum load: 200 kg. Supplied with mouth/nose hole plug, paper roll holder, 2 adjustable armrests and 4 adjustable feet.

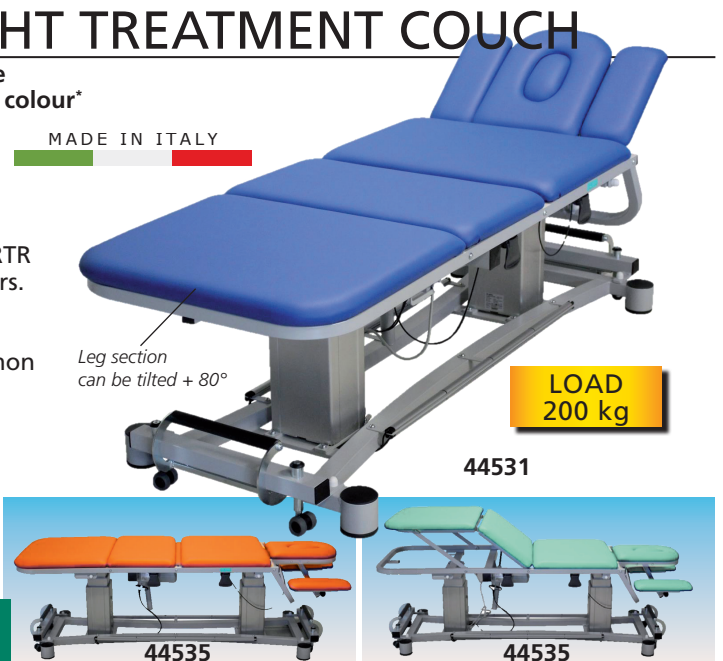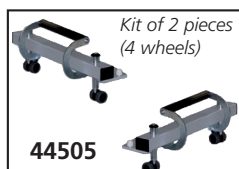

- 44505 RETRACTABLE/REMOVABLE WHEELS KIT - optional

