

AIR COMPRESSION THERAPY SOLUTIONS

Human body has a peripheral circulation of blood and lymph. Air compression therapy is an automatic drainage and circulation system which helps patients to recover from their disease by improving blood and lymph circulation and metabolism. It consists of air chambers and an air pump offering air compression massage to legs, arms or waist.


Medical indications (only 28441)

- DVT (Deep Vein Thrombosis)
- Lymph edema
- Post-paralytic
- Limb convulsion
- Phleboscrosis
- Chronic venous insufficiency (CVI)
- Post-mastectomy
- Post-trauma
- Spinal cord injury
- Varicose veins
- Sciatica
- Diabetes
- Sprains
- Arterial insufficiency

Home care indications (only 28441)


- Lymph edema
- Diseases caused by interrupted blood circulation
- Cellulitis fat dissolution
- Pregnant woman's limb treat.
- Numb feet or hands
- Cold feet or hands
- Settlement of stress
- Exercise effects
- Healthier body and elastic skin
- Old age, weakness and fatigue
- Insomnia

**MK400L
HELPS BLOOD
CIRCULATION**

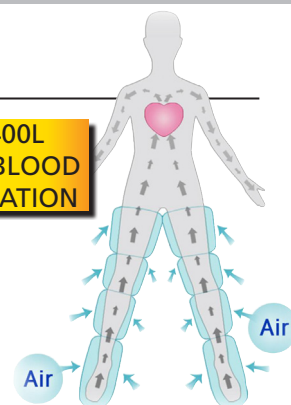

ANYCARE DIGITAL COMPRESSION MASSAGE

- **28415 ANYCARE DIGITAL COMPRESSION MASSAGE - LEG A**
 - **28416 ANYCARE DIGITAL COMPRESSION MASSAGE - LEG B**
- Pressure intensity and time can be adjusted according to user's needs. Individual chambers pressure testing. Portable and easy operation. For homecare use. Supplied with 2 legs A (74x62 cm) or B (79x75 cm) 4 chambers, 1 double hose and 2 plastic feet guides, multilingual manual (GB, FR, IT, ES).

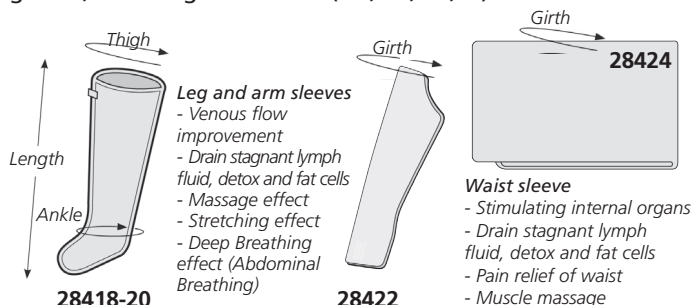

ACCESSORIES FOR ANYCARE AND MAYLEA						
New*	Old**	LEG SLEEVES	Patient height	Thigh	Ankle	Length
28418	-	Leg sleeve A	160-170 cm	62 cm	41 cm	74 cm
28419	-	Leg sleeve B	170-180 cm	75 cm	43 cm	79 cm
28420	-	Leg sleeve C	> 180 cm	76 cm	46 cm	94 cm
ARM - WAIST SLEEVES						
				Girth		Length
28422	28436	Arm sleeve		50 cm		88 cm
28424	28437	Waist sleeve		122 cm		38 cm
HOSES						
28426	-	Single hose to connect one sleeve				
28427	28440	Double hose to connect two sleeves				
28428	-	Distribution hose to connect three sleeves - needs 28426+28427				

*New sleeves can be used only with Anycare standard hose or codes 28426-8

**Old sleeves can be used only with Maylea standard hose or 28440

MK400L PROFESSIONAL COMPRESSION SYSTEM

- **28441 MK400L PROFESSIONAL COMPRESSION SYSTEM with 2 legs (XL) and 2 single hoses**

MK400L is a professional air compression therapy system with essential functions for different kinds of patients. This system provides 3 compression modes (wave, squeezing, double squeezing), reverse compression settings, skip functions to select directly part to be pressurized and compression speed controller. Touch key pad display. Manual: GB, IT, ES.

TECHNICAL SPECIFICATIONS

- Pressure range: 10-200 mmHg
- Use time: 1-90 minutes
- Automatic pressure detection
- Body size: 415x310x160 mm
- Operating voltage: 220-240 V, 50/60 Hz, 50 VA
- Body weight: 6.5 kg

You can use 2 cuffs at the same time, MK400 has 2 connection plugs



**PROFESSIONAL AIR
COMPRESSION
SYSTEM**

- **28442 LEG CUFF L - 6 CHAMBERS - optional**
- **28443 LEG CUFF XL - 6 CHAMBERS - spare**
- **28444 LEG EXTENSION ZIPPER XL - 10 cm extra extension for 28443**
- **28445 SINGLE HOSE - spare**
- **28446 ARM CUFF - 6 CHAMBERS - universal**
- **28447 WAIST CUFF - 6 CHAMBERS - min-max girth 94-125 cm**

A mode (Wave)

The selected air chamber inflates and deflates by one chamber sequentially from foot up to thigh.


B mode (Squeezing)

After all the selected air chambers inflate sequentially from foot up to thigh the inflated chambers deflate at once.

C mode (B mode two times + Reverse A mode 1 time)

The selected air chambers inflate and deflate by two chambers sequentially from foot up to thigh.

Reverse setting: it allows to change automatically the selected mode to reverse mode.

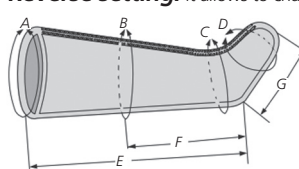

Cuff		Leg dimension (cm)						
Size		A	B	C	D	E	F	G
L		64,0	49,5	40,5	34,5	75,0	38,0	28,0
XL		73,5	55,0	42,0	35,5	85,0	43,0	31,0