

CHANGE-A-TIP™ CAUTERIES WITH HIGH TEMPERATURES

BOVIE CAUTERIES

- interchangeable sterile tips
- replaceable batteries
- reusable handles
- no sterilization necessary



GIMA code	INTERCHANGEABLE TIP CAUTERIES MINIMUM ORDER: 10 TIPS	Recommended for:
30590	Cautery handle "AA" (tip included)	All tips
30593	High temp loop tip - sterile	DER/ENT
30594	High temp 2" flex loop tip - sterile	GYN/ENT
30595	High temp 5" loop tip - long size - sterile	GYNAECOLOGY
30596	High temp fine tip - sterile	DER/CARDIO
30597	High temp elongated tip - sterile	VASECTOMY




REUSABLE HIGH TEMPERATURE CAUTERY WITH REPLACEABLE TIP

- CAUTERIES - 1,200°C (2,200°F) - made in Italy**
- interchangeable tips with nickel-chrome wire
 - replaceable batteries
 - reusable handles


Wide application in gynaecology, orthopedy, ear/nose/throat, dermatology, plastic and reconstruction surgery.

GIMA code	INTERCHANGEABLE TIP CAUTERIES MINIMUM ORDER: 10 TIPS
30598	Cautery handle "AA" (tip not included)
30599	High temp long fine tip - sterile - box of 10
30600	High temp large tip - sterile - box of 10
30601	High temp fine tip - sterile - box of 10

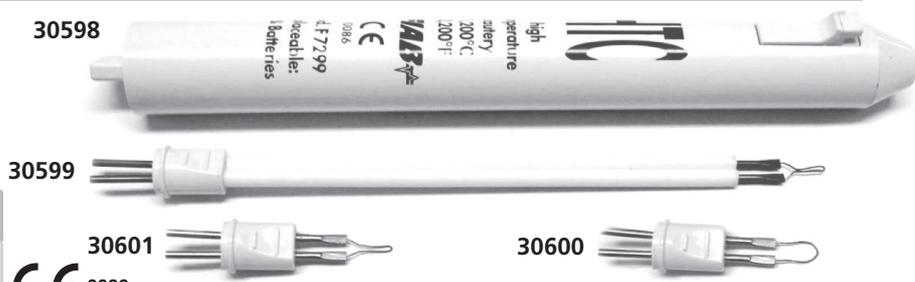

HIGH TEMPERATURE ELECTROCAUTERIES

- **30580 ELECTROCAUTERY - 1,200 °C - fine tip - made in Italy**
Electrocautery is a sterile bipolar battery operated disposable device for the cauterization of small vessels; the temperature of the tip is 1,200°C. With a fine tip required in most of the surgeries.
- **30582 GYNAECOLOGY CAUTERY - 1,200°C - fine tip**
Same as 30580 but with the tip fixed on a 12 cm long shaft.
- **30581 ELECTROCAUTERY - 600 °C - fine tip**
Same as 30580 but with lower temperature = 600°C. Especially designed for ophthalmology.
- **30583 EMERGENCY ELECTROCAUTERY - 800°C - Thick tip**
Same as 30580 but with "thick tip" for the perforation of nails in case of haematoma.

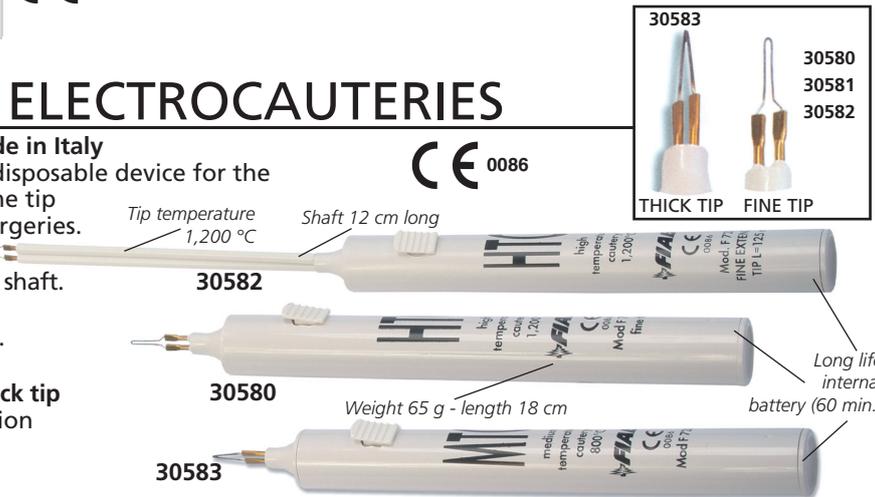

CARTS FOR ELECTROSURGICAL UNITS

Trolley designed to hold and carry electrosurgical units (including Argon units). Ergonomics has been particularly considered for the design, for example the height of 900mm with inclined top shelf allows easy reading of the unit's display. The trolley is provided with a second shelf with antiskid raised edges, two drawers placed at the base and a rear space with securing straps to hold two Argon cylinders. The upright stands holding the shelves are designed to take possible optional accessories. All parts of the trolley are made of epoxy-painted, tempered steel and aluminum. Available version with 4 socket group, fuses box, independent grounding plug, power cable 3m.

TECHNICAL SPECIFICATIONS		
Load capacity:	- Total load 200 kg	- Inclined shelf 50 kg
Inclined shelf:	42x45 cm (27864/6)	50x50 cm (27865/7)
Drawers and shelf:	42x31 cm (27864/6)	50x36 cm (27865/7)

GIMA code	CARTS FOR ELECTROSURGICAL UNITS
27864	DIATERMO CART - medium
27865	DIATERMO CART - large
27866	DIATERMO CART with power - medium
27867	DIATERMO CART with power - large

