



"KARMAN" AND "BERKELEY" INTRAUTERINE SUCTION CURETTES

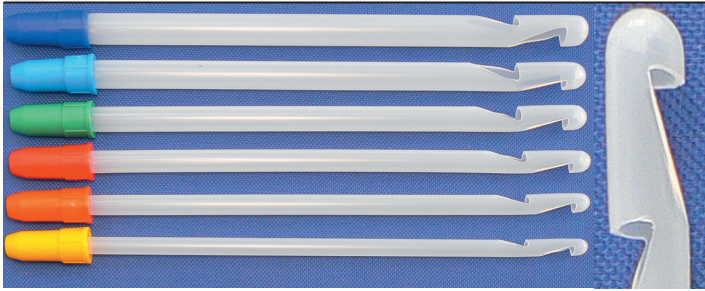

Length 24 cm with connector, 20.5 cm without

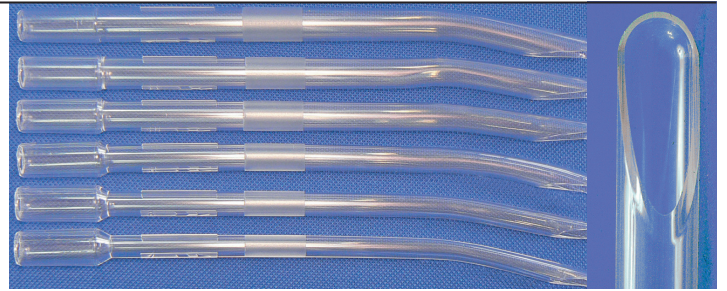

Length 20 cm with connector, 17 cm without


GIMA code	"KARMAN" STYLE CURETTES FLEXIBLE - SINGLE USE - STERILE	Minimum order
29858	Karman Curette N° 4 - 4 mm x 24 cm	box of 20
29859	Karman Curette N° 5 - 5 mm x 24 cm	box of 20
29860	Karman Curette N° 6 - 6 mm x 24 cm	box of 20
29861	Karman Curette N° 7 - 7 mm x 24 cm	box of 20
29862	Karman Curette N° 8 - 8 mm x 24 cm	box of 20
29863	Karman Curette N° 9 - 9 mm x 24 cm	box of 20
29864	Karman Curette N° 10 - 10 mm x 24 cm	box of 20
29865	Karman Curette N° 11 - 11 mm x 24 cm	box of 20
29866	Karman Curette N° 12 - 12 mm x 24 cm	box of 20

MADE IN ITALY



GIMA code	"BERKELEY" STYLE CURETTES CURVED - RIGID - SINGLE USE - STERILE	Minimum order
29870	Berkeley Curette N° 7 - Ø 7 mm x 20 cm	box of 20
29871	Berkeley Curette N° 8 - Ø 8 mm x 20 cm	box of 20
29872	Berkeley Curette N° 9 - Ø 9 mm x 20 cm	box of 20
29873	Berkeley Curette N° 10 - Ø 10 mm x 20 cm	box of 20
29874	Berkeley Curette N° 11 - Ø 11 mm x 20 cm	box of 20
29875	Berkeley Curette N° 12 - Ø 12 mm x 20 cm	box of 20

CANNULAS FOR ENDOMETRIAL SAMPLING

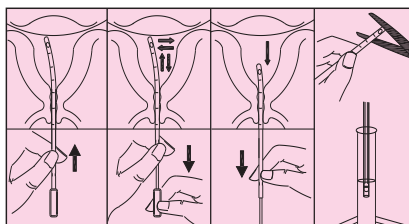

Rounded tip for easy and smooth insertion with 1 or 4 holes Endometrial sampler (mucus)

Graduation in cm for easy measurement of uterus length

Endometrial sampler for cytologic examination



29826

Asymmetrical V shape blades tip with protrusions for sampling



29826

Hysterometry can be measured by means of the reference points located at 3, 6 and 9 cm along the sheath

29940

- 29826 PIPELLA 1 hole - disposable, sterile - box of 100
- 29827 PIPELLA 4 holes - disposable, sterile - box of 100

For endometrial cytologic sampling.

It allows a painless and easy routine screening of the endometrium avoiding dilation and traction of the uterus neck.

Hole: Ø 2.0-2.5 mm

External diameter: 3.1 mm


Length of cannula: 23.5 cm

Length of handle: 26.5 cm

- 29940 ENDOGYN - disposable, sterile - box of 25

For cytologic endometrial sampling. It allows a painless and easy routine screening of the endometrium avoiding dilation and traction of the uterus neck.

An outer sheath covers a mandrel with asymmetrical "V" shaped blades on the top. Each arm features fine toothings and thin cutting edges which are able to remove and keep more endometrial tissue. The head of one of the two arms features a hemispherical head which, when it goes back inside the sheath, allows for previously removed tissue to be aspirated.

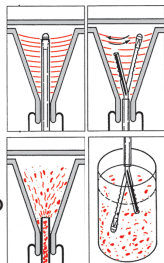

29940

MICROSCOPE SLIDES AND BOX

29775



29777

Twin frosted in both ends, ground edges, sandblasted band

29775

Single end frosted, sandblasted band, ground edges

29774

Plain side, ground edges

29776

BOX FOR SLIDES
Plastic, transparent box for storage of up to three slides.

Multilingual box:
GB, FR, IT, ES, PT, GR

MICROSCOPE SLIDES

Size 76x26 mm, 1.0-1.2 mm thick 45° corners

GIMA code	MICROSCOPE SLIDES box of 200 pcs. (4 boxes of 50 pcs.)
29774	Plain side, ground edges
29775	Single end frosted, sandblasted band, ground edges
29777	Twin frosted in both ends, ground edges, sandblasted band
BOX FOR SLIDES	
29776	Box for slides - box of 20 pcs

CITOLOGY FIXATION

- 29788 GIMAFIX SPRAY FOR CITOLOGY FIXATION - 100 ml - box of 12
- 29789 GIMAFIX SPRAY FOR CITOLOGY - 200 ml - box of 24

Fixative for cytology. Composition: isopropyl alcohol, polyethylene glycol.

Spray the product at a distance of 20 cm from the slide.

It dries quickly and forms a thin film to maintain over time the slide.

It complies with 2017/746 EU regulation.

Multilingual instructions:
GB, FR, IT, ES.



MADE IN ITALY

29788

29789

220