
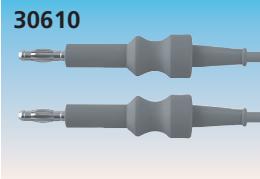


MONOPOLAR FORCEPS, SCISSORS AND CABLES

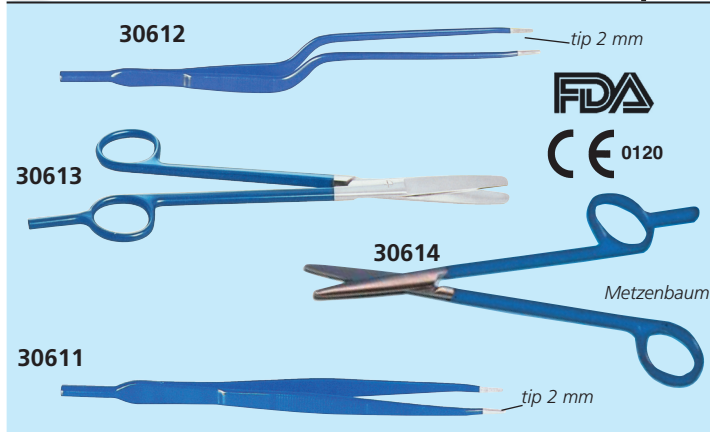

Monopolar forceps and scissors are provided with a 4 mm connection for standard monopolar cables. These can be used with all Gima, LED, Erbe and many other electrosurgical units in the market.

GIMA code	INSULATED MONOPOLAR INSTRUMENTS
30610	Monopolar cable - 4 mm pin - 3 m
30611	Mc INDOE straight monopolar forceps 15 cm (6")
30612	Bayonet monopolar forceps - 20 cm (7 1/2")
30613	Monopolar scissors - straight - 18 cm (7")
30614	Metzenbaum monopolar scissors - straight - 18 cm (7")
30615	Metzenbaum monopolar scissors - curved - 18 cm (7")

3M ELECTROSURGICAL GROUNDING PADS WITH SAFETY RING

- 30496 3M 9160 SPLIT GROUNDING PAD - box of 100
 - 30497 3M 9160F SPLIT GROUNDING PAD - box of 200 (40 pouches of 5 pieces)
 - 30498 3M 9165 SPLIT GROUNDING PAD - precord - box of 40
- The 3M Universal Grounding Pads all have a conductor area of 15 square inches and come with the patented 3M Green Safety Ring technology which allows them to perform the same as grounding pads that are up to 33% larger in conductor area. Green Ring allows adult and pediatric use.

They have fluid resistance, non-woven backing which allows heat to pass through up to 25% faster than foam, and they have no minimum age or weight limit.



30496, 30497
(3M 9160)



30498
(3M 9165)


HANDLES FOR ELECTROLYSIS NEEDLES


Needle handles for all GIMA Diatermo and electro-depilators

- 30552 HANDLE for electrolysis needles
- 30415 PROFESSIONAL HANDLE for electrolysis needles




BALLET STERILE DISPOSABLE ELECTROLYSIS NEEDLES






30400-30401



30405-30406



30410-30411

Stainless steel	Gold plated	Insulated	ELECTROLYSIS NEEDLES - STERILE	Minimum order
30400	30405	30410	K2 needle - 0.055 mm (0.002")	box of 50
30401	30406	30411	K3 needle - 0.075 mm (0.003")	box of 50

Ballet needles

Manufactured with a precise point and a mirror-like surface polish. The rounded tapered point of the needle assists the electrologist in locating the follicle opening and upon insertion the needle will slide smoothly and accurately to the base of the follicle.

Stainless steel needles

Ideal for patients with strong healthy skin who respond well to electrolysis treatment.

Insulated needles

Suggested for patients with high sensitivity to discomforts of electrolysis.

24 carat gold plated needles

For clients who suffer from allergies and sensitive skin.



ELECTRODEPILATOR

• 28340 ELECTRO-DEPILATOR 400

This device ensures a definitive depilation with thermolytic method. It delivers a current of radiofrequency to cause the heating and consequent coagulation of the papilla, reached through the introduction of a needle into the follicle of the hair. Provided with highly reliable electronic integrated components, with RF intensity adjustable with continuity and luminous working signal. It comes complete with footswitch, handpiece for depilation, depilation needle and power cable (2 m).

TECHNICAL SPECIFICATIONS

RF output frequency: 1.9 ±10% MHz
Modulation: sinusoidal
Maximum output power: 7 ±2Watt
Size: 260x265x110 mm Weight: 4 kg
Operating voltage: 115/230 V - 50/60 Hz switchable
Norms: CEI 62-39 Class: I-BF

- 30552 HANDPIECE FOR DEPILATION
- 30577 FOOTSWITCH

