
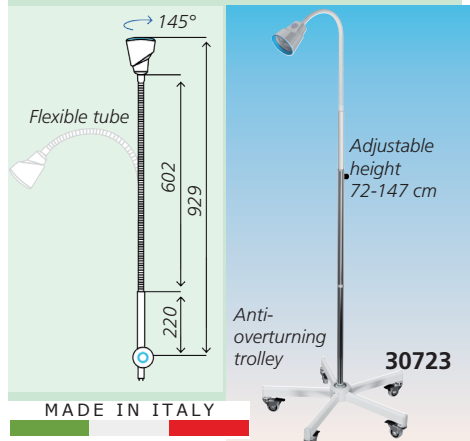


# HALOGEN AND LED LIGHTS FOR EXAMINATION AND MINOR SURGERY

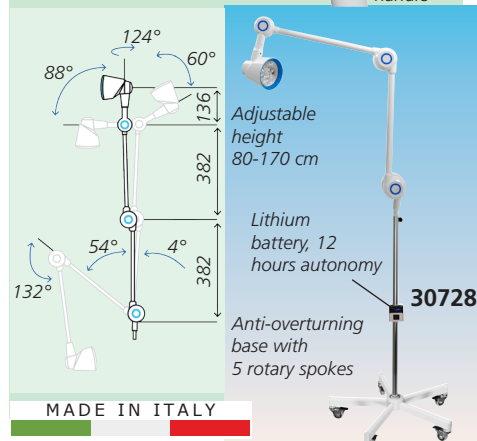

## ALFA-FIX AND ALFA-FLEX

Three LED light sources with coinciding lenses and polycarbonate protection shield. Thanks to three independent light sources, shadows can be reduced and provide a deep cylindrical light with highly reduced heat irradiation. Each LED integrates a resistance to ensure the continuous operation of the lamp even in the rare case of a LED fault.

### Design

New design makes the lamp extremely easy to use and manoeuvre, with an additional rotation of the reflector.

### Touch control

The light intensity can be adjusted through the touch button, a smooth touch of the control allows to regulate the light intensity from 4% to 100%.

### Interchangeable plugs

The transformer, integrated in the plug for a lighter and smooth movement



of the arm structure, is supplied with interchangeable plugs (UK, US, EU).

### Flexible arm (30723-5)

The flexible arm ensures easier light positioning and is covered with a smooth white shrink-wrap sheath for easier cleaning and disinfection.

### Battery kit

Included with 30728, available in 45 days for 30723 (to be purchased with light).

QR manual: GB, FR, IT, ES, PT.

GIMA code	ALFA-FIX AND ALFA-FLEX
30723	Alfa-flex LED - trolley
30724	Alfa-flex LED - wall
30725	Alfa-flex LED - table with clamp
30726	Alfa-fix LED - trolley
30727	Alfa-fix LED - wall
30728	Alfa-fix LED - trolley with battery

ACCESSORIES AND SPARE	
30766	Spare bulb Simplex - IRC 35 W for 30760-30763
30768	Spare trolley, base Ø 61 cm for Simplex and Alfaflex
30769	Spare trolley, base Ø 61 cm for Solesud2 and Alfafix
30775	Rail connection with clamp
30772	Spare glass Ø 8.5 cm for 30723-30727
30770	Spare glass Ø 12.3 cm for Solesud

Model	ALFA-FIX	ALFA-FLEX	LUXIFLEX	LED	LED Plus
Type of lighting	LED	LED	halogen	LED	LED Plus
Consumption (Watt)	8	8	35	6.5	10
Life of the bulbs (hours)	> 60,000	> 60,000	5,000	50,000	50,000
Light intensity at 50 cm (lux)	60,000	60,000	50,000	15,000	35,000
Light intensity at 1 m (lux)	20,000	20,000	-	-	-
Colour temperature ± 5% at 1 m (K)	4,500	4,500	4,000	4,500	4,500
Colour rendering index (Ra)%	94	94	95	80	97
Light field size at 50 cm (Ø mm)	150	150	80	80	80
External Ø of the reflector (mm)	95	95	80	80	80
Colour	white	white	white	white	white
Norms	IEC 60601-1	IEC 60601-1	EN 60601-1-2		
Marks	EU 2017/745	EU 2017/745	CE 93/42		
Weight (kg)	5	5			
Operating voltage V	100-240	100-240	230	110/230	
Operating Hz	50/60	50/60	50/60		

- 30754 LUXIFLEX on trolley - halogen
  - 30756 LUXIFLEX on trolley - LED
  - 30755 LUXIFLEX PLUS on trolley - LED
- Stylish and versatile light for examination and minor surgery procedures in health-centres, clinics and hospitals. Suitable also for veterinary examinations and procedures. Class II.
- Safe low voltage operation without infrared radiation. Electronic transformer lodged at the bottom of the lamp. Flexible arm and 5 castors trolley allow easy use and movability.
- 30757 SPARE BULB LUXIFLEX - halogen

