

VISIO 5 HQ WALL DIAGNOSTIC STATIONS 3.5 V - XENON OR LED

VISIO 5 HQ WALL DIAGNOSTIC STATIONS 3.5 V

Modern design diagnostic stations to satisfy all practitioners and clinics needs. Two configurations available:

- oto-ophthalmoscope with disposable ear tip dispenser
- oto-ophthalmoscope with ear tips dispenser, thermometer, B.P.M.

Two ophthalmalmo heads with 7 or 21 apertures
Two illumination technologies: LED or Xenon-Halogen.
Regulated output prolongs the bulb life and guarantees optimal light intensity

- ergonomic designed non-slip handles
- constant output transformer with integrated surge protection
- soft start technology to introduce voltage gently, eliminates cold-shock when turning ON, increasing bulb life
- adjustable light intensity, low power consumption
- automatic ON-OFF switch once the handle is removed or replaced in the cradle
- precise brightness control, between 3 and 100%
- spiral cord of 3 meters

Supplied with multilingual manual: GB, FR, IT, ES.



AVAILABLE WITH 5XE OR 5GE OPHTHALMO HEAD

7 APERTURES

- Small light ●
- Medium light ●
- Large light ●
- Grid ●
- Slit |
- Half light ▸
- Red free light* ●
- Blue light** ●

5XE: 21 APERTURES

5GE: 7 APERTURES

3 filters

31622-3

31620-1
31624-5

TECHNICAL SPECIFICATIONS OPHTHALMOSCOPES

	31620-1	31624-5	31622-3
Apertures	7	21**	21***
Focusing lenses	28	68	
Diopters range	-25 to +40	-30 to +38	

Visio 5XE HQ additional features (31622-3)
Coaxial optics for easier examination of undilated pupils
Detection of corneal abrasions with cobalt blue filter
Polarizing filter to virtually eliminate corneal reflection

3 FILTERS**

- Polarizing filter ●
- Neutral light ○
- Red free filter ●

**5XE (31622-3 only)

*** 7 apertures / 3 filters combination

VISIO 2000 F.O., XENON OTO, OPHTHALMO, ENT SETS - 3.5 V


Professional diagnostic sets with fiber optic illumination, Xenon light 3.5 V, rheostat to adjust the intensity of light. Made of resistant re-inforced plastic and durable easy locking system. Supplied in deluxe case, should be used with charging station 31542 and Li-ion rechargeable battery 31541. Only plug-in set (31537) has direct AC connection and included battery (31547).

- 31535 VISIO 2000 F.O. XENON OTOSCOPE - 3.5 V
- 31536 VISIO 2000 F.O. XENON OPHTHALMOSCOPE - 3.5 V
- 31537 VISIO 2000 F.O. XENON OTO-OPHTHALMOSCOPE PLUG-IN SET - 3.5 V - 2 pin plug
- 31538 VISIO 2000 F.O. XENON OTO-OPHTHALMOSCOPE - 3.5 V
- 31539 VISIO 2000 F.O. XENON DIAGNOSTIC SET - 3.5 V
- 31428 OPHTHALMOSCOPE BULB - spare
- 31544 OTOSCOPE BULB - spare
- 31541 LI-ION RECHARGEABLE BATTERY 3.5 V
- 31547 LI-ION RECHARGEABLE BATTERY 3.5 V for PLUG-IN HANDLE SET 31537- spare
- 31542 ZOOM CHARGING STATION

Supplied with EU 2 pin, UK and US plug. See page 377

COMPATIBLE MINI EAR SPECULA: 31485-8, 31491-2, 31495-9 page 266

Swivel type window with 3X magnification for sharp images and minimal reflection

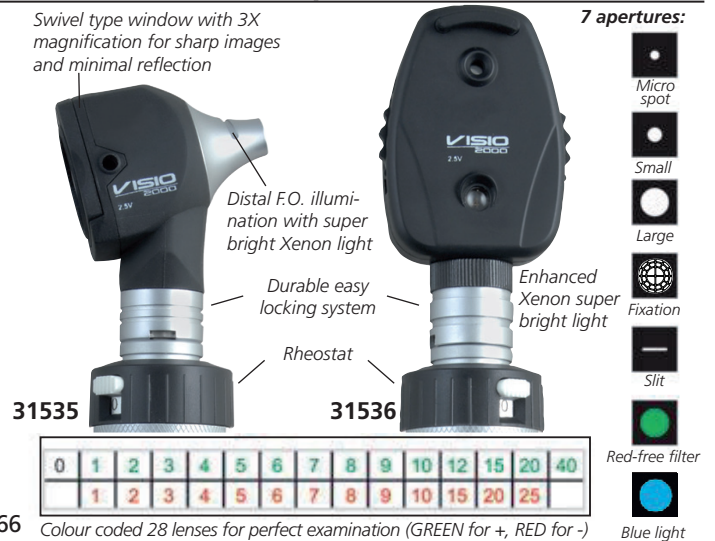

7 apertures:

- Micro spot ●
- Small ●
- Large ●
- Fixation ●
- Slit ●
- Red-free filter ●
- Blue light ●



F.O. Xenon otoscope set 3.5 V 31535

F.O. otoscope head, battery handle with bayonet lock, 1 set of 5 reusable eartips (2.5, 3.0, 4.0, 5.0 & 10.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare bulb, in a rigid case.

F.O. Xenon ophthalmoscope set 3.5 V

31536 F.O. ophthalmalmo head, battery handle with bayonet lock, spare bulb, in a rigid case.



F.O. Xenon oto/ophthalmo plug-in set 3.5 V 31537

F.O. otoscope head, F.O. ophthalmalmo head, battery handle with bayonet lock and plug-in, 1 set of 5 reusable eartips (2.5, 3.0, 4.0, 5.0 & 10.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto and ophthalmalmo bulbs, rechargeable Li-ion battery 31547, in a rigid case.



F.O. Xenon oto/ophthalmo set 3.5 V 31538

F.O. otoscope head, F.O. ophthalmalmo head, 2 battery handles with bayonet lock, nasal illuminator with F.O. light, 1 set of 5 reusable eartips (2.5, 3.0, 4.0, 5.0 & 10.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto bulb, suitable also for nasal illuminator, and ophthalmalmo bulb, in a rigid case.



F.O. Xenon ENT diagnostic set 3.5 V 31539

F.O. otoscope head, F.O. ophthalmalmo head, battery handle with bayonet lock, nasal illuminator with F.O. light, curved illuminator with bulb, plain mirror 20.0 mm, plastic tongue depressor, 1 set of 5 reusable eartips (2.5, 3.0, 4.0, 5.0 & 10.0 mm) tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto and ophthalmalmo bulbs, in a rigid case.