


ERKA FINESSE STETHOSCOPES

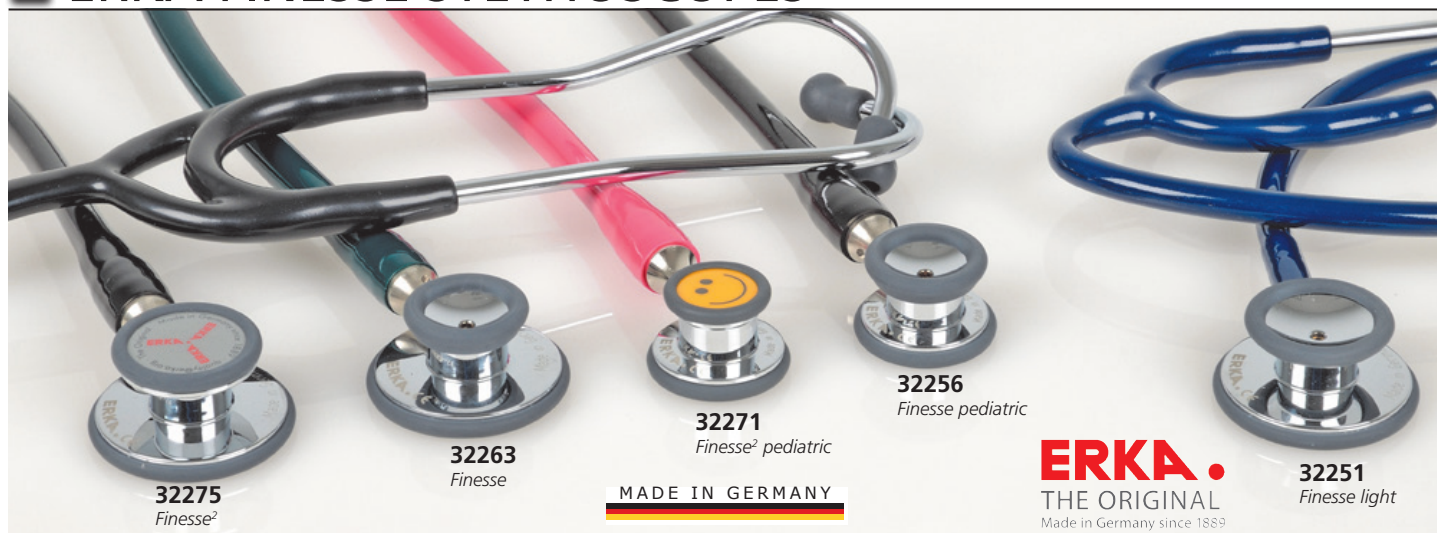


FINESSE STETHOSCOPES

With the reversible double chest-piece you'll be prepared for all eventualities in everyday medical practice. The membrane side guarantees optimal auscultation for the entire frequency spectrum and the open bell side is specifically designed for low frequencies.

FINESSE² STETHOSCOPES

The reversible, combined chest-piece has two membrane sides, each making the full frequency range perfectly audible. The large membrane is suitable for auscultation of adults and the small one for children or slim patients.

All Erka stethoscopes are supplied with carrying bag

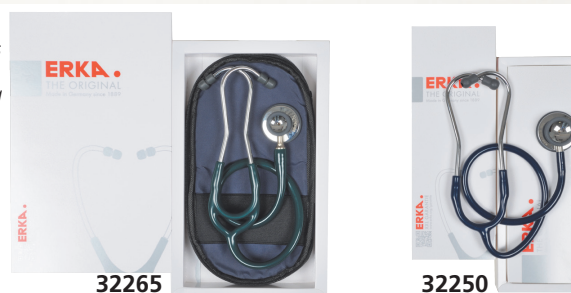

FEATURES	Finesse Light	Finesse Finesse²
Pre-tensioned ERKA. convex membrane provides excellent resonance and sound.	X	X
Non-chill ring ensures highest patient comfort.		X
Superior polish, chromium-plated, non-porous surface and brass chest-piece ensure perfect surface hygiene.	X	X
Two internal channels run parallel in the stethoscope tubing preventing any noise caused by friction of two separate tubes.		X
Through the separately integrated spring in the tube, there is no acoustic loss from the chest-piece to the ear.		X
Ergonomic design and super soft tips provide excellent acoustic seal against external noise interference and maximize user's comfort.	X	X
Standard 15° inclined binaurals can be individually adapted to the specific ear canal of the user.	X	X

GIMA code	ERKA code	ERKA STETHOSCOPES - tubing length - 55 cm
FINESSE LIGHT ADULT - head Ø 44.5/33.5 mm		
32250	520 000 10	■ Finesse light - black
32251	520 000 20	■ Finesse light - navy blue
32252	520 000 55	■ Finesse light - dark green
FINESSE PEDIATRIC - head Ø 31/26.5 mm		
32256	549 000 00	■ Finesse pediatric - black
32257	549 000 25	■ Finesse pediatric - light blue
32258	549 000 35	■ Finesse pediatric - pink
FINESSE ADULT - head Ø 44.5/33.5 mm		
32262	550 000 00	■ Finesse adult - black
32263	550 000 05	■ Finesse adult - dark grey
32264	550 000 20	■ Finesse adult - navy blue
32265	550 000 55	■ Finesse adult - dark green
32266	550 000 60	■ Finesse adult - burgundy
FINESSE² PEDIATRIC - head Ø 30/21.3 mm		
32270	536 000 00	■ Finesse² pediatric - black
32271	536 000 35	■ Finesse² pediatric - pink
FINESSE² ADULT - head Ø 47/30 mm		
32275	535 000 00	■ Finesse² adult - black
32276	535 000 05	■ Finesse² adult - dark grey
32277	535 000 20	■ Finesse² adult - navy blue

ERKAPHON STETHOSCOPES



Stethoscope with double chest piece is switchable from membrane to bell function. The membrane side is suitable for the auscultation of high and low frequency sounds while the bell side is designed for cardiac impulses.

ERKAPHON STETHOSCOPES:

High quality precision stethoscopes for excellent acoustics. All stethoscopes are equipped with 56 cm long Y tube, brass chrome-plated binaural and chest-piece (32350, 32358) or aluminium (32344-6), membrane Ø 44 mm, bell Ø 30 mm (32358 only).

GIMA code	ERKA code	ERKAPHON STETHOSCOPES
ERKAPHON FLAT CHESTPIECE - aluminium		
32344	544 000 10	Erkaphon flat chestpiece - black
32345	544 000 20	Erkaphon flat chestpiece - navy blue
32346	544 000 25	Erkaphon flat chestpiece - light blue
ERKAPHON FLAT CHESTPIECE - chromed		
32350	543 000 10	Erkaphon flat chestpiece - black
ERKAPHON DOUBLE CHESTPIECE - chromed		
32358	541 000 10	Erkaphon double chestpiece - black