

CHECKME™ LITE HEALTH MONITOR WITH BLUETOOTH

• 35185 CHECKME™ LITE HEALTH MONITOR with bluetooth

Innovative all-in-one device with multiple functionalities into one compact and portable unit. It allows to measure simultaneously ECG, heart rate, QRS duration, SpO₂, pulse rate and pulse index. With colour touch screen, well organized user interface and audio instructions, is very intuitive and easy to use. Clinical approved accuracy.

- supports both cable and cable-free measurement of ECG and SpO₂ (optional cable)
- checks, captures and identifies abnormal ECG in 20 sec.
- professional ECG algorithm provides heart rate, QRS duration, ST segment and rhythm analysis
- check blood oxygenation with built-in SpO₂ sensor in 20 sec
- high quality ECG waveform can be saved and replayed
- professional PC data management software for review and generating reports (compatible with windows OS) MyBeats Hub App.

• 35190 EXTERNAL SpO₂ ADULT SENSOR - reusable - optional

**DAILY CHECK
JUST TAKES
20 SECONDS**

CE 0197


**ECG + HR + SpO₂ +
QRS + Pulse Rate + PI**

Manual: GB, IT, ES, DE, FR
Audio instructions:
GB, PL, DE, IT, FR, ES, HU
PC software: GB



Download on the
App Store

GET IT ON
Google play

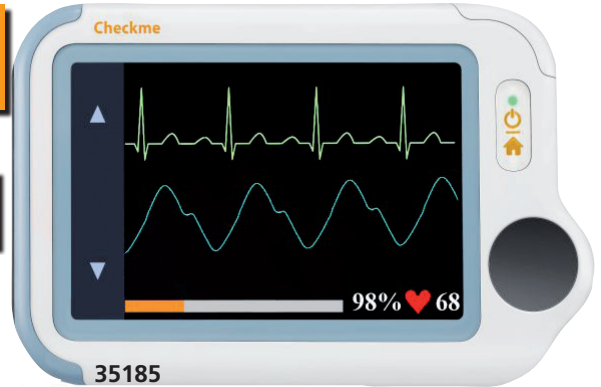

Size: 88x56x13 mm Weight: 64 g

STANDARD ACCESSORIES

PC software: GB Neck strap
ECG cable Charging cable, micro D
ECG electrodes

CHECKME™ PRO VITAL SIGNS MONITORS - WIRELESS WITH BLUETOOTH

• 35186 CHECKME™ PRO VITAL SIGNS MONITOR - with Bluetooth • 35187 CHECKME™ PRO VITAL SIGNS MONITOR WITH ECG HOLTER - with Bluetooth

All-in-one, compact, versatile and portable wireless vital signs monitors for checking your health anytime anywhere. Wide range of application scenarios, including: outpatient, emergency, surgeries, telemedicine. MyBeats Hub App.

• 35188 AIR BP BLOOD PRESSURE UNIT* - optional • 35189 AIR BP BLOOD PRESSURE UNIT* WITH DISPLAY - optional

Ultra-portable upper arm blood pressure device. See below.

*Includes main unit, cuff and charge cable

• 35190 EXTERNAL SpO₂ ADULT SENSOR - reusable - optional



Download on the
App Store

GET IT ON
Google play

Size: 88x56x13 mm
Weight: 80 g

**ECG + SpO₂ +
BP (needs 35188-9)
+ Sleep Apnea +
Temp + Pedometer
+ 24 hours ECG
holter (only 35187)**

CHECKME™ PRO FEATURES

ECG

- supports both cable and cable-free operation
- 500 Hz sampling rate and wide bandwidth
- 30/60s waveform can be replayed
- professional ECG algorithm to provide Heart Rate, QRS duration, ST segment, QT/QTc interval, and rhythm analysis ECG holter (only 35187)
- easy to wear and use
- continuously record ECG for 24 hours
- check report and ECG waveform on APP or PC
- the report includes heart rate, rhythm and HRV (heart rate variability)
- holter report can be stored in Cloud, printed or saved as PDF

SpO₂

- supports both cable and cable-free operation
- measurement of oxygen saturation, pulse rate and pulse strength

Sleep apnea screening

- SpO₂ monitor for OSA (Obstructive Sleep Apnea)
- records and tracks SpO₂ and Pulse Rate up to 10 hours
- mobile APP to check SpO₂ change in details

Blood Pressure via Bluetooth (needs code 35188 or 35189)

- manage all vital signs data together - see below

Temperature - Integrated infrared temperature sensor

- quick measurement with high accuracy
- safe and comfortable, no need to wake up child while sleeping

Pedometer

- accurately calculate steps, calories and amount of fat burned
- built-in high accuracy 3D smart sensor

APP - support iOS and Android

- built-in Bluetooth, all data can be exported from device to phone/tablet

PC browser


- medical data management software
- user data, user profile, readings and wave form can be exported from Checkme to the browser by USB

Manual: GB, ES, IT, DE, FR

Audio instructions: GB, ES, IT, DE, FR, PL, RU, HU

Pc software: GB

App: GB, DE, FR, ES

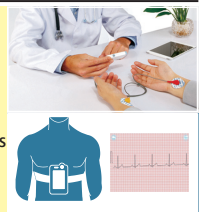

STANDARD ACCESSORIES

ECG cable, SpO₂ sensor, ECG electrodes, charge cable, neck strip, user manual, protective cover, wrist band, Y adaptor and sport chest band (only 35187)

AIR BP ULTRA-PORTABLE SMART BLOOD PRESSURE MONITOR

• 35188 AIR BP BLOOD PRESSURE UNIT

• 35189 AIR BP2 BLOOD PRESSURE UNIT WITH DISPLAY

Ultra-portable smart upper arm blood pressure device.

- one battery charge for 1000 measurements

- easy to use, robust and durable

- professional measurement

- readings:

- systolic, Diastolic, Pulse Rate;
- pulse Pressure, MAP (Mean Arterial Pressure)
- irregular heartbeat

- work with Smartphone

- real-time Virtual Mercury Column in App

- during measurement

- real-time Pulse Waveform



Smartphone
not included



35189

- multi-patient Management

- unlimited recording

- works with Checkme Pro 35186, 35187.

- Supplied with cuff, charge cable and

- multilingual manual (GB, FR, IT, ES,

- DE).



Download on the
App Store

GET IT ON
Google play



TECHNICAL SPECIFICATIONS

Cuff:

Wireless:

Charge input:

Battery life:

Pressure measurement range:

Pulse rate:

Mobile App for iOS:

Mobile App for Android:

Size - weight:

Clinical accuracy:

Upper arm circumference: 22-42cm

Bluetooth 4.0 BLE

Micro USB, DC 5V

1000 measurements, rechargeable

0-300 mmHg ±3 mmHg

40 to 200 bpm ±2 bpm

iOS 9.0 or above, iPhone 4s/iPad 3 or above

Android 4.4 or above, with Bluetooth 4.0 BLE

68x20 mm (main unit) - 30g (main unit)

meet IEC 80601-2-30