


DALLAS SPHYGMO: WALL AND TROLLEY



- **32749 DALLAS SPHYGMOMANOMETER - wall model**
Provided with rear plastic basket for cuff storage.
- **32747 DALLAS SPHYGMOMANOMETER - trolley model**
Provided with rear plastic basket and stainless steel trolley adjustable in height, with 5 plastic castors \varnothing 5 cm. All DALLAS sphygmomanometers are provided with extra large dial plate (\varnothing 150 mm) for easy reading, black nylon calibrated cuff and latex quality inflation system and bladder. Scale 0-300 mmHg. Multilingual manual: GB, FR, IT, ES, PT, DE, GR. Multilingual box (32749 only): GB, FR, IT, ES, PT, DE, GR, Arabic.

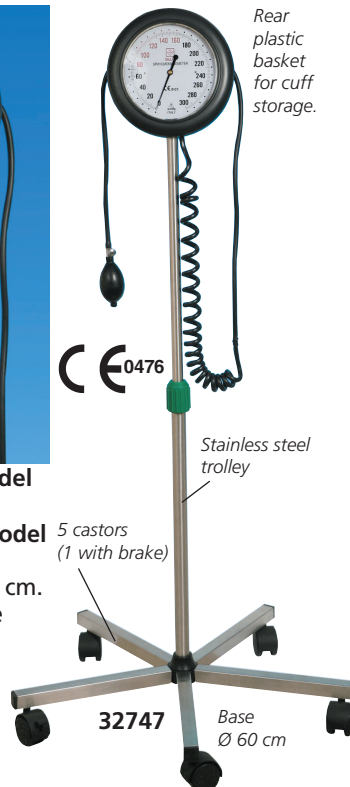

DENVER SPHYGMO

- **32744 DENVER SPHYGMOMANOMETER - trolley - adult cuff**
Trolley sphygmomanometer, light but with high stability and 360° swivelling castors. Adjustable height (90 to 135 cm) with metal basket. Square-shaped-dial 145x145 mm. Nylon cuff, 2-tube adult TPU bladder, PVC bulb, dual use air-release valve with spring, spiral PVC tube (1.8 m). Multilingual manual: GB, FR, IT, ES, PT, DE, GR.
- CALIBRATED NYLON CUFFS**
- **32839 CHILD CUFF - black**
- **32850 ADULT CUFF - black - spare**
- **32841 ADULT LARGE CUFF - black**
All cuffs with 2 tubes TPU bladder.



LCD DISPLAY, MERCURY-FREE SPHYGMOMANOMETERS

MERCURY FREE SPHYGMOMANOMETERS

LCD display mercury-free sphygmomanometers provide reliable and precise results.

Manual inflation and air release. Two measuring methods: button AUSC (Auscultation), for use with a stethoscope, button OSCI (Oscillometric), to get the result automatically, without need of a stethoscope, just inflate the cuff not less than 180 mmHg. Suitable for pediatric use with auscultation method only (need 32825). Supplied with adult cuff. Multilingual box and manual: GB, FR, IT, ES, PT, DE, GR, Arabic.



32800



32801

AND PROFESSIONAL MERCURY FREE SPHYGMOMANOMETERS

Professional highly advanced digital blood pressure monitor

- auscultatory blood pressure measurement
- digital readout for pulse rate
- cuff holder and grip for carrying
- durable and chemical resistant body
- full ABS housing lightweight body
- adjustable display column for easier analysis of results
- anti-glare display minimizes reflected light



32775

32776

- **32800 GREEN MERCURY FREE SPHYGMO - desk**
- **32801 GREEN MERCURY FREE SPHYGMO - on trolley with basket**
- **32825 PEDIATRIC CUFF 17-22 cm - optional**
- **32827 ADULT CUFF 22-32 cm - spare**

- **32775 UM-102A MERCURY FREE SPHYGMOMANOMETER - desk**
- **32776 UM-102B MERCURY FREE SPHYGMOMANOMETER - trolley**

TECHNICAL SPECIFICATIONS (32800, 32801)
Size: 360x96x66 mm (32800) Power supply: battery
Measuring range: 20-280 mmHg Accuracy: pressure: ± 3 mmHg

TECHNICAL SPECIFICATIONS
Size: 360x96x66 mm (32775)
Measuring range - pressure: 0-300 mmHg ± 3 mmHg
- pulse: 40-180 beats/minute $\pm 5\%$ (32775-6)
Power supply: 2x1.5V AA alkaline batteries