

## RESPIRA OXYGEN CONCENTRATOR 5 LITRES

### • 34582 RESPIRA OXYGEN CONCENTRATOR 5 I

Small and portable oxygen concentrator composed of air compressor, absorption tower of molecular sieve, sensor and flowrate controller. The machine adopts PSA (pressure swing absorbers) method to directly concentrate medical oxygen with high concentrator from the air. Lightweight, easy to use and safe. Stable performance and low noise level. Wide application in clinic, hospital, oxygen bar, gym center...

Supplied with humidifier bottle, power cord, nasal cannula and user manual (GB, FR\*, IT, ES\*, DE\*) \*only on [www.gimaitaly.com](http://www.gimaitaly.com)

Size: 508x260x h 530 mm. Weight: 21 kg

#### • 34583 PAPER FILTER - box of 10

#### • 34584 HUMIDIFIER BOTTLE - 1 I

#### • 34166 OXYGEN MASK with tube - adult


#### • 34167 OXYGEN MASK with tube - pediatric

#### • 34172 OXY GLASSES - nasal cannulas

#### TECHNICAL SPECIFICATIONS

Max. flowrate:	5l/min
Flow range:	0.5~5l/min
Oxygen concentration:	93% ±3%
Output pressure:	20kPa~50kPa
Input power:	450VA
Noise level:	≤55dB(A)
Operating voltage:	AC220V ± 10%, 50Hz ±1 Hz
Type of protection:	Class II, Type BF, IP21

Oxygen concentrator control panel

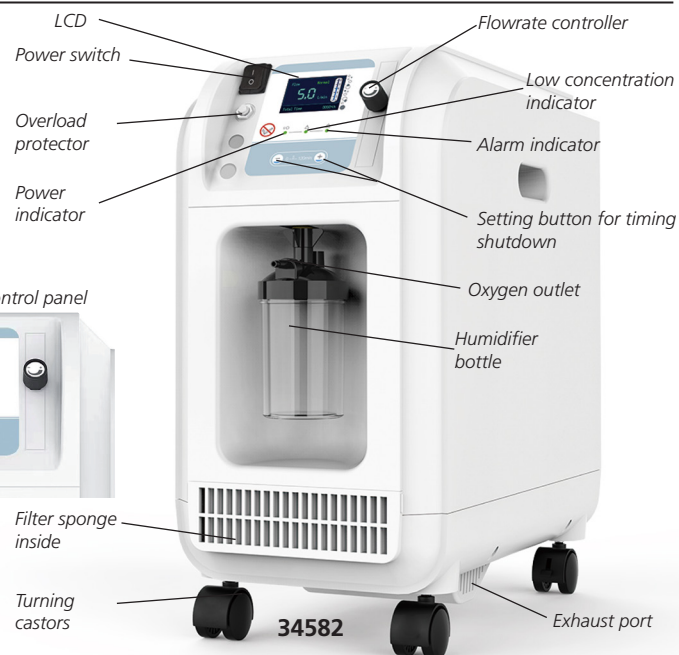

## SMART OXYGEN CONCENTRATORS 3 AND 5 LITRES

### • 34603 SMART OXYGEN CONCENTRATOR 3 I

### • 34598 SMART OXYGEN CONCENTRATOR 5 I

Medical grade, silent, small oxygen concentrators for home use. Provide high concentration of oxygen through nasal cannula. Work continuously. Supplied with filter and humidifier bottle. Class II, Type BF

#### • 34599 PAPER FILTER - box of 10

#### • 34608 HUMIDIFIER BOTTLE - spare

#### • 34167 OXYGEN MASK with tube - pediatric

#### • 34166 OXYGEN MASK with tube - adult

#### • 34172 OXY GLASSES - nasal cannulas

#### TECHNICAL SPECIFICATIONS

	34603	34598
Flow range:	0.5 ~ 3 l/min	0.5 ~ 5 l/min
Oxygen concentration:	93% ± 3%	95.5 ~ 87 %
Output pressure:	20 ~ 50 kPa	40 ~ 70 kPa
Input power:	230 W	350 W
Noise level (front - whole unit):	42 - 46 dBA	≤48 - ≤52 dBA
Size - weight:	39x24.5xh 50 cm - 11.8 kg	39x24xh 50 cm - 16 kg
Operating voltage:	230 V - 50 Hz	230 V - 50 Hz
	(110 V - 60 Hz on request)	
Minimum operating time*	30 min	

\*Time necessary to reach optimal oxygen concentration

**yuwell**  
Top manufacturer  
over 200,000 pcs/years



## OXYGEN CONCENTRATOR 5 LITRES

### • 34605 OXYGEN CONCENTRATOR 5 I

Bionic and streamlined fashion design, which won the national outlook top award-Red Star Award.

Low noise achieved by using advanced round airway and multi-level noise decreasing design.

Large LCD screen for clearer figures.

Cumulative time metering and timing off function.

Fault self-detection system.

Four alarm functions:

- power failure alarm
- pressure cycle failure alarm
- low oxygen concentration alarm
- compressor failure alarm

Supplied with filter, nasal oxygen cannula, nebulizer set and manual (GB, IT)

#### • 34607 PAPER FILTER - box of 10

#### (5 internal + 5 external)

#### • 34619 HUMIDIFIER BOTTLE

#### • 34166 OXYGEN MASK with tube - adult

#### • 34167 OXYGEN MASK with tube - pediatric

#### • 34172 OXY GLASSES - nasal cannulas

#### TECHNICAL SPECIFICATIONS

Flow range:	0.5 - 5 l/min
Oxygen concentration:	95.5% ~ 87%
Output pressure:	40 ~ 70 kPa
Operation noise level:	≤ 53 dB (A)
Input power:	400 VA
Net/gross weight:	25/27 kg
Size:	28.5x54.5xh71.2 cm
Packing size:	62x37xh 71 cm
Operating voltage:	230 V - 50 Hz
	(on request 110 V - 60 Hz)

**yuwell**  
Top manufacturer  
over 200,000 pcs/years

CE 0123

