

VACUUM MATTRESS STRETCHERS

- **34102 VACUUM MATTRESS STRETCHER**
- **34075 VACUUM MATTRESS STRETCHER PLUS**

Vacuum mattresses used to transport injured patients, to prevent compression and mitigate the risk of additional injuries.

High quality and durable TPU material, filled with harmless EPS particulate material and shaped to conform to the human body. They are compatible with X-rays and CT scan.

Healthcare professionals can use the air pump to generate vacuum or negative pressure, facilitating easy and rapid adjustment of softness and hardness based on severity of injury.

Lightweight, foldable after deflation and easy to carry.

Both equipped with four handles on both sides and supplied with air pump and carrying bag.

Max. load capacity: 160 kg (34102), 200 kg (34075).

Unfolded size:

204x104x7 cm (34102)

204x100x5 cm (34075)

Weight: 6 kg



34102



34075

SEMIRIGID SPLINT KIT

- **34119 SEMIRIGID SPLINT KIT**

Kit of 3 splints (arm, wrist, leg) in nylon carrying bag.

All splints are made of polyester + PVC, with lightweight stays providing high strength support. Adjustable velcro for easy use.

Multilingual instructions: GB, FR, IT, RO.



34119

SPINAL IMMOBILIZER

- **34032 SPINAL IMMOBILIZER - adult M (150-175 cm height)**
- **34062 SPINAL IMMOBILIZER - adult L (170-185 cm height)**

Designed to secure the patient, the surrounding is suitable for individuals with rib injuries.

Strips, available in different colours, can be quickly fixed.


It is equipped with nylon buckles and safety belts.

It is recommended that when used in conjunction with a neck fixation, the patient should be in sitting position, before being transferred to a PE stretcher.

Supplied with forehead strap for head fixation and soft strap for chin fixation. Max. load capacity: 175 kg.

Weight: 3 kg. Unfolded size: 79x82 cm. Folded size: 84x25x9 cm


Weight: 4.5 kg. Unfolded size: 88x86 cm. Folded size: 88x28x9 cm

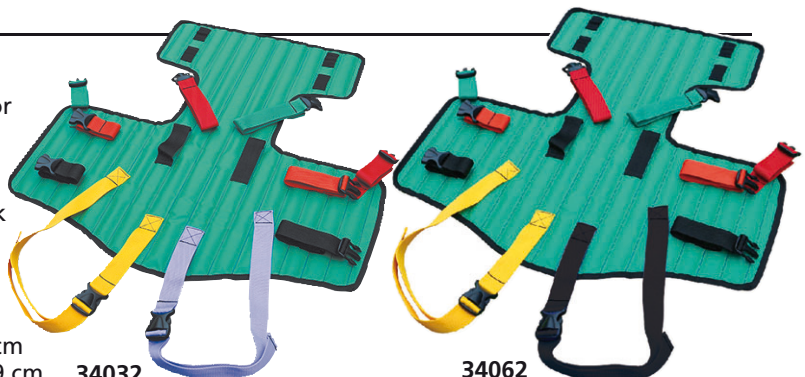

34032

34062

BODY BAGS

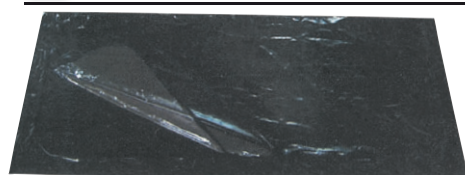

- **34116 BODY BAG - PVC - black**

Made of PVC (Inner transparent PVC, over black), provided with two long slide fasteners.

Size: 210x90 cm

Loading capacity: 150 kg.


- **34124 BODY BAG - PVC - black**

Made of PVC with one "U" zipper and 4 handles.

Size: 210x90 cm;

0.02 cm thickness

Loading capacity: 150 kg

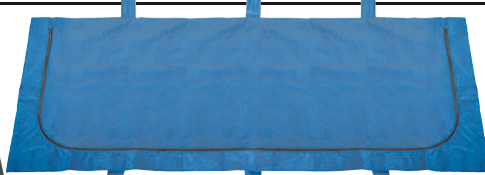

- **34117 BODY BAG - vinyl - blue**

Made of vinyl covered nylon conform to standard EN 1865. Provided with a long slide fastener.

Size: 210x90 cm.

Loading capacity: 150 kg.


- **34118 BODY BAG - vinyl - white**

Made of vinyl with one "U" zipper and 4 handles.

Size: 210x90 cm;

0.02 cm thickness

Loading capacity: 150 kg

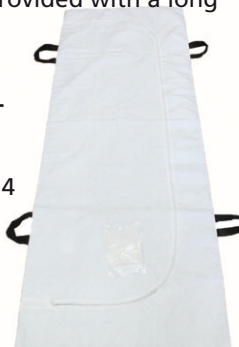

CARRYING COT

34090



With 8 handles for a secure patient transport

- **34090 CARRYING COT**

Made of polyester, PVC coated.

Anti-putridity, bacterial and fungus proof, waterproof, easy to wash and disinfect, X-ray translucent material.

Size 190x75 cm.

Loading capacity: 150 kg.

Certified EN1865-1-2010 + A1 2015.

Supplied in a carrying bag with multilingual manual: GB, FR, IT, ES, PT, DE

Multilingual instructions:
GB, FR, IT, ES, PT, DE, PL, DK,
SE, FI, NO, NL, CZ, SK, HU, HR, SI, RO, BG, GR, LT, LV, EE, Arabic.

Multilingual instructions:
GB, FR, IT, ES, PT, DE, PL, DK,
SE, FI, NO, NL, CZ, SK, HU, HR, SI, RO, BG, GR, LT, LV, EE, Arabic.