

AUTOMATIC LOADING STRETCHER




34055

• 34055 AUTOMATIC LOADING STRETCHER

Automatic stretcher to load and unload patients from emergency vehicles, can be moved by only one person. Specially designed with roll-in system to easily load and unload patients from emergency vehicle without lifting them. Legs fold under cot automatically as stretcher is rolled into the vehicle. Lock-in operating position as cot is removed from vehicle, left-hand and right-hand control levers unlock front and back legs respectively so as to be folded when cot is loaded. Polyethylene covered legs allow cot to roll smoothly in and out of vehicle. Mattress, belts, side rails, IV stand and fixing devices are included.


TECHNICAL SPECIFICATIONS

Maximum size (LxWxH):	193x60x88 cm
Minimum size (LxWxH):	193x60x21 cm
Packing size:	195x61x26 cm
Weight:	40 kg
Load:	159 kg

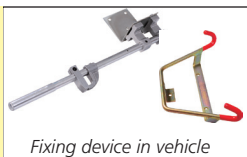

Fixing device in vehicle

MULTIPOSITION STRETCHER

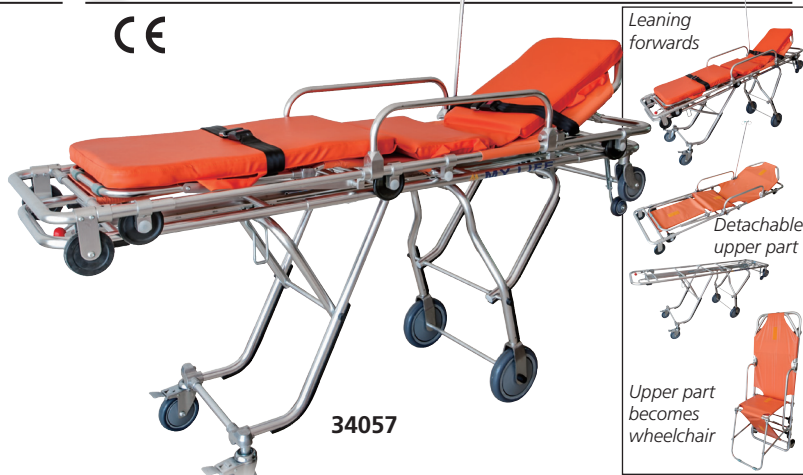


34057

• 34057 MULTIPOSITION AUTO STRETCHER - WHEELCHAIR

Unique design multifunctional automatic stretcher, easily adjustable at different positions for emergency treatment, also in difficult situations. Front and back fixation to lock the stretcher into the ambulance at a desired position. The upper part can become a wheelchair. The stretcher can be adjusted to 6 levels, fold back and forth according to practical needs. When bent forward, the back of the stretcher is supported by fixed pneumatic springs with strapless controlling up to 60 degrees angle. Mattress, belts, side rails, I.V. stand and fixing device are included.

TECHNICAL SPECIFICATIONS

Maximum size (LxWxH):	198x56x89 cm
Minimum size (LxWxH):	198x56x33 cm
Packing size:	199x58x35 cm
Load:	200 kg
Net weight, Gross weight:	40 kg, 46 kg
Maximum angle of leaning forwards backwards:	0-90 degrees
The height scope of touching down (4 grades):	93 cm - 82 cm - 70 cm - 35 cm

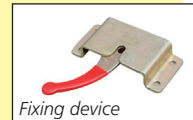

Fixing device

AUTOMATIC MULTIPOSITION STRETCHER



34072

• 34072 FULL AUTOMATIC MULTIPOSITION STRETCHER

Full automatic multiposition stretcher. Easily adjustable to 6 levels for any emergency situation. Backwards and forwards reclining, back and foot rest adjustable to different positions. Can be loaded into ambulance by one operator only and locked at desired position. Structure in high-strength aluminium alloy. Soft foam surface mattress, covered by special seamless welded PVC. Protective block side rust proof and easy to clean. Mattress, belt, side rails, IV stand and fixing device are included.

TECHNICAL SPECIFICATIONS

Maximum size (LxWxH):	198x55x93 cm
Minimum size (LxWxH):	198x55x30 cm
Packing size:	220x60x35 cm
Load:	159 kg
Net weight, Gross weight:	39 kg, 45 kg



WHEELCHAIR STRETCHER




34070

• 34070 WHEELCHAIR STRETCHER

Wheelchair stretcher with varied positions. Lightweight emergency cot quickly convert to any of four chair positions for maneuvering in narrow hallways, elevators and other confined spaces. Stretcher is made of high quality aluminium alloy. Soft foam mattress and protective block side, rust proof and easy to clean. Mattress, belts, side rails, I.V. stand and fixing device are included.

TECHNICAL SPECIFICATIONS

Size (LxWxH):	190x57x85 cm
Packing size:	198x60x25 cm
Load:	159 kg
Net weight, Gross weight:	32 kg, 40 kg

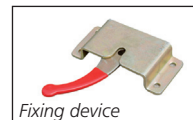

Fixing device

