

VETERINARY RI-SCOPE® OPERATION OTOSCOPE AND SETS

• 80480 RI-SCOPE® OPERATION VET OTOSCOPE HEAD 2.5V - 10541 needs 80488

- swivel lens with 2.5 X fold magnification
- moveable, slotted speculum holder can be locked in any position

- open design for unhindered use during examination and operation

- bayonet fitting for fast and secure attachment to the handle

- all-metal operation otoscope for long service life
- fast exchange of lamp at the front of the otoscope

• 80483 REUSABLE SPECULUM Ø 4 mm - 83 mm - 10881

• 80485 REUSABLE SPECULUM Ø 6 mm - 82 mm - 10882

• 80487 REUSABLE SPECULUM Ø 9.4 mm - 92 mm - 10884

- high quality, made of chrome-plated, corrosion-free metal
- matt inner surfaces reduce the risk of reflections

- slotted model (80487 only): optimum instrument guide provided by special edges on opening

• 80488 HANDLE TYPE AA for 80480 - 10642

Suitable for 2 AA batteries (not included)

RI-SCOPE® L OPHTHALMOSCOPE (only 80492, 80494)

- focusing wheel with corrective lenses from 0 to + 40 and 0 to - 35 diopters

- easy-to-operate aperture hand-wheel with semi-circle, small/medium/large circle fixation star, slit and red-free filter with contrast-enhancing effect

All sets supplied in glass-fiber reinforced and sturdy casing.

80480+
80487+
80488

• 80490 VET HALOGEN OTOSCOPE 2.5V AA handle - 3853

Ri-scope® operation otoscope

Metal specula, unslotted, sterilisable:

- Ø 4 mm, 83 mm long

- Ø 5 mm, 62 mm long

- Ø 6 mm, 82 mm long

Metal specula, slotted, sterilisable:

- Ø 7.5 mm, 68 mm long

- Ø 9.4 mm, 92 mm long

1 spare lamp

swivel lens with 2.5 X magnification



80490

• 80492 VET HALOGEN OTO-OPHTHALMOSCOPE SET 2.5V AA handle - 3861

Ri-scope® operation otoscope, Ri-scope® L1 Ophthalmoscope

Metal specula, unslotted, sterilisable:

- Ø 4 mm, 83 mm long

- Ø 5 mm, 62 mm long

- Ø 6 mm, 82 mm long


Metal specula, slotted, sterilisable:

- Ø 7.5 mm, 68 mm long

- Ø 9.4 mm, 92 mm long

2 spare lamps

swivel lens with 2.5 X magnification

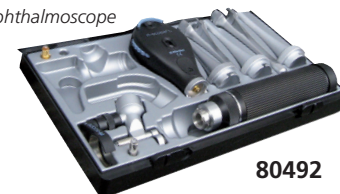

80492

• 80494 VET DELUXE OTO-OPHTHALMOSCOPE SET 2.5V AA handle - 3869

Ri-scope® operation otoscope, Ri-scope® L1 ophthalmoscope

Ri-scope® L tongue depressor, Ri-scope® L nasal speculum

Metal specula, unslotted, sterilisable:

- Ø 4 mm, 83 mm long

- Ø 5 mm, 62 mm long

- Ø 6 mm, 82 mm long

Metal specula, slotted, sterilisable:

- Ø 7.5 mm, 68 mm long

- Ø 9.4 mm, 92 mm long

2 spare lamps

swivel lens with 2.5 X magnification



80494

DIGITAL THERMOMETER



80450

°C and °F display

Battery cover

• 80450 VET DIGITAL THERMOMETER

- flexible tip, waterproof

LCD size 15x7 mm

Multilingual box and

manual: GB, FR, IT, ES, PT, GR.

TECHNICAL SPECIFICATIONS

Measurement range: 32~42.9 °C

Accuracy: ± 0.1 °C

Dimension: 138x21x10 mm

Battery: 1.5V LR41

Measurement time: 60 sec.

VET STETHOSCOPE

• 32545 REGALITE VET - black

Top quality single head stethoscope.

Chromed plated chestpiece with

"SUPERFLEX" inner

spring binaural and 1 piece

moulded PVC,

tubing 90 cm long.

Supplied in deluxe box,

with spare soft eartips

and one spare diaphragm.



90 CM
LONG Y TUBE

CE

32545

DIGITAL THERMOMETER

• 25588 VET RECTAL THERMOMETER

Designed to be ergonomically correct for animals.

Setting for pets, pigs, cattle and other animals.

Features:

- instant temperature in 16 seconds

- beeps indicator

- LCD display and memory

- automatic shut-off after 10 minutes

- display temperature in both Fahrenheit or Celsius scale

- reliable, safe and easy to use

Multilingual box and manual: GB, FR, IT, ES, PT, GR.

16 SECOND MEASUREMENT



25588

TECHNICAL SPECIFICATIONS

Measurement range: 32~42.9°C (89.6~109.2 °F)

Accuracy: ± 0.1 °C (35~39°C)

Battery life is approx. 5,000 readouts

Dimension: 190x110x35 mm

Battery: 2x1.5V AAA

Measurement time: 16 sec.

VISIOFOCUS® VET - BLUETOOTH

• 80457 VISIOFOCUS® VET - Bluetooth

VisioFocus Vet can measure body temperature of animals at a distance in few seconds without touching them.

Thanks to the optical aiming system, temperature is taken from the eyes through a LED system that indicates the correct reading distance. The temperature can also be taken from rectum area, auricle or gum. Built-in cap protects from dust and dirt.

Automatic adjustment to the environments: the software automatically applies a correction factor according to ambient conditions, just leaving it for at least 3-5 minutes in the room. Alternatively you can run a manual calibration (**Manual Quick Calibration System, MQCS**) in only 3 seconds.

Suitable for all animals

Multipurpose: used to take also the temperature of objects, food, liquids or other surfaces.

Bluetooth system to connect with smartphone:

VisioFocus Vet App helps to choose appropriate animal and record a large number of data.

Supplied with 4 AAA batteries.

80457



TECHNICAL SPECIFICATIONS

Body measurement range: 32-43.5°C (89.6-110.3°F)

General measurements range: 1-80°C (33.8-176 °F)

Resolution (body - general): 0.1°C (0.1°F)

Dimension: 144x43.5x21.5 mm (5.66x1.71x0.85") including cap

Battery: 4 AAA - LR03

Battery life is approx. 3 years or 30,000 measurements

Measurement time: immediate

Distance from the animal*: about 6 cm (2.4 inches)

*Calculated by an optical aiming system