

MULTICAREIN MULTIPARAMETER DIAGNOSTIC DEVICE

- 24147 MULTICAREIN* - Italian
- 24148 MULTICAREIN* - English

*Supplied with nylon bag, 5 cholesterol strips, 5 safety lancets, rapid guide (GB, FR, IT, ES, PT, GR, PL, HU) and IT or GB manual.

- 24149 MULTICAREIN STARTER KIT** - GB, IT

- 24153 MULTICAREIN STARTER KIT BLUETOOTH**

**Supplied in cardboard box with a nylon bag, 10 glucose, 5 cholesterol, 5 triglycerides strips, 25 safety lancets, user manual (GB, IT), rapid guide and QR code with complete manual in 8 languages (GB, IT, FR, ES, PT, GR, PL, HU).

Multiparameter diagnostic device for the periodical self-monitoring of risk factors: cholesterol, triglycerides and glucose. Portable device for easy test procedure and quick results.

Accurate with Amperometric Technology for glucose and Reflectometric Technology for cholesterol and triglycerides.

Switchable mg/dL - mmol/L

Multilingual box: GB, FR, IT, ES, GR, PT, PL, HU, RU and selectable internal software: GB, IT.

- 23982 SOFTWARE+CABLE FOR MULTICAREIN - GB, IT, - optional
- Operating system: Windows XP, 7, 8, 10, Vista, 16 Mb Ram.

The first and only multi-parameter device with double technology
AMPEROMETRIC SUCTION STRIP - glucose *Amperometric*

- Electro-chemical technology for a practical use
- Requires very little blood samples (0.5 µl)
- Short time for the results (5 sec.)



REFLECTOMETRIC STRIP - cholesterol and triglycerides
- It allows a more functional self-monitoring of triglycerides and cholesterol.




Tube of 5	Tube of 25	Tube of 50	STRIPS FOR MULTICAREIN*** 9 LANGUAGE BOX - spare
-	24156	24157	Glucose strips biosensor + 1 chip
24158	24159	-	Cholesterol strips + 1 chip
24160	24161	-	Triglycerides strips + 1 chip

Strips are sold in a 9 languages box: GB, IT, FR, ES, PT, PL, GR, HU, Arabic.

GIMA code	CONTROL SOLUTION FOR MULTICAREIN***
23978	Glucose control solution
23979	Cholesterol control solution
23980	Triglycerides control solution
LANCETS FOR MULTICAREIN***	
23917	Lancets 28G - sterile - bag of 25 pcs
23916	Lancets 28G - sterile - box of 100

***All spares are compatible with old Multicare versions

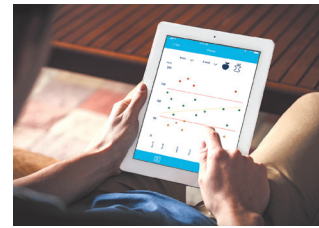

MULTICAREIN



3 IN 1:
CHOLESTEROL,
TRIGLYCERIDES,
GLUCOSE

SELF CONTROL

CE 0344

ISO
15197:2015



24147,
24148



Cholesterol

Glucose

Triglycerides

TECHNICAL SPECIFICATIONS

Measuring range

	Measuring range	Measuring time	Blood sample
Glucose electro-chemical:	10-600 mg/dL - 0.6-33.3 mmol/L	5"	0.5 µl
Cholesterol	130-400 mg/dL - 3-10.2 mmol/L	30"-60"	10 µl
Triglycerides	50-500 mg dL - 0.56 - 5.6 mmol/L	30"-60"	10 µl
Memory:	500 test results with date and time		
Display:	LCD		
Batteries:	2 lithium batteries Cr 2032 3 Volts (32212) suitable for 1 year or 1,000 measurements		
Size - weight:	97x49x20.5 mm - 65 g (battery included)		

MISSION® CHOLESTEROL METER + LIPID PANEL TEST

- 23932 MISSION® CHOLESTEROL METER

Provides high quality and reliable results for complete Lipid Panel testing in ≤2 min.

3 in 1 complete Lipid Panel device tests for Total Cholesterol (CHOL), High Density Lipoprotein (HDL) and Triglycerides (TRIG)

Also includes calculated tests: Low Density Lipoprotein (LDL), CHOL/HDL ratio and Cardiac Risk Assessment

Reliable

- can test whole blood, serum or plasma
- wide hematocrit range of 0-50%
- small printer available for immediate printing of results

Easy to use

- large LCD display for easy reading
- colour coded test devices with engraved test name for easy identification

Supplied with 2 control devices,

4 batteries, carrying case and manual (GB, IT).

Available FR, ES, DE, PT manual on www.gimaitaly.com

Multilingual box: GB, FR, IT, ES, DE, PT, PL, SE, FI, GR, NL, RO, BG, RU, SA.



23932

TECHNICAL SPECIFICATIONS

Detection principle	Reflectance Photometry
Time to results	≤ 2 minutes
Memory with date/time	200 results
Type of test	Quantitative
Tests	CHOL, HDL, TRIG
Calculated tests	LDL, CHOL/HDL ratio, Cardiac Risk Assessment
Specimen volume	3-1 Lipid Panel Device: 35 µL CHOL Individual Test Device: 10 µL
Measurement units	mg/dL or mmol/L
Measurement range	CHOL: 100-500 mg/dL (2.59-12.93 mmol/L) HDL: 15-100 mg/dL (0.39-2.59 mmol/L) TRIG: 45-650 mg/dL (0.51-7.34 mmol/L)
Pc interface for data transfer	Mini USB Port
Calibration	Code Chip (included with 23933-5)
Power source	4 AAA or AC Adaptor (Mini USB, 5V dc, 50 mA)
Battery life	≥ 1,000 tests
Automatic shut off	5 minutes
Meter dimensions - weight	137x79x26 mm (5.4" x 3.1" x 1.0") - 145 g (5.11 oz)
Display dimensions	50x50 mm (2.0" x 2.0")

GIMA code	ACCESSORIES AND SPARES	Minimum order
23935	3 in 1 Lipid panel tests + capil. tubes 35 ul + lancets*	box of 5
23933	3 in 1 Lipid panel tests + capillary tubes 35 ul	box of 25
23934	Total Cholesterol tests + capillary tubes 10 ul	box of 25
24024	Capillary transfer tubes 10 ul - plastic	box of 50
23917	Gima lancets 28G	box of 25
23916	Gima lancets 28G	box of 100
23919	Auto lancet device - AST method	-
23936	Control solution	1 kit
23929	Mission printer**	-
35674	Thermal paper roll 57 mm x 25 m	box of 20

*Safety lancets

**Includes 1 paper roll, cable, adaptor