GIMA

PROFESSIONAL MEDICAL PRODUCTS

TELO PORTAFERITI 190X75 CM
CARRYING COT 190X75 CM
CIVIÈRE DE TRANSPORT 190X75 CM
PORTA-HERIDOS 190X75 CM
MACA DE LONA PARA TRANSPORTAR FERIDOS 190X75 CM
TRAGETUCH FÜR VERLETZTE 190X75 CM

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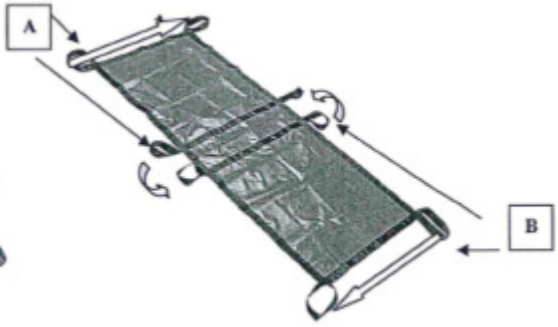
Fig. 2

ENGLISH

INDICATIONS: the flexible carrying cot moves patient securely in restricted areas where standard stretchers have difficulties to be used. **It should not be used if there is a spine, backbone or pelvis suspected injury.**
SPECIFICATIONS: strong polyester fabric, PVC coated, with four secured handles per side to enable a comfortable patient carriage. The materials are waterproof and are easily washable and disinfectable, bacterial and fungus proof, non putrescent. Transparent to X rays.
INSTRUCTIONS: **should be used by skilled staff**. After laying the patient on the carrying cot,
four rescuers (as per fig. 1) should pull up contemporary on the lifting handles and bring the patient to destination. However if four are not available, two rescuers (as per fig. 2).
CLEANING INSTRUCTIONS: after using the carrying sheet, it should be accurately washed with water and soap. It can also be disinfected by a bacterial liquid solution.
STORAGE: no special instructions are requested. Better if kept in a dry place, out of dust and direct sunlight.
MAXIMUM LOAD: 150 kg.
Made according EN 1865-1:2010+A1:2015

GIMA WARRANTY TERMS
The Gima 12-month standard B2B warranty applies.

REF	**IT -** Codice prodotto **GB -** Product code **FR -** Code produit **ES -** Código producto **PT -** Código produto **DE -** Erzeugniscode	LOT	**IT -** Numero di lotto **GB -** Lot number **FR -** Numéro de lot **ES -** Número de lote **PT -** Número de lote **DE -** Chargennummer
	**IT -** Conservare al riparo dalla luce solare **GB -** Keep away from sunlight **FR -** Á conserver à l'abri de la lumière du soleil **ES -** Conservar al amparo de la luz solar **PT -** Guardar ao abrigo da luz solar **DE -** Vor Sonneneinstrahlung geschützt lagern		**IT -** Conservare in luogo fresco ed asciutto **GB -** Keep in a cool, dry place **FR -** Á conserver dans un endroit frais et sec **ES -** Conservar en un lugar fresco y seco **PT -** Armazenar em local fresco e seco **DE -** An einem kühlen und trockenen Ort lagern
CE	**IT -** Dispositivo medico conforme al regolamento (UE) 2017/745 **GB -** Medical Device compliant with Regulation (EU) 2017/745 **FR -** Dispositif médical conforme au règlement (UE) 2017/745 **ES -** Producto sanitario conforme con el reglamento (UE) 2017/745 **PT -** Dispositivo médico em conformidade com a regulamento (UE) 2017/745 **DE -** Medizinprodukt im Sinne der Verordnung (EU) 2017/745		**IT -** Leggere le istruzioni per l'uso **GB -** Consult instructions for use **FR -** Consulter les instructions d'utilisation **ES -** Consultar las instrucciones de uso **PT -** Consulte as instruções de uso **DE -** Gebrauchsanweisung beachten
	**IT -** Fabbricante **GB -** Manufacturer **FR -** Fabricant **ES -** Fabricante **PT -** Fabricante **DE -** Hersteller	MD	**IT -** Dispositivo medico **GB -** Medical Device **FR -** Dispositif médical **ES -** Producto sanitario **PT -** Dispositivo médico **DE -** Medizinprodukt
	**IT -** Data di fabbricazione **GB -** Date of manufacture **FR -** Date de fabrication **ES -** Fecha de fabricación **PT -** Data de fabrico **DE -** Herstellungsdatum		**IT -** Non sterile **GB -** Non-sterile **FR -** Pas stérile **ES -** No estéril **PT -** Não estéril **DE -** Nicht steril
UDI	**IT -** Identificatore univoco del dispositivo **GB -** Unique device identifier **FR -** Identifiant unique de l'appareil **ES -** Identificador de dispositivo único **PT -** Identificador exclusivo do dispositivo **DE -** Eindeutige Kennung des Geräts		**IT -** Attenzione: Leggere e seguire attentamente le istruzioni (avvertenze) per l'uso **GB -** Caution: read instructions (warnings) carefully **FR -** Attention: lisez attentivement les instructions (avertissements) **ES -** Precaución: lea las instrucciones (advertencias) cuidadosamente **PT -** Cuidado: leia as instruções (avisos) cuidadosamente **DE -** Achtung: Anweisungen (Warnungen) sorgfältig lesen
	**IT -** Non utilizzare se l'imballaggio è danneggiato **GB -** Don't use if package is damaged **FR -** Ne pas utiliser si le colis est endommagé **ES -** No usar si el paquete está dañado **PT -** Não use se o pacote estiver danificado **DE -** Nicht verwenden, wenn das Paket beschädigt ist	CH REP	**IT -** Rappresentante Autorizzato in Svizzera **GB -** Authorized Representative in Switzerland **FR -** Représentant autorisé en Suisse **ES -** Representante autorizado en Suiza **PT -** Representante autorizado na Suíça **DE -** Bevollmächtigter Vertreter in der Schweiz

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Gima S.p.A.