


INFUSION SETS

INFUSION SETS

Infusion set includes:
ABS spike with air vent,
6 cm PVC chamber, 150 cm
PVC tube, ABS roller,
ABS luer lock terminal
without needle, latex free injection site.
Sterilized by Ethylene oxide.

- 26062 Y INFUSION SET - polybag - box of 500 (20 bags of 25)
- 26069 HOSE INFUSION SET - polybag - box of 500
- 26064 HOSE INFUSION SET - blister - box of 360 (18 bags of 20)

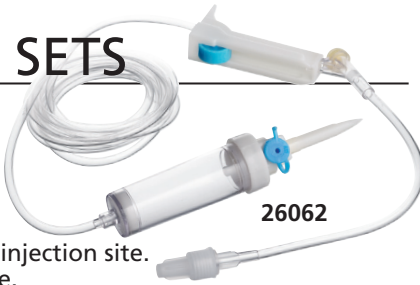

26062

NEEDLE-LESS VALVE

- 26059 NEEDLE-LESS VALVE - sterile - box of 500

Device for I.V. infusion and/or withdrawal,
needle-less closed circuit.

For protection of vascular access and of all
infusion systems equipped with male luer-lock
standard connections.

This closed and self-sealing protection
system allows to infuse or inject fluid,
make withdrawals eliminating the risk
of accidental puncture to operators
and reducing the risk of contamination
to patients.

Medical Device Class I.

Latex free
PVC free

26059



INFUSION SET FOR PUMP

- 26073 ARIES I.V. INFUSION SET - box of 150 (6 bags of 25)

Infusion set for the
administration of
solution into the
human body with
pump.

Tested with infusion pump
Ampall Argus 707 and Optima
PT.



BD SPINAL NEEDLES

Single use latex free sterile BD quincke
point spinal anesthesia needles.



23707

GIMA code	BD code	BD QUINCKE POINT SPINAL ANESTHESIA NEEDLE - STERILE	Minimum order
23705	405248	18G 3.50 IN 1.2x90 mm - pink	box of 25
23707	405253	20G 3.50 IN 0.9x90 mm - yellow	box of 25
23709	405256	22G 3.50 IN 0.7x90 mm - black	box of 25
23710	405240	23G 3.50 IN 0.64x90 mm - turquoise	box of 25
23711	405257	25G 3.50 IN 0.5x90 mm - orange	box of 25

MULTI-VENOUS IV TRAINING ARMS



26048

Length 72.4 cm Weight 1.7 kg



Laerdal
helping save lives

26049

Length 73.7 cm Weight 1.79 kg

- 26048 LAERDAL MALE MULTI-VENOUS IV TRAINING KIT
- 26049 LAERDAL FEMALE MULTI-VENOUS IV TRAINING KIT

Lifelike arm reproduction with multi-vein system
designed for peripheral intravenous therapy.
Kit includes: arm, replacement skin and multi
vein system, simulated blood, blood bag with
tubing and connector, clamp and hook, 5 syringes,
manikin lubricant, carry case and directions for use.

POCKET SHARP CONTAINER

- 25780 POCKET SHARP CONTAINER 0.4 l - 12.3x6.2xh 16.5 cm - box of 70

Single use container for both sharp
and cutting waste disposal.
The permanent and inviolable
closure that prevents waste leakage
and the high container resistance
guarantee top quality and safety.
Equipped with a patented double
opening system, it can be used
by just one hand.

The double opening includes one side
with needles and insulin pen removal,
and the other side with a large
opening for entire syringe disposal.

Double opening system

Temporary closure

Definitive closure



25780

MADE IN ITALY

FEATURES
25780, 25781
Pocket size
Ergonomic shape (25780)
Practical opening system
Needles removal
Insulin pen removal device
Maximum filling level
Shock resistant
Puncture resistant
Sturdy handle
Possible side fastening
Autoclavable
at 134°C for 18'
No emission of any
harmful gas during the
incineration process

DISPO SHARP CONTAINER

- 25781 DISPO SHARP CONTAINER 0.5 l - 8.3x h 15.1 cm - box of 187

Single use conical shape puncture
resistant container for the disposal
of sharp and cutting waste.

The permanent and inviolable closure
that prevents waste leakage and the
high container resistance guarantee
top quality and safety.

Temporary and permanent closures with
a wide opening for the disposal
of larger materials, and three indents for
the needles and insulin pen removal.



Easy lid opening



Temporary closure



Permanent closure

25781

COMPARTMENT TRAY

- 37725 COMPARTMENT TRAY

35x24,5xh 4,4 cm - box of 24
Compartment tray suitable for
carrying small containers and
minor hospital waste.

Tray fits PBS containers 25790-91-92
and CS container 25782, see page 74.



37725

+ 25790

MADE IN ITALY

PRO SPIKE FOR USED NEEDLES

- 27149 PRO SPIKE CONTAINER

27149

Small disposable plastic container for the separate
collection of spiky and cutting wastes (needles,
scalpels) on the operation field. It has a clever lid
with a safe closing system and a clip for belts and
backpacks attachment.

It can be disposed of with the contents.

Size: 90x54x185 mm
Weight: 88 g

