

CHLAMYDIA RAPID TEST


• **29117 CHLAMYDIA TEST - box of 20 tests**
Provide accurate results in 10 minutes, giving doctors immediate diagnosis. High sensitivity and specificity. Each box includes 20 tests, reagents, female cervical swabs, tubes, droppers, pipette and user manual. For professional use.



Male/Female swab specimen
- Add reagent A and B to specimen
- Add 3 drops of extracted specimen
- Read results after 10 minutes




Male urine specimen
- Add reagent B and A to specimen
- Add 3 drops of extracted specimen
- Read results after 10 minutes

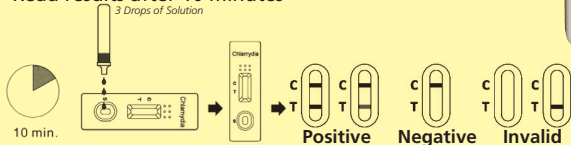

VAGINITIS TEST KIT

The test includes:
1 test + vaginal swab
+ fluid specimen
(buffer solution)
in single box +
English/Italian
manual

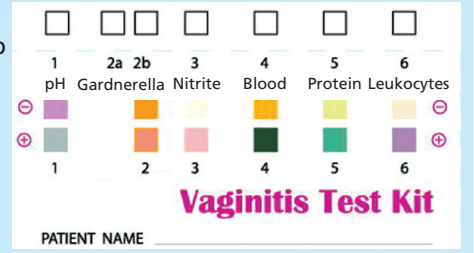


• **29111 VAGINITIS TEST KIT**

A unique, easy-to-use and revolutionary test for gynecologists and operators. Finds the causes of the inflammation/vaginal infection in just 10 minutes. It replaces traditional bacterial tests that take days and costs only a tenth. When necessary dilute vaginal fluid specimen in a buffer solution and then transfer to reagent pad. The instruction manual contains a result guide and can be used to recognize: Trichomoniasis, negative gram bacteria, Gardnerella, Chlamydia, serious bacterial vaginitis. For professional use.

GYN&CUBE PESSARY

• **29887 GYN&CUBE PESSARY Ø 28/38 mm, h 38 mm - standard**
• **29888 GYN&CUBE PESSARY Ø 32/44 mm, h 44 mm - large**
Silicone, autoclavable, latex free pessary with pyramidal shape and draining holes. Kevlar reinforced lanyard. Multilingual instructions: GB, FR, IT, ES, DE.

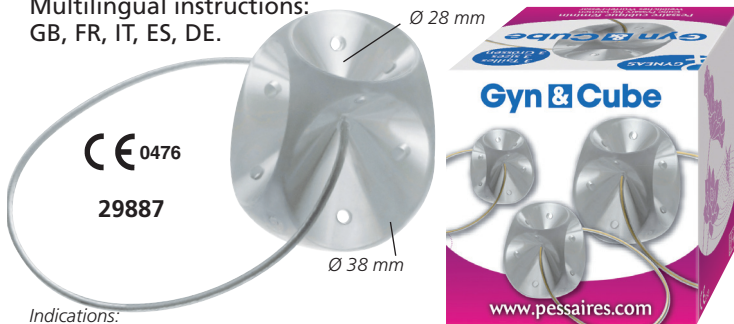


Indications:

- stress urinary incontinence during sustained exertion (aerobics, tennis, jogging)
- moderate or exteriorized prolapse - moderate cystocele and rectocele
- after childbirth, transient prolapse or incontinence.

VAGINAL PH TEST

The kit contains:
10 cotton swabs +
10 reagent strips in
an individual box
with user manual
(English + Italian)

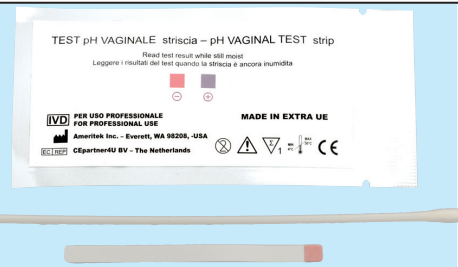

• **29110 VAGINAL PH TEST - box of 10 pcs.**

Quick test to determine vaginal PH.
An easy routine test to check vaginal PH in just 1 minute to make sure it is within the physiological limits (< 4.7). An increase in the Ph level may indicate a change in the vaginal environment. This could mean there is a possible vaginitis. The test is ideal for pregnant women and for women with particular diseases such as diabetes. For professional use.

RUBBER AND SILICONE PESSARY FOR UTERINE PROLAPSE - LATEX FREE



Individually packed in cardboard box. Multilingual box and instructions: GB, FR, IT, ES, PT, DE, GR.



GIMA code	RUBBER PESSARY - autoclavable*
30010	Rubber Pessary Ø 55 mm
30011	Rubber Pessary Ø 60 mm
30012	Rubber Pessary Ø 65 mm
30013	Rubber Pessary Ø 70 mm
30014	Rubber Pessary Ø 75 mm
30015	Rubber Pessary Ø 80 mm
30016	Rubber Pessary Ø 85 mm
30017	Rubber Pessary Ø 90 mm
30018	Rubber Pessary Ø 95 mm
30019	Rubber Pessary Ø 100 mm
SILICONE PESSARY - autoclavable*	
29892	Silicone Pessary Ø 55 mm
29893	Silicone Pessary Ø 60 mm
29894	Silicone Pessary Ø 65 mm
29895	Silicone Pessary Ø 70 mm
29896	Silicone Pessary Ø 75 mm
29897	Silicone Pessary Ø 80 mm
29898	Silicone Pessary Ø 85 mm
29899	Silicone Pessary Ø 90 mm
29900	Silicone Pessary Ø 95 mm
29901	Silicone Pessary Ø 100 mm
BLUE SILICONE PESSARY** - autoclavable*	
29904	Blue silicone Pessary Ø 65 mm
29905	Blue silicone Pessary Ø 70 mm
29906	Blue silicone Pessary Ø 75 mm
29907	Blue silicone Pessary Ø 80 mm
29908	Blue silicone Pessary Ø 85 mm
29909	Blue silicone Pessary Ø 90 mm

*Autoclavable at 121°C/15 min or 134°C/10 min for up to 3 times
**Sterile