

■ AMBU® NEUROLINE™ CONCENTRIC - SINGLE USE **Ambu**

Designed to obtain and transfer signals reliably during EMG procedures. The needle tip cut is designed with a small angle to reduce the penetration force and reduce pain for the patient during needle insertion.

The cannula is made of stainless steel, to provide balance between stiffness and flexibility.

- Silver sensor for suitable impedance level
- stable impedance level throughout the examination
- silicone coating reduces friction when inserting and moving the needle in the target area

All needles feature a colour-coded hub with a tactile index mark to indicate the recording position of the bevel.



AMBU® NEUROLINE™ CONCENTRIC NEUROLOGY NEEDLES					
GIMA code	Ambu® code	Length mm	Ø diam.	Gauge	Minimum order
37520	74025-30/25	25	0.30	30G	box of 25
37523	74038-45/25	38	0.45	26G	box of 25
37524	74050-45/25	50	0.45	26G	box of 25

CONCENTRIC CABLE* FOR 37520-37524

37526	1741	Concentric cable - 100 cm	1
-------	------	---------------------------	---

*Needed for the use of neurology needle electrodes

The hub is ribbed and slightly contoured so that it can be solidly held between the thumb and index finger.

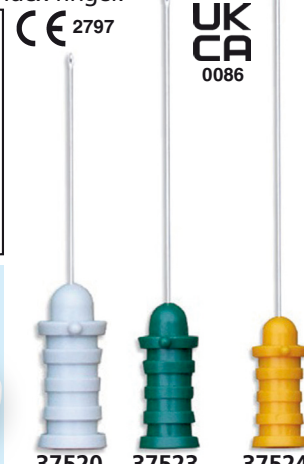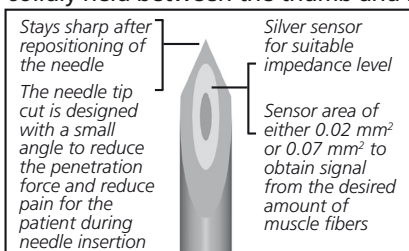

■ AMBU® NEUROLINE™ 700

Small size, multiple purpose neurology electrodes that can be used for all neurological examinations where the electrode is placed outside the hair zone.

Ambu Neuroline solid gel electrode is self-adhesive, and it is easy to reposition, the special "duck foot" design ensures excellent adhesion.

The solid gel electrode is made of a very soft and flexible material, which ensures patient comfort and trace quality even during sleep studies.

Ag/AgCl sensor for best trace quality.

Ideal for examinations of the face, infants and children. Single-patient.



AMBU® NEUROLINE™ 700 ELECTRODES				
GIMA code	Ambu® code	Size mm	Cable length cm	Minimum order
37533	700 10-K/12	20x15	100	box of 120
37536	700 15-K/C/12	20x15	150	box of 120


■ AMBU® NEUROLINE™ 715

Ambu® Neuroline™ 715 is a self-adhesive and repositionable surface electrode. The electrode features a skin friendly solid gel that allows operator to reposition the electrode several times without losing its adhesive power. A transparent backing allows visually guided placement of sensor over the target area. The silver/silver chloride sensor ensures stable signal quality. The electrode is available with a thin and flexible colour-coded lead wire for placement identification. The lead wire has a length of 80 cm and comes with a 1.5 mm touch proof connector. Ideal for Neurophysiological examinations of muscles and/or nerves.



AMBU® NEUROLINE™ 715 ELECTRODES				
GIMA code	Ambu® code	Size mm	Cable length cm	Minimum order
37541	71508-K/C/12	30x22	80	box of 120

■ FIXA DISPOSABLE SKIN STAPLER

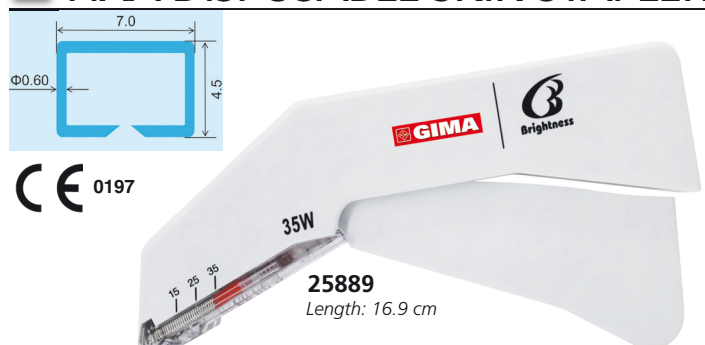

- **25889 FIXA DISPOSABLE SKIN STAPLER - sterile - box of 10**
- Designed for closing external skin under high tension, these skin staplers securely hold tissue in place until sufficient wound healing allows for their removal.
- User-friendly with superior closure performance.
- The device features a clearly marked cartridge scale that allows the operator to monitor the exact number of remaining staples.
- The triangular design of the staple head ensures that skin edges and layers are automatically and neatly aligned during application.
- Each unit is intended for single use only to effectively prevent cross-infection. It includes 35 staples, W size.

NEWS - 10

■ BD Q-SYTE™ NEEDLE-FREE CONNECTOR

- **26058 BD Q-SYTE™ NEEDLE-FREE CONNECTOR - box of 50**
- The BD Q-Syte™ Needle-Free Connector has a split septum design that complies with the guidelines for reducing catheter-related bloodstream infection (CRBSI) rates. Straight, open fluid path design provides a high flow rate. Compact and lightweight design, free of complex mechanical valve systems. Clear housing allows healthcare professionals to see the effectiveness of their own technique and to promptly correct less effective flushing. DEHP free. Made with no natural rubber or latex. Class IIa.

Multilingual manual: GB, FR, IT, ES, PT, DE, NL, DK, SE, NO, FI, LT, LV, EE, PL, CZ, SK, HU, SI, HR, RO, BG, GR, UA, TR, RU.



■ SKIN STAPLE REMOVERS

- **25894 GIMA DISPOSABLE SKIN STAPLE REMOVER - sterile - box of 5**
- Designed to safely remove skin staples by using a pincer-like mechanism. Ergonomic grip. Plastic handle and stainless-steel tips.

