

## PATIENT TRANSFER CHAIRS

- 44840 ISCHIA PATIENT TRANSFER CHAIR - blue
- 44842 ISCHIA PATIENT TRANSFER CHAIR - any colour

Transfer chair for patient transport within hospital or community environment, suited for admission, accident & emergency department (A&E) and ward transfer purposes. It is also well suited for external transfers or when there is a need to jump over obstacles or uneven ground smoothly. Solid front wheels with Ø 30 cm and 2 swivelling rear castors Ø 12.5 cm.

The carer can easily move the chair with patient on board, even in extremely small spaces. With footrest fully touching the ground and backrest up, it enables the patient to safely perform any activity involving sitting or standing up. Armrests can be flipped back to allow side access.

The backrest can be reclined via gas spring control creating a seat tilt (-33°/+6°) that allows for a better position to rest and relax.

Equipped with pushing handle, back pocket and pedal for the activation of "central braking system" and "steering-lock system" on the two rear wheels.

Max weight capacity 160 kg. Class I M fire resistant.




Push handle and adjustable backrest lever

AVAILABLE IN 24 COLOURS


Weight: 42 kg - Seat size: 470x545 mm  
Footrest size: 243x380 mm

Armrest flip-back

Seat and backrest with padded surface for comfort

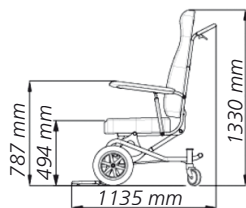

44840

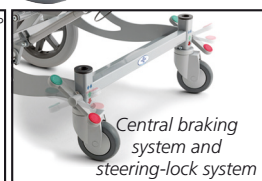

## HYDRAULIC TRANSFER CHAIR

- 43430 HYDRAULIC PATIENT TRANSFER CHAIR

Hydraulic transfer chair made of steel, ideal for transferring the patient from bed to toilet, shower or other sites.

It includes 3 interchangeable seat panels:

- fixed toilet panel that can be used with or without a removable toilet pan
- removable plastic shower panel
- removable soft cushion seat panel

Shower and soft cushion seat panels can be disassembled whenever necessary to accommodate varied conditions, leaving the area more hygienic and cleaner.

- 180° easy opening and closing, flexible for movement
- hydraulic lifting adjustment up to 15 cm
- load capacity: 125 kg

- 4 directional castors (2 frontal with brake)

Packing size: 79x63x33 cm. Seat size: 48x37xh 42-57cm

Multilingual manual: GB, FR, IT, ES, PT, DE, SE, PL, HU, HR, RO, GR.

3 IN 1  
TRANSFER CHAIR  
SHOWER CHAIR  
TOILET CHAIR



43430



## ENT CHAIRS - HYDRAULIC

- 27558 OTO P/I ENT CHAIR - metal blue
- 27559 OTO P/I ENT CHAIR - green Toronto
- 27560 OTO P/I ENT CHAIR - colour on request

ENT patient chair equipped with hydraulic system allowing the lowering and raising of the chair.

The lifting mechanism is operated via a foot control pedal located at the rear of the base. The robust, zinc-coated painted metal structure ensures maximum stability and prevents tipping. The upholstered sections feature high-density, non-deformable foam covered in bi-elastic, tear-resistant, and antimicrobial Skai® synthetic leather, making it exceptionally easy to clean and sanitize. The wide, stable metal base is finished in a sleek black coating.

Adjustable detachable armrests

LIFTING CAPACITY:  
150 kg

Chair weight: 60 kg

Hydraulic pump push: 3000N  
Back rest push: 800N

Fixed headrest



## COMMODE WHEELCHAIR

- 43202 FOLDING COMMODE WHEELCHAIR

- one-button foldable commode
- EVA seat and backrest with soft surface and odour free
- detachable seat and backrest pads
- rust-proof frame, made of aluminum alloy
- equipped with removable toilet pan

FOLDABLE

Packing size: 87x23.5x63 cm

Max. load: 125 kg

Weight: 9.4 kg

Multilingual manual:

GB, FR, IT, ES, PT, DE, SE, PL, HR, RO, GR.



43202